

SEDEO PRO+



© 2015 Handicare

All rights reserved.

The information provided herein may not be reproduced and/or published in any form, by print, photo print, microfilm or any other means whatsoever (electronically or mechanically) without the prior written authorisation of Handicare.

The information provided is based on general data concerning the constructions known at the time of the publication of this manual. Handicare executes a policy of continuous improvement and reserves the right to changes and modifications.

The information provided is valid for the product in its standard version. Handicare cannot be held liable for possible damage resulting from specifications of the product deviating from the standard configuration.

The available information has been prepared with all possible diligence, but Handicare cannot be held liable for possible errors in the information or the consequences thereof.

Handicare accepts no liability for loss resulting from work executed by third parties.

Names, trade names, etc. used by Handicare may not, as per the legislation concerning the protection of trade names, be considered as being available.

Table of content

1	Introduction.....	5
1.1	This manual	5
1.2	Identification of the product	5
1.3	Symbols used in this manual	6
2	Safety.....	7
2.1	Personnel qualification	7
2.2	Cautions and warning statements	7
2.3	Used decals on the seating system	8
3	Tools, parts and components	9
3.1	Tools	9
4	Spare parts.....	10
4.1	Use of the parts lists	10
4.2	Sedeo Pro+	11
1	Seat frame	12
2	Backrest frame	13
3	Seat cushion	14
4	Backrest cushion	23
5	Mounting bracket for legrest	25
6	Standard legrest	26
7	Manually elevating legrest	28
8	Electrically elevating legrest	30
9	Central legrest	32
10	Electrically elevating central legrest	33
11	Flip up armrest	35
12	Headrest	36
13	Joystick bracket	37
14	Lighting	38
15	Abduction support	39
16	Mirror - Side plate	40
17	Armrest bag - Bag hooks - Hemiplegia armrest	41
18	Positioning belt	42
19	Cane holder - Thoracic support	43
20	Worktops power wheelchair	44
5	Service instructions	45
5.1	Maintenance plan / Reimbursement plan	45
5.2	Assembly, replacement and adjustment instructions Sedeo Pro+	46
5.2.1	Adjusting or mounting the seat cushion fixation set	46
5.2.2	Replacing the seat cushion cover	47
5.2.3	Mounting the backrest cushion fixation set	49
5.2.4	Replacing the backrest cushion cover	50
5.2.5	Mounting the thoracic supports	51
5.2.6	Mounting the standard legrest	52
5.2.7	Mounting the manually elevating legrest	54
5.2.8	Mounting the electrically elevating legrest	57
5.2.9	Mounting the central legrest	60
5.2.10	Mounting the electrically central legrest	61
5.2.11	Mounting the standard or swing-away controller bracket	63
5.2.12	Mounting the lights on the seat frame	64
5.2.13	Mounting the lights on the central legrest frame	65
5.2.14	Mounting the side support	66
5.2.15	Mounting the abduction support	67
5.2.16	Mounting the armrest pads	69

5.2.17	Mounting the hemiplegia armrest	70
5.2.18	Mounting the mirror	72
5.2.19	Mounting the positioning belt	74
5.2.20	Mounting the cane holder	75
5.2.21	Mounting the table	76
6	Trouble shooting.....	77
7	Specifications	78
7.1	CE declaration of conformity	78
7.2	Technical information Sedeo Pro+	78
8	Warranty	80
8.1	Definitions of terms	80
8.2	Warranty period table	80

1 Introduction

1.1 This manual

This manual contains the instructions for repairs and general maintenance of the Sedeo Pro+ seating system. Mechanics who do repairs on this seating system must be well trained and familiar with the repair methods and the maintenance of the seating system.

Always make sure that the work is carried out safely, particularly with respect to procedures requiring the seating system to be lifted up.

We advise that you contact our service department before doing repair work on a seating system that has been involved in an accident.

The following specifications are important when ordering parts:

- Model
- Year of manufacture
- Identification number
- Part number
- Name of the part concerned

This information is provided on the identification plate. See 'Identification of the product'.

Available documentation

The following technical documentation is available / required to service this seating system:

- User manual
- Service manual

Service and technical support


For information concerning specific settings, maintenance or repair works please contact your supplier. He is always prepared to help you.

Ensure you have at hand:

- Model
- Year of manufacture
- Identification number

This information is provided on the identification plate. See 'Identification of the product'.

1.2 Identification of the product

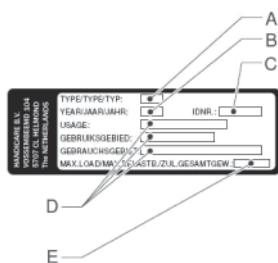

The identification plate contains the following data:

- A. Model
- B. Year of manufacture
- C. Identification number
- D. Use area indoors or outdoors
- E. Maximum load in kg

1.3 Symbols used in this manual

Warning symbol



Follow the instructions next to this symbol closely.

Not paying careful attention to these instructions could result in physical injury or damage to the seating system or the environment.

Reference symbol



The symbol refers to a separate user manual. This reference will indicate the specific user manual and the section to which is being referred.

2 Safety

2.1 Personnel qualification

Service technicians:

Repairs may only be carried out by trained and authorised service technicians.

During the execution of their work, they are at all times fully responsible for the fulfilment of locally applicable safety guidelines and standards.

Temporary employees and persons in training may only carry out repair and replacement work under the supervision of an authorised service technician.

2.2 Cautions and warning statements

⚠ Safety


Safety information is indicated with the warning symbol.

- Follow the instructions carefully next to these warning symbols! Not paying careful attention to these instructions could result in physical injury or damage to the seating system or the environment. Where ever possible, safety information is provided in the relevant chapter.

⚠ Temperature

- Avoid physical contact with the seating system's motors at all times. Motors are continuously in motion during use and can reach high temperatures. After use, the motors will cool down slowly. Physical contact could cause burns.
- If you do not use the seating system, ensure that it is not exposed to direct sunlight for lengthy periods of time. Certain parts of the seating system, such as the seat, the back and the armrests can become hot if they have been exposed to full sunlight for too long. This may cause burns or allergic reactions to the skin.

⚠ Moving parts

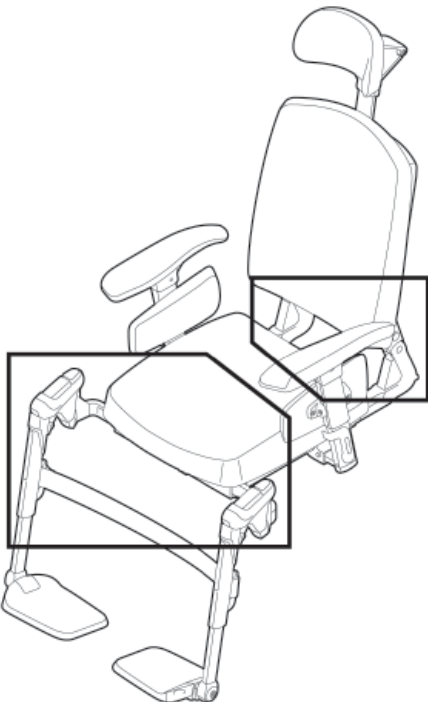

- A seating system has moving and rotating parts. Contact with moving parts may result in serious physical injury or damage to the seating system. Contact with the moving parts of the seating system should be avoided.
- Electrical backrest adjustment
- Electrical elevating legrests

Decals and instructions on the seating system

Decals and instructions on the seating system

- The signs, symbols and instructions affixed to the seating system comprise part of the safety facilities. They must never be covered or removed. They must remain present and clearly legible throughout the entire lifespan of the seating system.
- Replace or repair all illegible or damaged signs, symbols and instructions immediately. Please contact your supplier for assistance.

2.3 Used decals on the seating system

Check manual before using



The seating system must not be lifted up by the armrest assembly, complete and/or legrests.



The stability of the seating system is influenced by the adjustment of the seating system's centre of gravity. Such adjustments should be made by professionals only.

3 Tools, parts and components

3.1 Tools

The tools below are needed for various mechanical settings and maintenance:

Quantity	Description	Size (mm)
1	Screwdriver, medium	-
1	Screwdriver, crosshead	-
1	Hammer (plastic)	-
1	Chaser	-
2	Open ended spanner	17
2	Ring spanner	13
1	Torx key	T25, T30
1	Loctite 243	-
#	Tie wraps	

Remark Size: the EU key width (of the relevant tool).



Use only high-quality tools for the adjustment(s) described.

Make sure that the hexagon of the spanner is attached securely to the hexagon of the fastening article. This prevents the 'rotation' of hexagons of fastening articles and spanners that may affect proper adjustment.

4 Spare parts

4.1 Use of the parts lists

This document is meant as a reference book to be used to order parts for the seating systems that are shown on the front cover.

How to order:

When ordering parts, please specify:

- Serial number (see the identification plate)
- Group (to which the relevant part belongs)
- Article number
- Number of parts required
- Description (in the relevant language)
- Dimensions (if applicable)

Remark:

- If a part does not have a position number, it means that the part concerned cannot be purchased separately. The part concerned is part of the assembly shown. This assembly must be ordered as one piece. It has to be replaced in its entirety.
- Boxed position numbers refer to the relevant drawing.

Order address:

Please mail or fax your orders to your supplier.

Service technicians:

Repairs may only be carried out by trained and authorised service technicians.

During the execution of their work they are at all times fully responsible for the fulfilment of locally applicable safety guidelines and standards.


Temporary employees and persons in training may only carry out repair and replacement work under the supervision of an authorised service technician.

4.2 Sedeo Pro+




Pos	Description	Page
1	Seat frame	12
2	Backrest frame	13
3	Seat cushion	14
4	Backrest cushion	23
5	Mounting bracket for legrest	25
6	Standard legrest	26
7	Manually elevating legrest	28
8	Electrically elevating legrest	30
9	Central legrest	32
10	Electrically elevating central legrest	33
11	Flip up armrest	35
12	Headrest	36
13	Joystick bracket	37
14	Lighting	38
15	Abduction support	39
16	Mirror - Side plate	40
17	Armrest bag - Bag hooks - Hemiplegia armrest	41
18	Positioning belt	42
19	Cane holder - Thoracic support	43
20	Worktops power wheelchair	44

1 Seat frame

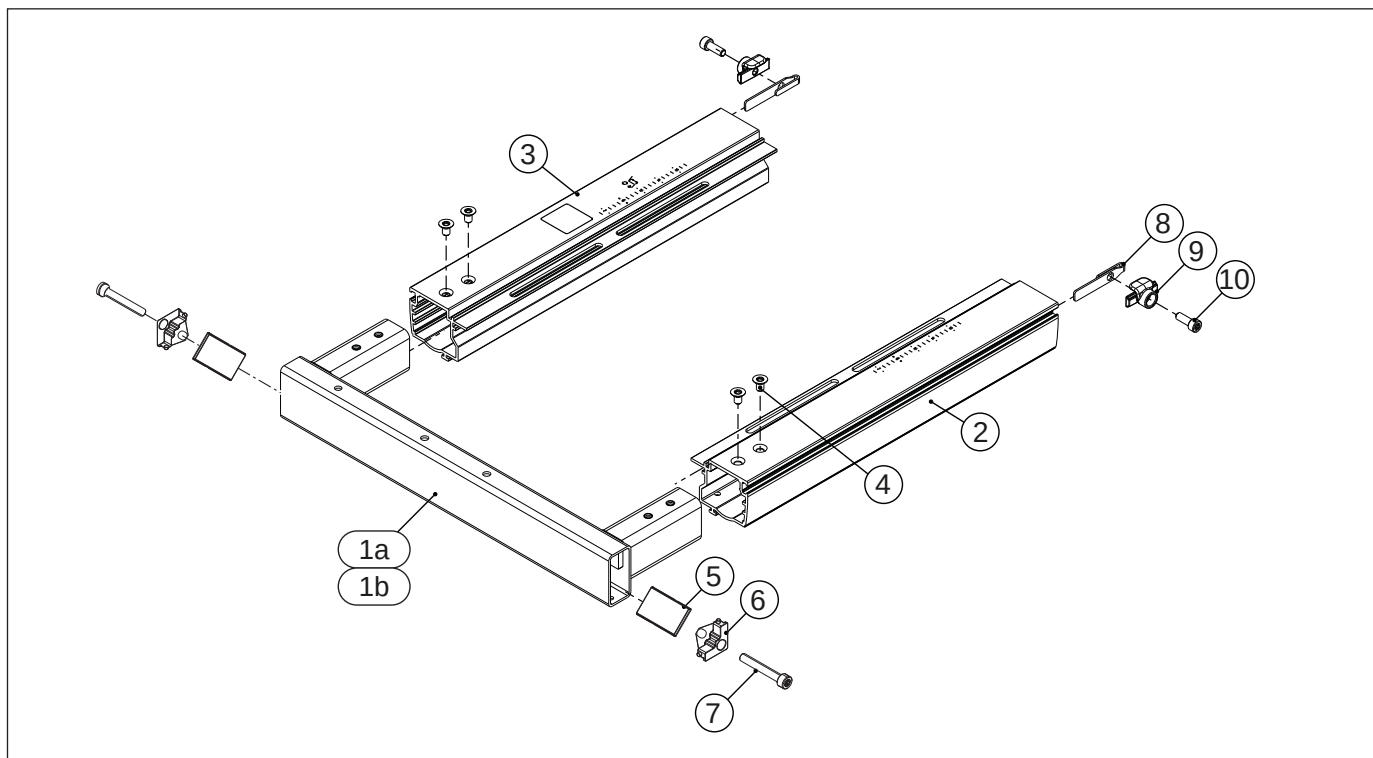


Pos	Qty	Article number	Description
1a	1	9003096	Legrest bracket
1b	-	-	Central legrest bracket, see 'Electrically elevating central legrest'
2	1	9003100	Seat frame profile, left
3	1	9003099	Seat frame profile, right
4	4	1016778	Torx countersunk socket head screw M6x10 iso10642 (100 pieces)
5-7	1	1010873	Legrest clamp, set
7	2	1015338	Torx socket head screw M6x50 iso14579 (100 pieces)
8-10	1	1010876	Seat depth fixation (2 pieces)
10	2	1015333	Torx socket head screw M6x16 iso14579 Precote 85 (100 pieces)

2 Backrest frame

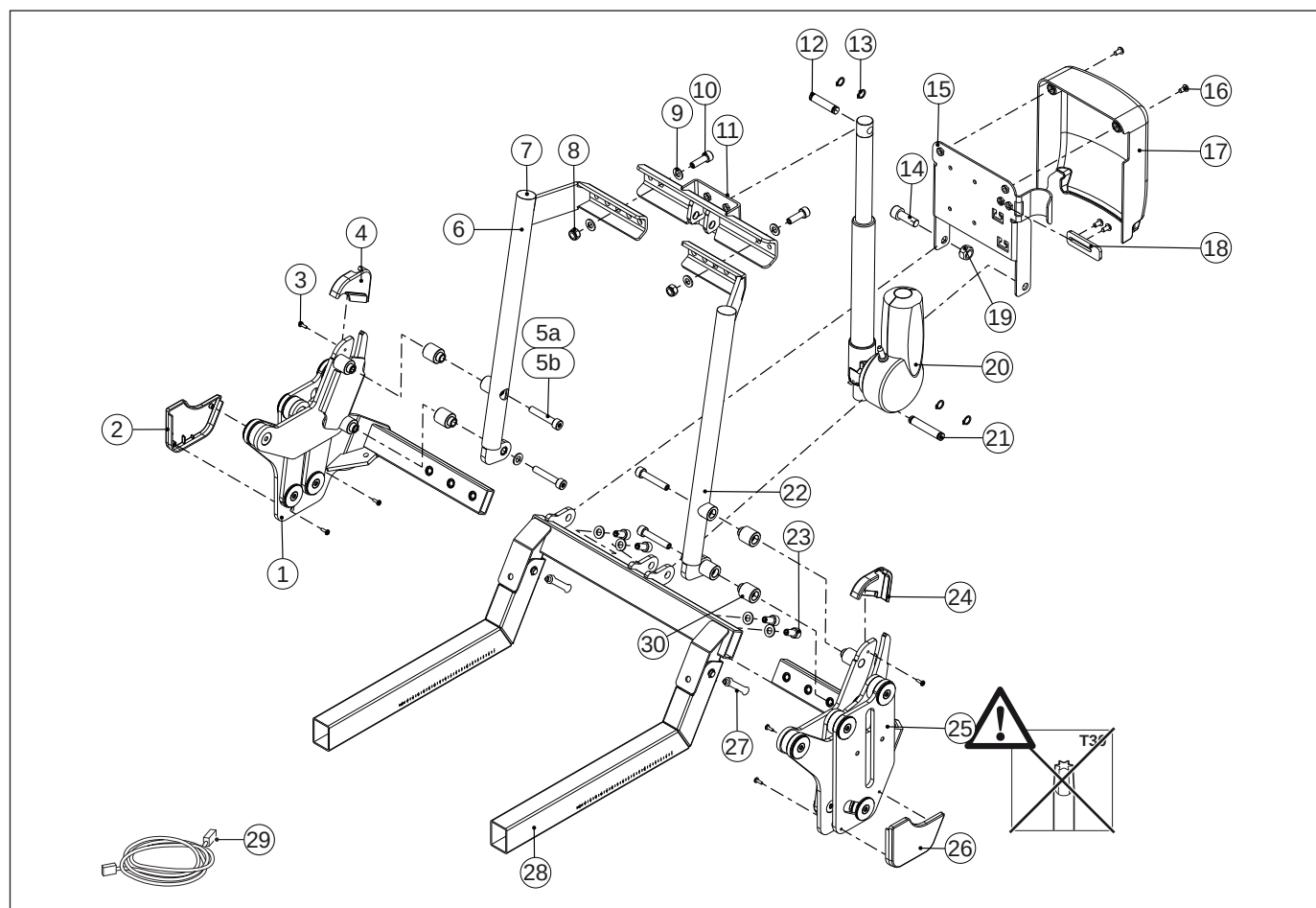

Pos	Qty	Article number	Description
1,10,24	1	1016807	Backrest mechanism, right
2-4	1	1016809	Backrest adjustment covers, right
3	6	-	Screw Ejot, K30x10
3,25,27	1	1016810	Backrest adjustment covers, left
5a	4	1015334	Torx socket head screw M6x25 iso14579 (100 pieces)
5b	4	00000.4032	Socket head screw M8x60 DIN912
5a-12,23	1	1016811	Backrest frame
5b,31	1	1016818	Spacer set seat width 58 cm
8	2	-	Lock nut M8 DIN985/8/VZ
9	10	-	Washer flat M8 DIN125A/VZ
9,23,25	1	1016808	Backrest mechanism, left
10	2	1016851	Torx socket head screw M8x25 iso14579 (100 pieces)
12,13,20,21	1	1016813	Backrest actuator
14	1	-	Socket head screw M10x30
14-19	1	1016814	Electronics box R-net
16	4	-	Buttonhead M5x8 iso7380
16,17	1	1016815	Electronics box cover
19	1	-	Lock nut M10 D985/8/VZ
23	4	1016850	Torx socket head screw M8x16 iso14579 (100 pieces)
27	2	9010103	Clamping lever variable M6x20
27,28	1	1016816	Backrest sliders with cross bar
29	1	1016817	Actuator cable R-net

3 Seat cushion

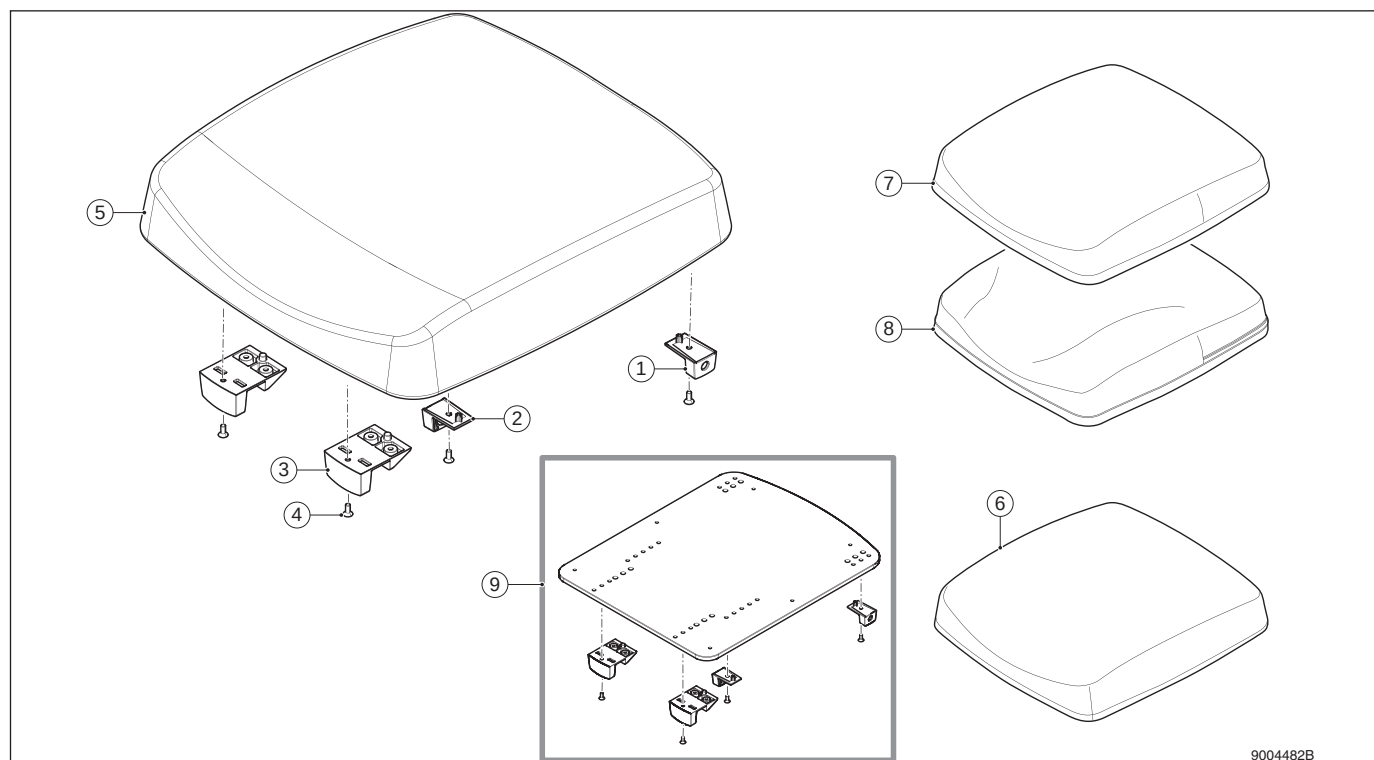

Pos	Qty	Article number	Description
1-4	1	1010953	Seat cushion mounting set
4	4	1015346	Torx socket head screw M5x16 T30 (100 pieces)
9	1	1012267	Interface board Sedeo original
1-6	1	1011726	Seat cushion padded 42x50 Stretch Rubber
1-6	1	1011727	Seat cushion padded 42x54 Stretch Rubber
1-6	1	1011728	Seat cushion padded 42x58 Stretch Rubber
1-6	1	1011729	Seat cushion padded 46x50 Stretch Rubber
1-6	1	1011730	Seat cushion padded 46x54 Stretch Rubber
1-6	1	1011731	Seat cushion padded 46x58 Stretch Rubber
1-6	1	1011732	Seat cushion padded 50x50 Stretch Rubber
1-6	1	1011733	Seat cushion padded 50x54 Stretch Rubber
1-6	1	1011734	Seat cushion padded 50x58 Stretch Rubber
1-6	1	1011735	Seat cushion padded 54x50 Stretch Rubber
1-6	1	1011736	Seat cushion padded 54x54 Stretch Rubber
1-6	1	1011737	Seat cushion padded 54x58 Stretch Rubber
6	1	9004873	Cover seat cushion padded 42x50 Stretch Rubber
6	1	9004874	Cover seat cushion padded 42x54 Stretch Rubber
6	1	9004875	Cover seat cushion padded 42x58 Stretch Rubber
6	1	9004876	Cover seat cushion padded 46x50 Stretch Rubber
6	1	9004877	Cover seat cushion padded 46x54 Stretch Rubber
6	1	9004878	Cover seat cushion padded 46x58 Stretch Rubber
6	1	9004879	Cover seat cushion padded 50x50 Stretch Rubber
6	1	9004880	Cover seat cushion padded 50x54 Stretch Rubber
6	1	9004881	Cover seat cushion padded 50x58 Stretch Rubber
6	1	9004882	Cover seat cushion padded 54x50 Stretch Rubber
6	1	9004883	Cover seat cushion padded 54x54 Stretch Rubber
6	1	9004884	Cover seat cushion padded 54x58 Stretch Rubber
1-5,7,8	1	1011660	Seat cushion padded 42x50 Comfort
1-5,7,8	1	1011661	Seat cushion padded 42x54 Comfort
1-5,7,8	1	1011662	Seat cushion padded 42x58 Comfort
1-5,7,8	1	1011663	Seat cushion padded 46x50 Comfort
1-5,7,8	1	1011664	Seat cushion padded 46x54 Comfort
1-5,7,8	1	1011665	Seat cushion padded 46x58 Comfort

Continued from **Seat cushion**

Pos	Qty	Article number	Description
1-5,7,8	1	1011666	Seat cushion padded 50x50 Comfort
1-5,7,8	1	1011667	Seat cushion padded 50x54 Comfort
1-5,7,8	1	1011668	Seat cushion padded 50x58 Comfort
1-5,7,8	1	1011669	Seat cushion padded 54x50 Comfort
1-5,7,8	1	1011670	Seat cushion padded 54x54 Comfort
1-5,7,8	1	1011671	Seat cushion padded 54x58 Comfort
7-8	1	9004861	Cover seat cushion padded 42x50 Comfort
7-8	1	9004862	Cover seat cushion padded 42x54 Comfort
7-8	1	9004863	Cover seat cushion padded 42x58 Comfort
7-8	1	9004864	Cover seat cushion padded 46x50 Comfort
7-8	1	9004865	Cover seat cushion padded 46x54 Comfort
7-8	1	9004866	Cover seat cushion padded 46x58 Comfort
7-8	1	9004867	Cover seat cushion padded 50x50 Comfort
7-8	1	9004868	Cover seat cushion padded 50x54 Comfort
7-8	1	9004869	Cover seat cushion padded 50x58 Comfort
7-8	1	9004870	Cover seat cushion padded 54x50 Comfort
7-8	1	9004871	Cover seat cushion padded 54x54 Comfort
7-8	1	9004872	Cover seat cushion padded 54x58 Comfort
7	1	9005515	Top cover seat cushion padded 42x50 Comfort
7	1	9005516	Top cover seat cushion padded 42x54 Comfort
7	1	9005517	Top cover seat cushion padded 42x58 Comfort
7	1	9005518	Top cover seat cushion padded 46x50 Comfort
7	1	9005519	Top cover seat cushion padded 46x54 Comfort
7	1	9005520	Top cover seat cushion padded 46x58 Comfort
7	1	9005521	Top cover seat cushion padded 50x50 Comfort
7	1	9005522	Top cover seat cushion padded 50x54 Comfort
7	1	9005523	Top cover seat cushion padded 50x58 Comfort
7	1	9005524	Top cover seat cushion padded 54x50 Comfort
7	1	9005525	Top cover seat cushion padded 54x54 Comfort
7	1	9005526	Top cover seat cushion padded 54x58 Comfort
1-6	1	1011738	Seat cushion memory foam 42x50 Stretch Rubber
1-6	1	1011739	Seat cushion memory foam 42x54 Stretch Rubber
1-6	1	1011740	Seat cushion memory foam 42x58 Stretch Rubber
1-6	1	1011741	Seat cushion memory foam 46x50 Stretch Rubber
1-6	1	1011742	Seat cushion memory foam 46x54 Stretch Rubber
1-6	1	1011743	Seat cushion memory foam 46x58 Stretch Rubber
1-6	1	1011744	Seat cushion memory foam 50x50 Stretch Rubber
1-6	1	1011745	Seat cushion memory foam 50x54 Stretch Rubber
1-6	1	1011746	Seat cushion memory foam 50x58 Stretch Rubber
1-6	1	1011747	Seat cushion memory foam 54x50 Stretch Rubber
1-6	1	1011748	Seat cushion memory foam 54x54 Stretch Rubber
1-6	1	1011749	Seat cushion memory foam 54x58 Stretch Rubber
6	1	9004897	Cover seat cushion memory foam 42x50 Stretch Rubber
6	1	9004898	Cover seat cushion memory foam 42x54 Stretch Rubber
6	1	9004899	Cover seat cushion memory foam 42x58 Stretch Rubber
6	1	9004900	Cover seat cushion memory foam 46x50 Stretch Rubber
6	1	9004901	Cover seat cushion memory foam 46x54 Stretch Rubber
6	1	9004902	Cover seat cushion memory foam 46x58 Stretch Rubber
6	1	9004903	Cover seat cushion memory foam 50x50 Stretch Rubber
6	1	9004904	Cover seat cushion memory foam 50x54 Stretch Rubber
6	1	9004905	Cover seat cushion memory foam 50x58 Stretch Rubber
6	1	9004906	Cover seat cushion memory foam 54x50 Stretch Rubber
6	1	9004907	Cover seat cushion memory foam 54x54 Stretch Rubber
6	1	9004908	Cover seat cushion memory foam 54x58 Stretch Rubber
1-5,7,8	1	1011672	Seat cushion memory foam 42x50 Comfort
1-5,7,8	1	1011674	Seat cushion memory foam 42x54 Comfort
1-5,7,8	1	1011675	Seat cushion memory foam 42x58 Comfort
1-5,7,8	1	1011676	Seat cushion memory foam 46x50 Comfort

Continued from **Seat cushion**

Pos	Qty	Article number	Description
1-5,7,8	1	1011677	Seat cushion memory foam 46x54 Comfort
1-5,7,8	1	1011678	Seat cushion memory foam 46x58 Comfort
1-5,7,8	1	1011679	Seat cushion memory foam 50x50 Comfort
1-5,7,8	1	1011680	Seat cushion memory foam 50x54 Comfort
1-5,7,8	1	1011681	Seat cushion memory foam 50x58 Comfort
1-5,7,8	1	1011682	Seat cushion memory foam 54x50 Comfort
1-5,7,8	1	1011683	Seat cushion memory foam 54x54 Comfort
1-5,7,8	1	1011684	Seat cushion memory foam 54x58 Comfort
7-8	1	9004885	Cover seat cushion memory foam 42x50 Comfort
7-8	1	9004886	Cover seat cushion memory foam 42x54 Comfort
7-8	1	9004887	Cover seat cushion memory foam 42x58 Comfort
7-8	1	9004888	Cover seat cushion memory foam 46x50 Comfort
7-8	1	9004889	Cover seat cushion memory foam 46x54 Comfort
7-8	1	9004890	Cover seat cushion memory foam 46x58 Comfort
7-8	1	9004891	Cover seat cushion memory foam 50x50 Comfort
7-8	1	9004892	Cover seat cushion memory foam 50x54 Comfort
7-8	1	9004893	Cover seat cushion memory foam 50x58 Comfort
7-8	1	9004894	Cover seat cushion memory foam 54x50 Comfort
7-8	1	9004895	Cover seat cushion memory foam 54x54 Comfort
7-8	1	9004896	Cover seat cushion memory foam 54x58 Comfort
7	1	9005527	Top cover seat cushion memory foam 42x50 Comfort
7	1	9005528	Top cover seat cushion memory foam 42x54 Comfort
7	1	9005529	Top cover seat cushion memory foam 42x58 Comfort
7	1	9005530	Top cover seat cushion memory foam 46x50 Comfort
7	1	9005531	Top cover seat cushion memory foam 46x54 Comfort
7	1	9005532	Top cover seat cushion memory foam 46x58 Comfort
7	1	9005533	Top cover seat cushion memory foam 50x50 Comfort
7	1	9005534	Top cover seat cushion memory foam 50x54 Comfort
7	1	9005535	Top cover seat cushion memory foam 50x58 Comfort
7	1	9005536	Top cover seat cushion memory foam 54x50 Comfort
7	1	9005537	Top cover seat cushion memory foam 54x54 Comfort
7	1	9005538	Top cover seat cushion memory foam 54x58 Comfort
1-6	1	1011751	Seat cushion support 42x50 Stretch Rubber
1-6	1	1011752	Seat cushion support 42x54 Stretch Rubber
1-6	1	1011753	Seat cushion support 42x58 Stretch Rubber
1-6	1	1011754	Seat cushion support 46x50 Stretch Rubber
1-6	1	1011755	Seat cushion support 46x54 Stretch Rubber
1-6	1	1011756	Seat cushion support 46x58 Stretch Rubber
1-6	1	1011757	Seat cushion support 50x50 Stretch Rubber
1-6	1	1011758	Seat cushion support 50x54 Stretch Rubber
1-6	1	1011759	Seat cushion support 50x58 Stretch Rubber
1-6	1	1011760	Seat cushion support 54x50 Stretch Rubber
1-6	1	1011761	Seat cushion support 54x54 Stretch Rubber
1-6	1	1011762	Seat cushion support 54x58 Stretch Rubber
6	1	9004921	Cover seat cushion support 42x50 Stretch Rubber
6	1	9004922	Cover seat cushion support 42x54 Stretch Rubber
6	1	9004923	Cover seat cushion support 42x58 Stretch Rubber
6	1	9004924	Cover seat cushion support 46x50 Stretch Rubber
6	1	9004925	Cover seat cushion support 46x54 Stretch Rubber
6	1	9004926	Cover seat cushion support 46x58 Stretch Rubber
6	1	9004927	Cover seat cushion support 50x50 Stretch Rubber
6	1	9004928	Cover seat cushion support 50x54 Stretch Rubber
6	1	9004929	Cover seat cushion support 50x58 Stretch Rubber
6	1	9004930	Cover seat cushion support 54x50 Stretch Rubber
6	1	9004931	Cover seat cushion support 54x54 Stretch Rubber
6	1	9004932	Cover seat cushion support 54x58 Stretch Rubber
1-5,7,8	1	1011685	Seat cushion support 42x50 Comfort
1-5,7,8	1	1011686	Seat cushion support 42x54 Comfort

Continued from **Seat cushion**

Pos	Qty	Article number	Description
1-5,7,8	1	1011687	Seat cushion support 42x58 Comfort
1-5,7,8	1	1011688	Seat cushion support 46x50 Comfort
1-5,7,8	1	1011689	Seat cushion support 46x54 Comfort
1-5,7,8	1	1011690	Seat cushion support 46x58 Comfort
1-5,7,8	1	1011691	Seat cushion support 50x50 Comfort
1-5,7,8	1	1011692	Seat cushion support 50x54 Comfort
1-5,7,8	1	1011693	Seat cushion support 50x58 Comfort
1-5,7,8	1	1011694	Seat cushion support 54x50 Comfort
1-5,7,8	1	1011695	Seat cushion support 54x54 Comfort
1-5,7,8	1	1011696	Seat cushion support 54x58 Comfort
7-8	1	9004909	Cover seat cushion support 42x50 Comfort
7-8	1	9004910	Cover seat cushion support 42x54 Comfort
7-8	1	9004911	Cover seat cushion support 42x58 Comfort
7-8	1	9004912	Cover seat cushion support 46x50 Comfort
7-8	1	9004913	Cover seat cushion support 46x54 Comfort
7-8	1	9004914	Cover seat cushion support 46x58 Comfort
7-8	1	9004915	Cover seat cushion support 50x50 Comfort
7-8	1	9004916	Cover seat cushion support 50x54 Comfort
7-8	1	9004917	Cover seat cushion support 50x58 Comfort
7-8	1	9004918	Cover seat cushion support 54x50 Comfort
7-8	1	9004919	Cover seat cushion support 54x54 Comfort
7-8	1	9004920	Cover seat cushion support 54x58 Comfort
7	1	9005539	Top cover seat cushion support 42x50 Comfort
7	1	9005540	Top cover seat cushion support 42x54 Comfort
7	1	9005541	Top cover seat cushion support 42x58 Comfort
7	1	9005542	Top cover seat cushion support 46x50 Comfort
7	1	9005543	Top cover seat cushion support 46x54 Comfort
7	1	9005544	Top cover seat cushion support 46x58 Comfort
7	1	9005545	Top cover seat cushion support 50x50 Comfort
7	1	9005546	Top cover seat cushion support 50x54 Comfort
7	1	9005547	Top cover seat cushion support 50x58 Comfort
7	1	9005548	Top cover seat cushion support 54x50 Comfort
7	1	9005549	Top cover seat cushion support 54x54 Comfort
7	1	9005550	Top cover seat cushion support 54x58 Comfort
1-6	1	1011763	Seat cushion support TR with memory foam 42x50 Stretch Rubber
1-6	1	1011764	Seat cushion support TR with memory foam 42x54 Stretch Rubber
1-6	1	1011765	Seat cushion support TR with memory foam 42x58 Stretch Rubber
1-6	1	1011766	Seat cushion support TR with memory foam 46x50 Stretch Rubber
1-6	1	1011767	Seat cushion support TR with memory foam 46x54 Stretch Rubber
1-6	1	1011768	Seat cushion support TR with memory foam 46x58 Stretch Rubber
1-6	1	1011769	Seat cushion support TR with memory foam 50x50 Stretch Rubber
1-6	1	1011770	Seat cushion support TR with memory foam 50x54 Stretch Rubber
1-6	1	1011771	Seat cushion support TR with memory foam 50x58 Stretch Rubber
1-6	1	1011772	Seat cushion support TR with memory foam 54x50 Stretch Rubber
1-6	1	1011773	Seat cushion support TR with memory foam 54x54 Stretch Rubber
1-6	1	1011774	Seat cushion support TR with memory foam 54x58 Stretch Rubber
6	1	9005007	Cover seat cushion support TR with memory foam 42x50 Stretch Rubber
6	1	9005008	Cover seat cushion support TR with memory foam 42x54 Stretch Rubber
6	1	9005009	Cover seat cushion support TR with memory foam 42x58 Stretch Rubber
6	1	9005010	Cover seat cushion support TR with memory foam 46x50 Stretch Rubber
6	1	9005011	Cover seat cushion support TR with memory foam 46x54 Stretch Rubber
6	1	9005012	Cover seat cushion support TR with memory foam 46x58 Stretch Rubber
6	1	9005013	Cover seat cushion support TR with memory foam 50x50 Stretch Rubber
6	1	9005014	Cover seat cushion support TR with memory foam 50x54 Stretch Rubber
6	1	9005015	Cover seat cushion support TR with memory foam 50x58 Stretch Rubber
6	1	9005016	Cover seat cushion support TR with memory foam 54x50 Stretch Rubber
6	1	9005017	Cover seat cushion support TR with memory foam 54x54 Stretch Rubber
6	1	9005018	Cover seat cushion support TR with memory foam 54x58 Stretch Rubber

Continued from **Seat cushion**

Pos	Qty	Article number	Description
1-5,7,8	1	1011697	Seat cushion support TR with memory foam 42x50 Comfort
1-5,7,8	1	1011698	Seat cushion support TR with memory foam 42x54 Comfort
1-5,7,8	1	1011699	Seat cushion support TR with memory foam 42x58 Comfort
1-5,7,8	1	1011700	Seat cushion support TR with memory foam 46x50 Comfort
1-5,7,8	1	1011701	Seat cushion support TR with memory foam 46x54 Comfort
1-5,7,8	1	1011702	Seat cushion support TR with memory foam 46x58 Comfort
1-5,7,8	1	1011703	Seat cushion support TR with memory foam 50x50 Comfort
1-5,7,8	1	1011704	Seat cushion support TR with memory foam 50x54 Comfort
1-5,7,8	1	1011705	Seat cushion support TR with memory foam 50x58 Comfort
1-5,7,8	1	1011706	Seat cushion support TR with memory foam 54x50 Comfort
1-5,7,8	1	1011707	Seat cushion support TR with memory foam 54x54 Comfort
1-5,7,8	1	1011750	Seat cushion support TR with memory foam 54x58 Comfort
7-8	1	9004995	Cover seat cushion support TR with memory foam 42x50 Comfort
7-8	1	9004996	Cover seat cushion support TR with memory foam 42x54 Comfort
7-8	1	9004997	Cover seat cushion support TR with memory foam 42x58 Comfort
7-8	1	9004998	Cover seat cushion support TR with memory foam 46x50 Comfort
7-8	1	9004999	Cover seat cushion support TR with memory foam 46x54 Comfort
7-8	1	9005000	Cover seat cushion support TR with memory foam 46x58 Comfort
7-8	1	9005001	Cover seat cushion support TR with memory foam 50x50 Comfort
7-8	1	9005002	Cover seat cushion support TR with memory foam 50x54 Comfort
7-8	1	9005003	Cover seat cushion support TR with memory foam 50x58 Comfort
7-8	1	9005004	Cover seat cushion support TR with memory foam 54x50 Comfort
7-8	1	9005005	Cover seat cushion support TR with memory foam 54x54 Comfort
7-8	1	9005006	Cover seat cushion support TR with memory foam 54x58 Comfort
7	1	9005551	Top cover seat cushion support TR with memory foam 42x50 Comfort
7	1	9005552	Top cover seat cushion support TR with memory foam 42x54 Comfort
7	1	9005553	Top cover seat cushion support TR with memory foam 42x58 Comfort
7	1	9005554	Top cover seat cushion support TR with memory foam 46x50 Comfort
7	1	9005555	Top cover seat cushion support TR with memory foam 46x54 Comfort
7	1	9005556	Top cover seat cushion support TR with memory foam 46x58 Comfort
7	1	9005557	Top cover seat cushion support TR with memory foam 50x50 Comfort
7	1	9005558	Top cover seat cushion support TR with memory foam 50x54 Comfort
7	1	9005559	Top cover seat cushion support TR with memory foam 50x58 Comfort
7	1	9005560	Top cover seat cushion support TR with memory foam 54x50 Comfort
7	1	9005561	Top cover seat cushion support TR with memory foam 54x54 Comfort
7	1	9005562	Top cover seat cushion support TR with memory foam 54x58 Comfort
1-6	1	9008325	Form seat 42x50 Stretch rubber
1-6	1	9008326	Form seat 42x54 Stretch rubber
1-6	1	9008327	Form seat 42x58 Stretch rubber
1-6	1	9008328	Form seat 46x50 Stretch rubber
1-6	1	9008329	Form seat 46x54 Stretch rubber
1-6	1	9008330	Form seat 46x58 Stretch rubber
1-6	1	9008331	Form seat 50x50 Stretch rubber
1-6	1	9008332	Form seat 50x54 Stretch rubber
1-6	1	9008333	Form seat 50x58 Stretch rubber
1-6	1	9008334	Form seat 54x50 Stretch rubber
1-6	1	9008335	Form seat 54x54 Stretch rubber
1-6	1	9008336	Form seat 54x58 Stretch rubber
6	1	9008138	Cover form seat 42x50 Stretch rubber
6	1	9008139	Cover form seat 42x54 Stretch rubber
6	1	9008140	Cover form seat 42x58 Stretch rubber
6	1	9008141	Cover form seat 46x50 Stretch rubber
6	1	9008142	Cover form seat 46x54 Stretch rubber
6	1	9008143	Cover form seat 46x58 Stretch rubber
6	1	9008144	Cover form seat 50x50 Stretch rubber
6	1	9008145	Cover form seat 50x54 Stretch rubber
6	1	9008146	Cover form seat 50x58 Stretch rubber
6	1	9008147	Cover form seat 54x50 Stretch rubber
6	1	9008148	Cover form seat 54x54 Stretch rubber

Continued from **Seat cushion**

Pos	Qty	Article number	Description
6	1	9008149	Cover form seat 54x58 Stretch rubber
1-6	1	1011775	Seat cushion balance 42x50 Stretch Rubber
1-6	1	1011776	Seat cushion balance 42x54 Stretch Rubber
1-6	1	1011777	Seat cushion balance 42x58 Stretch Rubber
1-6	1	1011778	Seat cushion balance 46x50 Stretch Rubber
1-6	1	1011779	Seat cushion balance 46x54 Stretch Rubber
1-6	1	1011780	Seat cushion balance 46x58 Stretch Rubber
1-6	1	1011781	Seat cushion balance 50x50 Stretch Rubber
1-6	1	1011782	Seat cushion balance 50x54 Stretch Rubber
1-6	1	1011783	Seat cushion balance 50x58 Stretch Rubber
1-6	1	1011784	Seat cushion balance 54x50 Stretch Rubber
1-6	1	1011785	Seat cushion balance 54x54 Stretch Rubber
1-6	1	1011786	Seat cushion balance 54x58 Stretch Rubber
6	1	9004945	Cover seat cushion balance 42x50 Stretch Rubber
6	1	9004946	Cover seat cushion balance 42x54 Stretch Rubber
6	1	9004947	Cover seat cushion balance 42x58 Stretch Rubber
6	1	9004948	Cover seat cushion balance 46x50 Stretch Rubber
6	1	9004949	Cover seat cushion balance 46x54 Stretch Rubber
6	1	9004950	Cover seat cushion balance 46x58 Stretch Rubber
6	1	9004951	Cover seat cushion balance 50x50 Stretch Rubber
6	1	9004952	Cover seat cushion balance 50x54 Stretch Rubber
6	1	9004953	Cover seat cushion balance 50x58 Stretch Rubber
6	1	9004954	Cover seat cushion balance 54x50 Stretch Rubber
6	1	9004955	Cover seat cushion balance 54x54 Stretch Rubber
6	1	9004956	Cover seat cushion balance 54x58 Stretch Rubber
1-5,7,8	1	1011708	Seat cushion balance 42x50 Comfort
1-5,7,8	1	1011709	Seat cushion balance 42x54 Comfort
1-5,7,8	1	1011710	Seat cushion balance 42x58 Comfort
1-5,7,8	1	1011711	Seat cushion balance 46x50 Comfort
1-5,7,8	1	1011712	Seat cushion balance 46x54 Comfort
1-5,7,8	1	1011713	Seat cushion balance 46x58 Comfort
1-5,7,8	1	1011714	Seat cushion balance 50x50 Comfort
1-5,7,8	1	1011715	Seat cushion balance 50x54 Comfort
1-5,7,8	1	1011716	Seat cushion balance 50x58 Comfort
1-5,7,8	1	1011717	Seat cushion balance 54x50 Comfort
1-5,7,8	1	1011718	Seat cushion balance 54x54 Comfort
1-5,7,8	1	1011719	Seat cushion balance 54x58 Comfort
7-8	1	9004933	Cover seat cushion balance 42x50 Comfort
7-8	1	9004934	Cover seat cushion balance 42x54 Comfort
7-8	1	9004935	Cover seat cushion balance 42x58 Comfort
7-8	1	9004936	Cover seat cushion balance 46x50 Comfort
7-8	1	9004937	Cover seat cushion balance 46x54 Comfort
7-8	1	9004938	Cover seat cushion balance 46x58 Comfort
7-8	1	9004939	Cover seat cushion balance 50x50 Comfort
7-8	1	9004940	Cover seat cushion balance 50x54 Comfort
7-8	1	9004941	Cover seat cushion balance 50x58 Comfort
7-8	1	9004942	Cover seat cushion balance 54x50 Comfort
7-8	1	9004943	Cover seat cushion balance 54x54 Comfort
7-8	1	9004944	Cover seat cushion balance 54x58 Comfort
7	1	9005563	Top cover seat cushion balance 42x50 Comfort
7	1	9005564	Top cover seat cushion balance 42x54 Comfort
7	1	9005565	Top cover seat cushion balance 42x58 Comfort
7	1	9005566	Top cover seat cushion balance 46x50 Comfort
7	1	9005567	Top cover seat cushion balance 46x54 Comfort
7	1	9005568	Top cover seat cushion balance 46x58 Comfort
7	1	9005569	Top cover seat cushion balance 50x50 Comfort
7	1	9005570	Top cover seat cushion balance 50x54 Comfort
7	1	9005571	Top cover seat cushion balance 50x58 Comfort

Continued from **Seat cushion**

Pos	Qty	Article number	Description
7	1	9005572	Top cover seat cushion balance 54x50 Comfort
7	1	9005573	Top cover seat cushion balance 54x54 Comfort
7	1	9005574	Top cover seat cushion balance 54x58 Comfort
1-6	1	1011787	Seat cushion fluid 42x50 Stretch Rubber
1-6	1	9008469	Seat cushion fluid 42x54 Stretch Rubber
1-6	1	9008470	Seat cushion fluid 42x58 Stretch Rubber
1-6	1	1011788	Seat cushion fluid 46x50 Stretch Rubber
1-6	1	9008471	Seat cushion fluid 46x54 Stretch Rubber
1-6	1	9008472	Seat cushion fluid 46x58 Stretch Rubber
1-6	1	1011789	Seat cushion fluid 50x50 Stretch Rubber
1-6	1	9008473	Seat cushion fluid 50x54 Stretch Rubber
1-6	1	9008474	Seat cushion fluid 50x58 Stretch Rubber
1-6	1	9008475	Seat cushion fluid 54x50 Stretch Rubber
1-6	1	9008476	Seat cushion fluid 54x54 Stretch Rubber
1-6	1	9008477	Seat cushion fluid 54x58 Stretch Rubber
6	1	9004960	Cover seat cushion fluid 42x50 Stretch Rubber
6	1	9008451	Cover seat cushion fluid 42x54 Stretch Rubber
6	1	9008452	Cover seat cushion fluid 42x58 Stretch Rubber
6	1	9004961	Cover seat cushion fluid 46x50 Stretch Rubber
6	1	9008453	Cover seat cushion fluid 46x54 Stretch Rubber
6	1	9008454	Cover seat cushion fluid 46x58 Stretch Rubber
6	1	9004962	Cover seat cushion fluid 50x50 Stretch Rubber
6	1	9008455	Cover seat cushion fluid 50x54 Stretch Rubber
6	1	9008456	Cover seat cushion fluid 50x58 Stretch Rubber
6	1	9008457	Cover seat cushion fluid 54x50 Stretch Rubber
6	1	9008458	Cover seat cushion fluid 54x54 Stretch Rubber
6	1	9008459	Cover seat cushion fluid 54x58 Stretch Rubber
1-5,7,8	1	1011720	Seat cushion fluid 42x50 Comfort
1-5,7,8	1	9008478	Seat cushion fluid 42x54 Comfort
1-5,7,8	1	9008479	Seat cushion fluid 42x58 Comfort
1-5,7,8	1	1011721	Seat cushion fluid 46x50 Comfort
1-5,7,8	1	9008480	Seat cushion fluid 46x54 Comfort
1-5,7,8	1	9008481	Seat cushion fluid 46x58 Comfort
1-5,7,8	1	1011722	Seat cushion fluid 50x50 Comfort
1-5,7,8	1	9008482	Seat cushion fluid 50x54 Comfort
1-5,7,8	1	9008483	Seat cushion fluid 50x58 Comfort
1-5,7,8	1	9008484	Seat cushion fluid 54x50 Comfort
1-5,7,8	1	9008485	Seat cushion fluid 54x54 Comfort
1-5,7,8	1	9008486	Seat cushion fluid 54x58 Comfort
7-8	1	9004957	Cover seat cushion fluid 42x50 Comfort
7-8	1	9008433	Cover seat cushion fluid 42x54 Comfort
7-8	1	9008434	Cover seat cushion fluid 42x58 Comfort
7-8	1	9004958	Cover seat cushion fluid 46x50 Comfort
7-8	1	9008435	Cover seat cushion fluid 46x54 Comfort
7-8	1	9008436	Cover seat cushion fluid 46x58 Comfort
7-8	1	9004959	Cover seat cushion fluid 50x50 Comfort
7-8	1	9008437	Cover seat cushion fluid 50x54 Comfort
7-8	1	9008438	Cover seat cushion fluid 50x58 Comfort
7-8	1	9008439	Cover seat cushion fluid 54x50 Comfort
7-8	1	9008440	Cover seat cushion fluid 54x54 Comfort
7-8	1	9008441	Cover seat cushion fluid 54x58 Comfort
7	1	9005575	Top cover seat cushion fluid 42x50 Comfort
7	1	9008442	Top cover seat cushion fluid 42x54 Comfort
7	1	9008443	Top cover seat cushion fluid 42x58 Comfort
7	1	9005576	Top cover seat cushion fluid 46x50 Comfort
7	1	9008444	Top cover seat cushion fluid 46x54 Comfort
7	1	9008445	Top cover seat cushion fluid 46x58 Comfort
7	1	9005577	Top cover seat cushion fluid 50x50 Comfort


Continued from **Seat cushion**

Pos	Qty	Article number	Description
7	1	9008446	Top cover seat cushion fluid 50x54 Comfort
7	1	9008447	Top cover seat cushion fluid 50x58 Comfort
7	1	9008448	Top cover seat cushion fluid 54x50 Comfort
7	1	9008449	Top cover seat cushion fluid 54x54 Comfort
7	1	9008450	Top cover seat cushion fluid 54x58 Comfort
1-6	1	1011790	Seat cushion Air contour seat ROHO 42x50 Stretch Rubber
1-6	1	9008393	Seat cushion Air contour seat ROHO 42x54 Stretch Rubber
1-6	1	9008394	Seat cushion Air contour seat ROHO 42x58 Stretch Rubber
1-6	1	1011791	Seat cushion Air contour seat ROHO 46x50 Stretch Rubber
1-6	1	9008395	Seat cushion Air contour seat ROHO 46x54 Stretch Rubber
1-6	1	9008396	Seat cushion Air contour seat ROHO 46x58 Stretch Rubber
1-6	1	1011792	Seat cushion Air contour seat ROHO 50x50 Stretch Rubber
1-6	1	9008397	Seat cushion Air contour seat ROHO 50x54 Stretch Rubber
1-6	1	9008398	Seat cushion Air contour seat ROHO 50x58 Stretch Rubber
1-6	1	9008399	Seat cushion Air contour seat ROHO 54x50 Stretch Rubber
1-6	1	9008400	Seat cushion Air contour seat ROHO 54x54 Stretch Rubber
1-6	1	9008401	Seat cushion Air contour seat ROHO 54x58 Stretch Rubber
6	1	1017330	Cover seat cushion Air contour seat ROHO 42x50 Stretch Rubber
6	1	1017332	Cover seat cushion Air contour seat ROHO 42x54 Stretch Rubber
6	1	1017333	Cover seat cushion Air contour seat ROHO 42x58 Stretch Rubber
6	1	1017338	Cover seat cushion Air contour seat ROHO 46x50 Stretch Rubber
6	1	1017340	Cover seat cushion Air contour seat ROHO 46x54 Stretch Rubber
6	1	1017341	Cover seat cushion Air contour seat ROHO 46x58 Stretch Rubber
6	1	1017346	Cover seat cushion Air contour seat ROHO 50x50 Stretch Rubber
6	1	1017348	Cover seat cushion Air contour seat ROHO 50x54 Stretch Rubber
6	1	1017349	Cover seat cushion Air contour seat ROHO 50x58 Stretch Rubber
6	1	1017354	Cover seat cushion Air contour seat ROHO 54x50 Stretch Rubber
6	1	1017356	Cover seat cushion Air contour seat ROHO 54x54 Stretch Rubber
6	1	1017357	Cover seat cushion Air contour seat ROHO 54x58 Stretch Rubber
1-6	1	1014300	Seat cushion pan 5 cm, 42x50 Sellaskin
1-6	1	1011863	Seat cushion pan 5 cm, 42x54 Sellaskin
1-6	1	1011928	Seat cushion pan 5 cm, 42x58 Sellaskin
1-6	1	1011929	Seat cushion pan 5 cm, 46x50 Sellaskin
1-6	1	1011930	Seat cushion pan 5 cm, 46x54 Sellaskin
1-6	1	1011932	Seat cushion pan 5 cm, 46x58 Sellaskin
1-6	1	1011933	Seat cushion pan 5 cm, 50x50 Sellaskin
1-6	1	1011934	Seat cushion pan 5 cm, 50x54 Sellaskin
1-6	1	1011935	Seat cushion pan 5 cm, 50x58 Sellaskin
1-6	1	1011936	Seat cushion pan 5 cm, 54x50 Sellaskin
1-6	1	1011937	Seat cushion pan 5 cm, 54x54 Sellaskin
1-6	1	1011938	Seat cushion pan 5 cm, 54x58 Sellaskin
1-6	1	9008369	Seat cushion pan 7,5 cm, 42x50 Sellaskin
1-6	1	9008370	Seat cushion pan 7,5 cm, 42x54 Sellaskin
1-6	1	9008371	Seat cushion pan 7,5 cm, 42x58 Sellaskin
1-6	1	9008372	Seat cushion pan 7,5 cm, 46x50 Sellaskin
1-6	1	9008373	Seat cushion pan 7,5 cm, 46x54 Sellaskin
1-6	1	9008374	Seat cushion pan 7,5 cm, 46x58 Sellaskin
1-6	1	9008375	Seat cushion pan 7,5 cm, 50x50 Sellaskin
1-6	1	9008376	Seat cushion pan 7,5 cm, 50x54 Sellaskin
1-6	1	9008377	Seat cushion pan 7,5 cm, 50x58 Sellaskin
1-6	1	9008378	Seat cushion pan 7,5 cm, 54x50 Sellaskin
1-6	1	9008379	Seat cushion pan 7,5 cm, 54x54 Sellaskin
1-6	1	9008380	Seat cushion pan 7,5 cm, 54x58 Sellaskin
1-6	1	9008381	Seat cushion pan 10 cm, 42x50 Sellaskin
1-6	1	9008382	Seat cushion pan 10 cm, 42x54 Sellaskin
1-6	1	9008383	Seat cushion pan 10 cm, 42x58 Sellaskin
1-6	1	9008384	Seat cushion pan 10 cm, 46x50 Sellaskin
1-6	1	9008385	Seat cushion pan 10 cm, 46x54 Sellaskin

Continued from **Seat cushion**

Pos	Qty	Article number	Description
1-6	1	9008386	Seat cushion pan 10 cm, 46x58 Sellaskin
1-6	1	9008387	Seat cushion pan 10 cm, 50x50 Sellaskin
1-6	1	9008388	Seat cushion pan 10 cm, 50x54 Sellaskin
1-6	1	9008389	Seat cushion pan 10 cm, 50x58 Sellaskin
1-6	1	9008390	Seat cushion pan 10 cm, 54x50 Sellaskin
1-6	1	9008391	Seat cushion pan 10 cm, 54x54 Sellaskin
1-6	1	9008392	Seat cushion pan 10 cm, 54x58 Sellaskin
1-6	1	1011952	Seat cushion board with rear pan 42x50 Sellaskin
1-6	1	1011953	Seat cushion board with rear pan 42x54 Sellaskin
1-6	1	1011954	Seat cushion board with rear pan 42x58 Sellaskin
1-6	1	1011955	Seat cushion board with rear pan 46x50 Sellaskin
1-6	1	1011956	Seat cushion board with rear pan 46x54 Sellaskin
1-6	1	1011957	Seat cushion board with rear pan 46x58 Sellaskin
1-6	1	1011958	Seat cushion board with rear pan 50x50 Sellaskin
1-6	1	1011959	Seat cushion board with rear pan 50x54 Sellaskin
1-6	1	1011960	Seat cushion board with rear pan 50x58 Sellaskin
1-6	1	1011961	Seat cushion board with rear pan 54x50 Sellaskin
1-6	1	1011962	Seat cushion board with rear pan 54x54 Sellaskin
1-6	1	1011963	Seat cushion board with rear pan 54x58 Sellaskin
1-6	1	1011939	Seat cushion board 42x50 Sellaskin
1-6	1	1011940	Seat cushion board 42x54 Sellaskin
1-6	1	1011941	Seat cushion board 42x58 Sellaskin
1-6	1	1011942	Seat cushion board 46x50 Sellaskin
1-6	1	1011943	Seat cushion board 46x54 Sellaskin
1-6	1	1011944	Seat cushion board 46x58 Sellaskin
1-6	1	1011945	Seat cushion board 50x50 Sellaskin
1-6	1	1011946	Seat cushion board 50x54 Sellaskin
1-6	1	1011948	Seat cushion board 50x58 Sellaskin
1-6	1	1011949	Seat cushion board 54x50 Sellaskin
1-6	1	1011950	Seat cushion board 54x54 Sellaskin
1-6	1	1011951	Seat cushion board 54x58 Sellaskin
-	-	1017108	Seat cushion SW 58 tailor made, please specify
-	-	9008228	Front pan 5 cm for pan seat, SW 42 cm
-	-	9008229	Front pan 5 cm for pan seat, SW 46 cm
-	-	9008230	Front pan 5 cm for pan seat, SW 50 cm
-	-	9008231	Front pan 5 cm for pan seat, SW 54 cm
-	-	9008260	Front pan 7,5 cm for pan seat, SW 42 cm
-	-	9008261	Front pan 7,5 cm for pan seat, SW 46 cm
-	-	9008262	Front pan 7,5 cm for pan seat, SW 50 cm
-	-	9008263	Front pan 7,5 cm for pan seat, SW 54 cm
-	-	9008268	Front pan 10 cm for pan seat, SW 42 cm
-	-	9008269	Front pan 10 cm for pan seat, SW 46 cm
-	-	9008270	Front pan 10 cm for pan seat, SW 50 cm
-	-	9008271	Front pan 10 cm for pan seat, SW 54 cm
-	-	9007089	Rear pan 5 cm for pan seat
-	-	9008274	Rear pan 7,5 cm for pan seat
-	-	1015093	Sellaskin Sedeo Pro (per meter)
-	-	1015091	Comfort fabric (per meter)
-	-	00004.0291	Black Stretch Rubber (per meter)

4 Backrest cushion

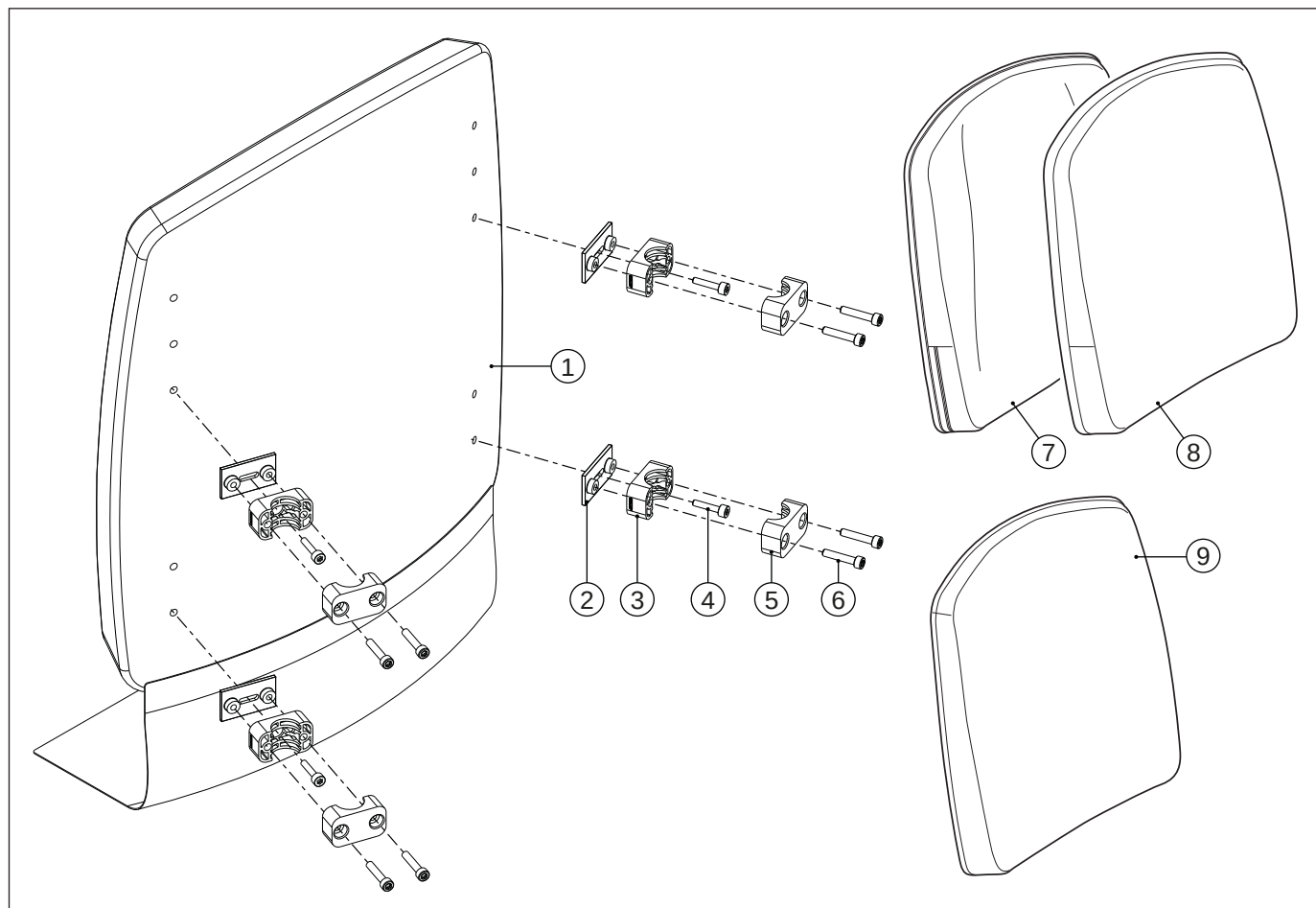


Pos	Qty	Article number	Description
2-6	1	1016821	Backrest mounting, set
4	4	1016852	Torx socket head screw M5x25 iso14579 (100 pieces)
6	4	1016780	Torx socket head screw M6x30 iso14579 (100 pieces)
1-6,9	1	1017058	Backrest flexible with adjustable straps active 42 x 44 Stretch rubber, Sedeo Pro+
1-6,9	1	1017059	Backrest flexible with adjustable straps active 42 x 50 Stretch rubber, Sedeo Pro+
1-6,9	1	1017060	Backrest flexible with adjustable straps active 42 x 56 Stretch rubber, Sedeo Pro+
1-6,9	1	1017061	Backrest flexible with adjustable straps active 46 x 44 Stretch rubber, Sedeo Pro+
1-6,9	1	1017062	Backrest flexible with adjustable straps active 46 x 50 Stretch rubber, Sedeo Pro+
1-6,9	1	1017063	Backrest flexible with adjustable straps active 46 x 56 Stretch rubber, Sedeo Pro+
1-6,9	1	1017064	Backrest flexible with adjustable straps active 50 x 44 Stretch rubber, Sedeo Pro+
1-6,9	1	1017065	Backrest flexible with adjustable straps active 50 x 50 Stretch rubber, Sedeo Pro+
1-6,9	1	1017066	Backrest flexible with adjustable straps active 50 x 56 Stretch rubber, Sedeo Pro+
1-6,9	1	1017067	Backrest flexible with adjustable straps active 54 x 44 Stretch rubber, Sedeo Pro+
1-6,9	1	1017068	Backrest flexible with adjustable straps active 54 x 50 Stretch rubber, Sedeo Pro+
1-6,9	1	1017069	Backrest flexible with adjustable straps active 54 x 56 Stretch rubber, Sedeo Pro+
9	1	1017070	Cover backrest flexible with adjustable straps active 42 x 44 Stretch rubber, Sedeo Pro+
9	1	1017071	Cover backrest flexible with adjustable straps active 42 x 50 Stretch rubber, Sedeo Pro+
9	1	1017072	Cover backrest flexible with adjustable straps active 42 x 56 Stretch rubber, Sedeo Pro+
9	1	1017073	Cover backrest flexible with adjustable straps active 46 x 44 Stretch rubber, Sedeo Pro+
9	1	1017074	Cover backrest flexible with adjustable straps active 46 x 50 Stretch rubber, Sedeo Pro+
9	1	1017075	Cover backrest flexible with adjustable straps active 46 x 56 Stretch rubber, Sedeo Pro+
9	1	1017076	Cover backrest flexible with adjustable straps active 50 x 44 Stretch rubber, Sedeo Pro+
9	1	1017077	Cover backrest flexible with adjustable straps active 50 x 50 Stretch rubber, Sedeo Pro+
9	1	1017078	Cover backrest flexible with adjustable straps active 50 x 56 Stretch rubber, Sedeo Pro+
9	1	1017079	Cover backrest flexible with adjustable straps active 54 x 44 Stretch rubber, Sedeo Pro+
9	1	1017080	Cover backrest flexible with adjustable straps active 54 x 50 Stretch rubber, Sedeo Pro+
9	1	1017081	Cover backrest flexible with adjustable straps active 54 x 56 Stretch rubber, Sedeo Pro+

Continued from **Backrest cushion**


Pos	Qty	Article number	Description
1-8	1	1017082	Backrest flexible with adjustable straps active 42 x 44 Comfort, Sedeo Pro+
1-8	1	1017083	Backrest flexible with adjustable straps active 42 x 50 Comfort, Sedeo Pro+
1-8	1	1017084	Backrest flexible with adjustable straps active 42 x 56 Comfort, Sedeo Pro+
1-8	1	1017085	Backrest flexible with adjustable straps active 46 x 44 Comfort, Sedeo Pro+
1-8	1	1017086	Backrest flexible with adjustable straps active 46 x 50 Comfort, Sedeo Pro+
1-8	1	1017087	Backrest flexible with adjustable straps active 46 x 56 Comfort, Sedeo Pro+
1-8	1	1017088	Backrest flexible with adjustable straps active 50 x 44 Comfort, Sedeo Pro+
1-8	1	1017089	Backrest flexible with adjustable straps active 50 x 50 Comfort, Sedeo Pro+
1-8	1	1017090	Backrest flexible with adjustable straps active 50 x 56 Comfort, Sedeo Pro+
1-8	1	1017091	Backrest flexible with adjustable straps active 54 x 44 Comfort, Sedeo Pro+
1-8	1	1017092	Backrest flexible with adjustable straps active 54 x 50 Comfort, Sedeo Pro+
1-8	1	1017093	Backrest flexible with adjustable straps active 54 x 56 Comfort, Sedeo Pro+
7,8	1	1017094	Cover backrest flexible with adjustable straps active 42 x 44 Comfort, Sedeo Pro+
7,8	1	1017095	Cover backrest flexible with adjustable straps active 42 x 50 Comfort, Sedeo Pro+
7,8	1	1017096	Cover backrest flexible with adjustable straps active 42 x 56 Comfort, Sedeo Pro+
7,8	1	1017097	Cover backrest flexible with adjustable straps active 46 x 44 Comfort, Sedeo Pro+
7,8	1	1017098	Cover backrest flexible with adjustable straps active 46 x 50 Comfort, Sedeo Pro+
7,8	1	1017099	Cover backrest flexible with adjustable straps active 46 x 56 Comfort, Sedeo Pro+
7,8	1	1017100	Cover backrest flexible with adjustable straps active 50 x 44 Comfort, Sedeo Pro+
7,8	1	1017101	Cover backrest flexible with adjustable straps active 50 x 50 Comfort, Sedeo Pro+
7,8	1	1017102	Cover backrest flexible with adjustable straps active 50 x 56 Comfort, Sedeo Pro+
7,8	1	1017103	Cover backrest flexible with adjustable straps active 54 x 44 Comfort, Sedeo Pro+
7,8	1	1017104	Cover backrest flexible with adjustable straps active 54 x 50 Comfort, Sedeo Pro+
7,8	1	1017105	Cover backrest flexible with adjustable straps active 54 x 56 Comfort, Sedeo Pro+
8	1	9008058	Top cover backrest flexible with adjustable straps active 42 x 44 Comfort
8	1	9005429	Top cover backrest flexible with adjustable straps active 42 x 50 Comfort
8	1	9005430	Top cover backrest flexible with adjustable straps active 42 x 56 Comfort
8	1	9008059	Top cover backrest flexible with adjustable straps active 46 x 44 Comfort
8	1	9005431	Top cover backrest flexible with adjustable straps active 46 x 50 Comfort
8	1	9005432	Top cover backrest flexible with adjustable straps active 46 x 56 Comfort
8	1	9008060	Top cover backrest flexible with adjustable straps active 50 x 44 Comfort
8	1	9005433	Top cover backrest flexible with adjustable straps active 50 x 50 Comfort
8	1	9005434	Top cover backrest flexible with adjustable straps active 50 x 56 Comfort
8	1	9008061	Top cover backrest flexible with adjustable straps active 54 x 44 Comfort
8	1	9005435	Top cover backrest flexible with adjustable straps active 54 x 50 Comfort
8	1	9005436	Top cover backrest flexible with adjustable straps active 54 x 56 Comfort
		1017109	Backrest cushion SW 58 tailor made, please specify

5 Mounting bracket for legrest

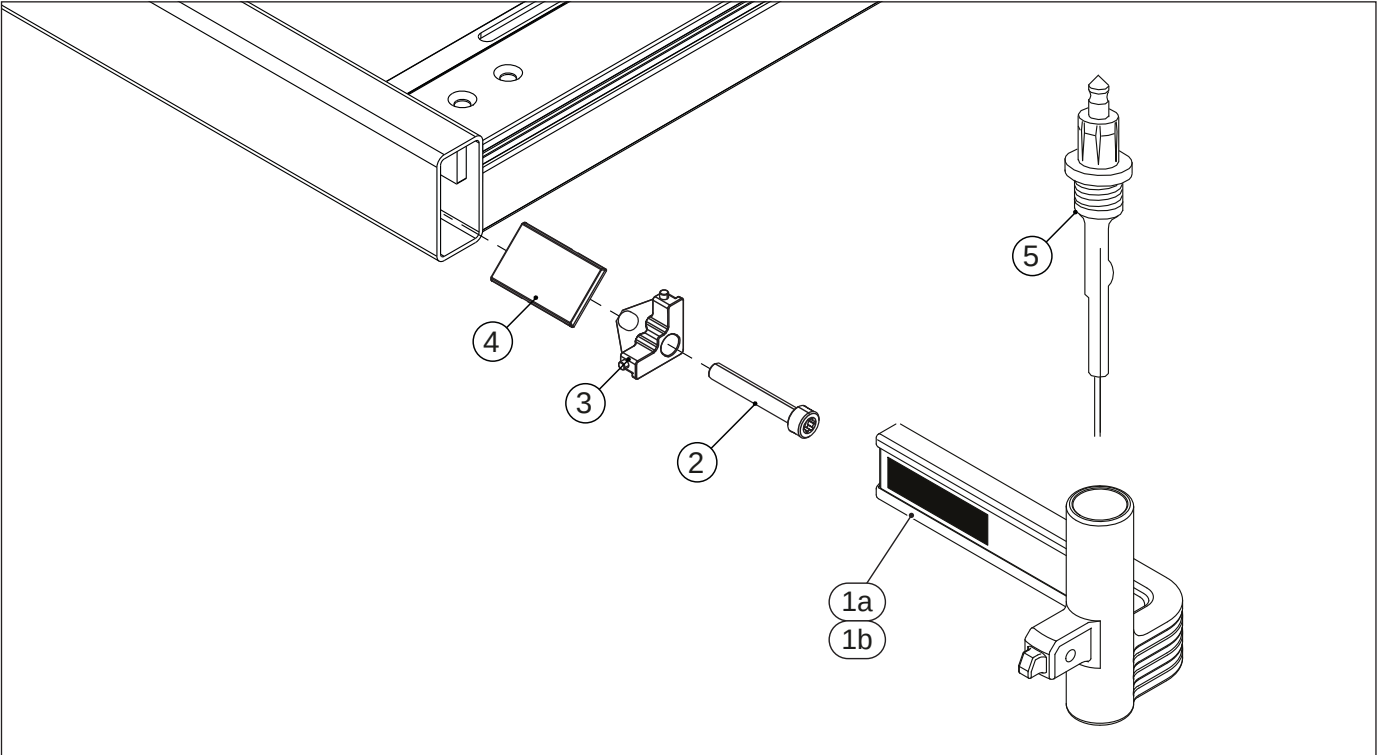

Pos	Qty	Article number	Description
1a	1	9001156	Legrest holder right
1b	1	9001155	Legrest holder left
2	2	1015338	Torx socket head screw M6x50 iso 14579 (100 pieces)
2-4,6	1	1010873	Legrest clamp, set
5	2	9002923	Legrest harness

6 Standard legrest

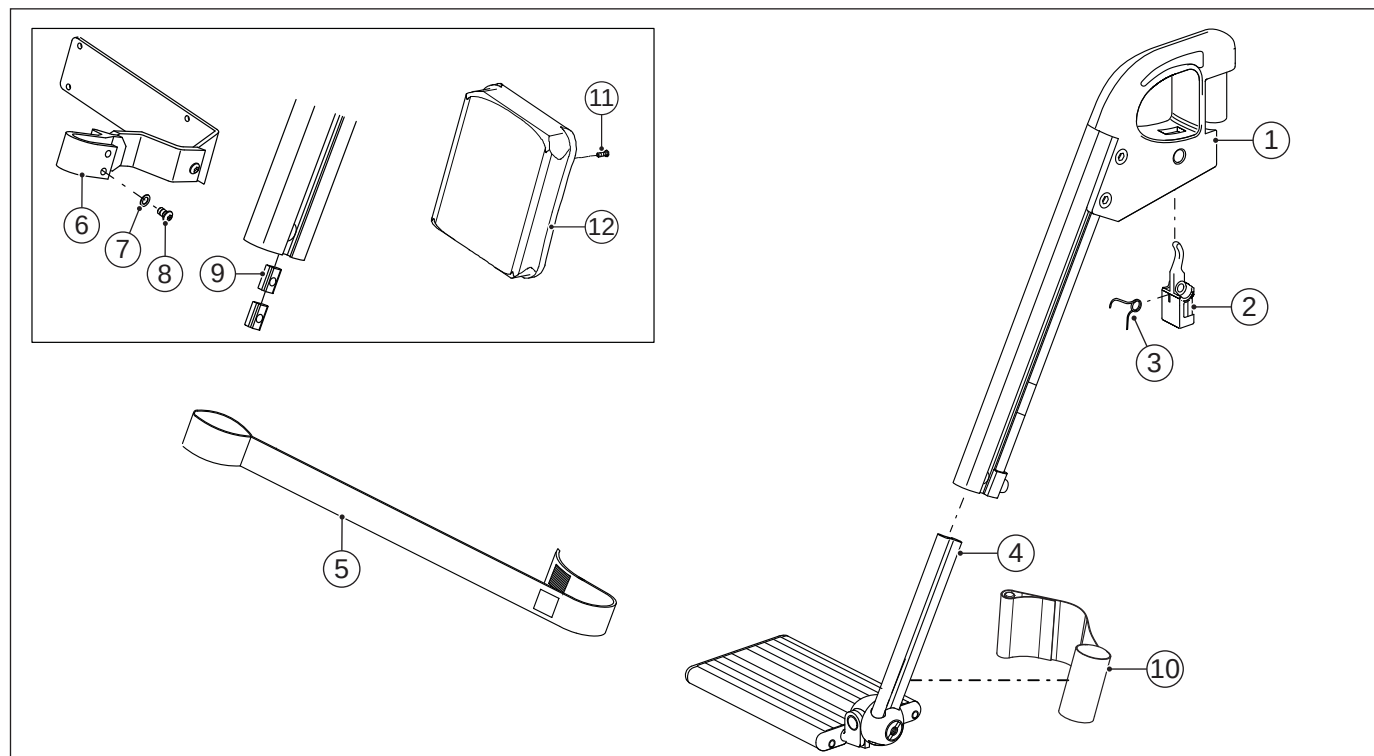


Pos	Qty	Article number	Description
1-3	1	1002514	Sedeo Original standard legrest upper part right
1-3	1	1002515	Sedeo Original standard legrest upper part left
1-4	1	1003563	Sedeo Original standard legrest for SW 42 cm (size 15) right
1-4	1	01500.9019	Sedeo Original standard legrest for SW 46 cm (size 17) right
1-4	1	1003565	Sedeo Original standard legrest for SW 50 cm (size 19) right
1-4	1	1003567	Sedeo Original standard legrest for SW 54 cm (size 21) right
1-4	1	1003562	Sedeo Original standard legrest for SW 42 cm (size 15) left
1-4	1	01500.9119	Sedeo Original standard legrest for SW 46 cm (size 17) left
1-4	1	1003564	Sedeo Original standard legrest for SW 50 cm (size 19) left
1-4	1	1003566	Sedeo Original standard legrest for SW 54 cm (size 21) left
2	2	01500.1112	Legrest release lever
3	2	01500.2212	Legrest release spring
4	1	140.00171.000	Sedeo Original footplate for SW 42 cm (size 15) right
4	1	140.00121.000	Sedeo Original footplate for SW 46 cm (size 17) right
4	1	140.00181.000	Sedeo Original footplate for SW 50 cm (size 19) right
4	1	1003220	Sedeo Original footplate for SW 54 cm (size 21) right
4	1	140.00201.000	Sedeo Original footplate for SW 42 cm (size 15) left
4	1	140.00131.000	Sedeo Original footplate for SW 46 cm (size 17) left
4	1	140.00211.000	Sedeo Original footplate for SW 50 cm (size 19) left
4	1	1003219	Sedeo Original footplate for SW 54 cm (size 21) left
5	1	01263.9000	Calf strap
6-9	2	1004189	Calf plate hinge standard
6-9	2	1004188	Calf plate hinge small
7	4	00000.2601	Toothed spring washer M6
8	4	00000.3725	Button head M6X10
9	4	01500.2110	Legrest clamping strip short
10	1	01273.9110	Canvas heel support left = right
11	4	00000.3206	Screw M5x16
11,12	1	220.00121.505	Calf plate standard
11,12	1	220.00141.505	Calf plate small
-	1	01503.8000	Sedeo Original standard legrest top complete right
-	1	01503.8100	Sedeo Original standard legrest top complete left
-	1	1015137	Sedeo Original standard legrests with calf strap for SW 42 cm, set (mounting brackets included)
-	1	1015138	Sedeo Original standard legrests with calf strap for SW 46 cm, set (mounting brackets included)

Continued from **Standard legrest**

Pos	Qty	Article number	Description
-	1	1015139	Sedeo Original standard legrests with calf strap for SW 50 cm, set (mounting brackets included)
-	1	1015140	Sedeo Original standard legrests with calf strap for SW 54 cm, set (mounting brackets included)
-	1	1015141	Sedeo Original standard legrests with calf support for SW 42 cm, set (mounting brackets included)
-	1	1015142	Sedeo Original standard legrests with calf support for SW 46 cm, set (mounting brackets included)
-	1	1015143	Sedeo Original standard legrests with calf support for SW 50 cm, set (mounting brackets included)
-	1	1015144	Sedeo Original standard legrests with calf support for SW 54 cm, set (mounting brackets included)
-	1	9006539	Torx set Sedeo Original legrest
-	2	1002599	Legrest inner extrusion 10 cm extended
-	2	1005119	Depth adjustable calf plate bracket (max. 4,5 cm)
-	1	1005160	Depth adjustable calf plate bracket (max. 4,5 cm), set

7 Manually elevating legrest

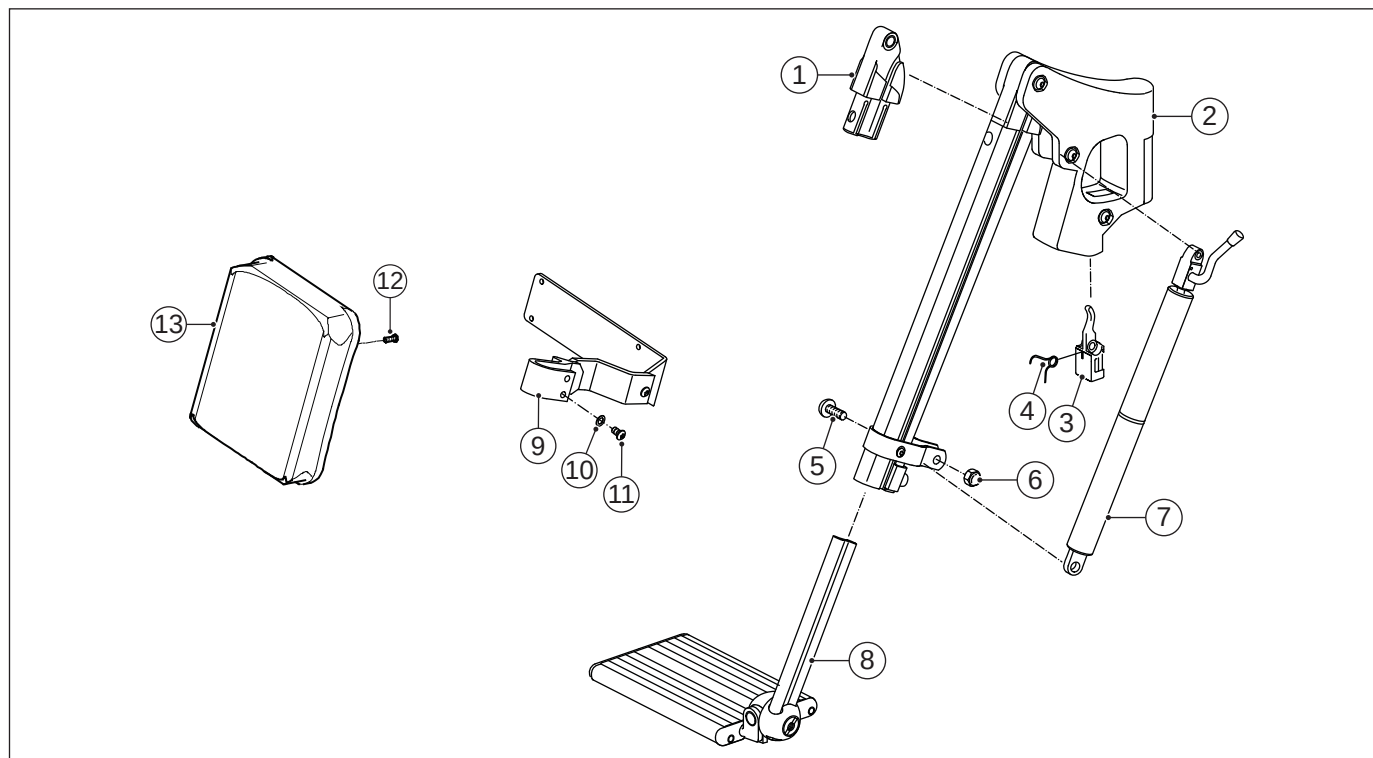


Pos	Qty	Article number	Description
1	2	1008580	Hinge Sedeo Original manually elevating legrest
1-7	1	1002516	Sedeo Original manually elevating legrest upper part right
1-7	1	1002517	Sedeo Original manually elevating legrest upper part left
1-11	1	140.00351.000	Sedeo Original manually elevating legrest for SW 42 cm (size 15) right
1-11	1	01500.9219	Sedeo Original manually elevating legrest for SW 46 cm (size 17) right
1-11	1	140.00361.000	Sedeo Original manually elevating legrest for SW 50 cm (size 19) right
1-11	1	1003601	Sedeo Original manually elevating legrest for SW 54 cm (size 21) right
1-11	1	140.00321.000	Sedeo Original manually elevating legrest for SW 42 cm (size 15) left
1-11	1	01500.9319	Sedeo Original manually elevating legrest for SW 46 cm (size 17) left
1-11	1	140.00331.000	Sedeo Original manually elevating legrest for SW 50 cm (size 19) left
1-11	1	1003598	Sedeo Original manually elevating legrest for SW 54 cm (size 21) left
3	2	01500.1112	Legrest release lever
4	2	01500.2212	Legrest release spring
5	2	00000.3706	Button head M8X25
6	2	00000.1802	Cap lock nut M8
7	1	1008077	Gaspring with gaspring control right
7	1	1008076	Gaspring with gaspring control left
-	1	1009968	Gaspring 200N
8	1	140.00171.000	Sedeo Original footplate for SW 42 cm (size 15) right
8	1	140.00121.000	Sedeo Original footplate for SW 46 cm (size 17) right
8	1	140.00181.000	Sedeo Original footplate for SW 50 cm (size 19) right
8	1	1003220	Sedeo Original footplate for SW 54 cm (size 21) right
8	1	140.00201.000	Sedeo Original footplate for SW 42 cm (size 15) left
8	1	140.00131.000	Sedeo Original footplate for SW 46 cm (size 17) left
8	1	140.00211.000	Sedeo Original footplate for SW 50 cm (size 19) left
8	1	1003219	Sedeo Original footplate for SW 54 cm (size 21) left
9-11	2	1004189	Calf plate hinge standard
9-11	2	1004188	Calf plate hinge small
10	4	00000.2601	Toothed spring washer M6
11	4	00000.3725	Button head M6X10
12	4	00000.3206	Screw M5x16
12,13	1	220.00121.505	Calf plate standard
12,13	1	220.00141.505	Calf plate small
-	1	01503.8200	Sedeo Original manually elevating legrest top complete right
-	1	01503.8300	Sedeo Original manually elevating legrest top complete left

Continued from **Manually elevating legrest**

Pos	Qty	Article number	Description
-	1	1015145	Sedeo Original manually elevating legrests with calf support for SW 42 cm, set (mounting brackets included)
-	1	1015146	Sedeo Original manually elevating legrests with calf support for SW 46 cm, set (mounting brackets included)
-	1	1015147	Sedeo Original manually elevating legrests with calf support for SW 50 cm, set (mounting brackets included)
-	1	1015148	Sedeo Original manually elevating legrests with calf support for SW 54 cm, set (mounting brackets included)
-	1	9006539	Torx set Sedeo Original legrest
-	2	1002599	Legrest inner extrusion 10 cm extended
-	2	1005119	Depth adjustable calf plate bracket (max. 4,5 cm)
-	1	1005160	Depth adjustable calf plate bracket (max. 4,5 cm), set

8 Electrically elevating legrest

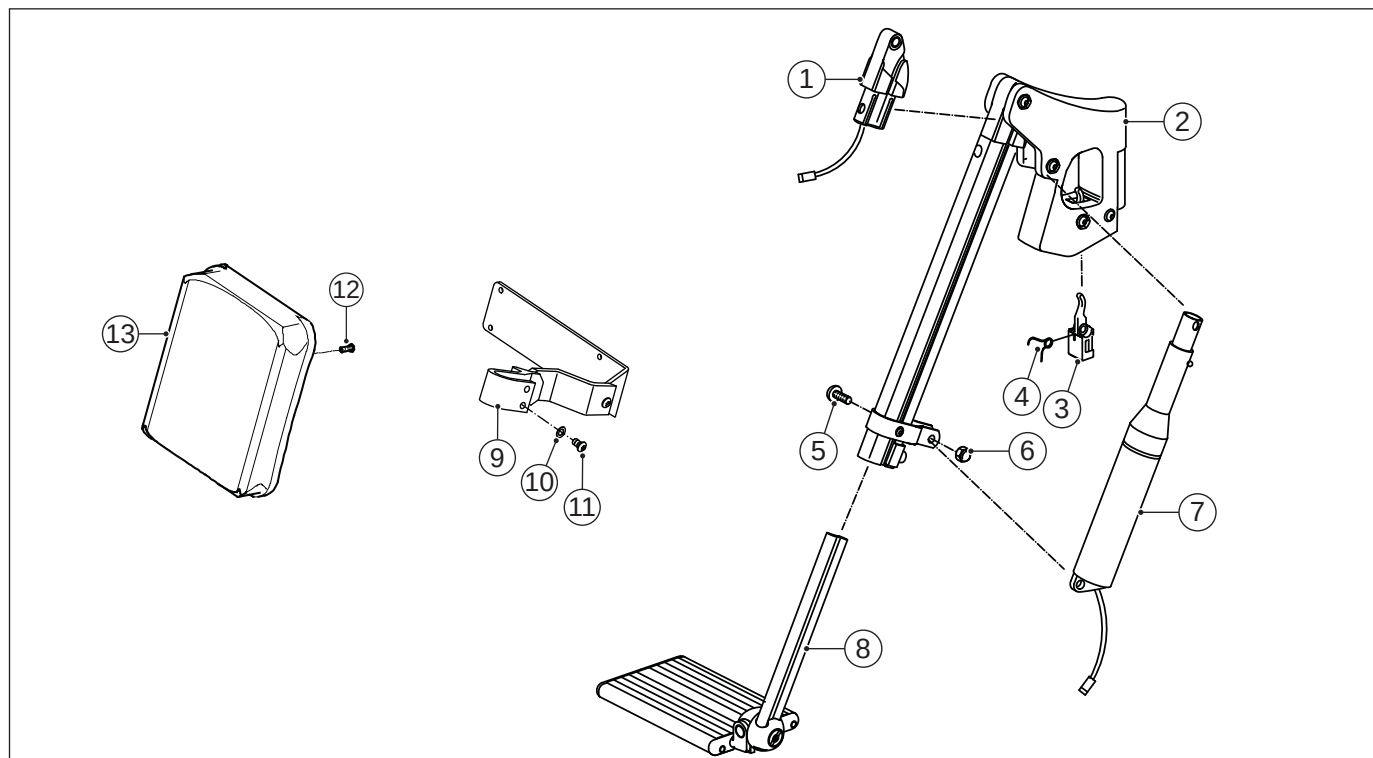


Pos	Qty	Article number	Description
1	2	1008580	Hinge Sedeo Original manually elevating legrest
1-7	1	1003028	Sedeo Original electrically elevating legrest upper part right
1-7	1	1003029	Sedeo Original electrically elevating legrest upper part left
1-11	1	1003618	Sedeo Original electrically elevating legrest for SW 42 cm (size 15) right
1-11	1	01502.9013	Sedeo Original electrically elevating legrest for SW 46 cm (size 17) right
1-11	1	1003619	Sedeo Original electrically elevating legrest for SW 50 cm (size 19) right
1-11	1	1003620	Sedeo Original electrically elevating legrest for SW 54 cm (size 21) right
1-11	1	1003613	Sedeo Original electrically elevating legrest for SW 42 cm (size 15) left
1-11	1	01502.9113	Sedeo Original electrically elevating legrest for SW 46 cm (size 17) left
1-11	1	1003614	Sedeo Original electrically elevating legrest for SW 50 cm (size 19) left
1-11	1	1003615	Sedeo Original electrically elevating legrest for SW 54 cm (size 21) left
3	2	01500.1112	Legrest release lever
4	2	01500.2212	Legrest release spring
5	2	00000.3706	Button head M8X25
6	2	00000.1802	Cap lock nut M8
7	2	9001232	Legrest actuator inc. kabel
8	1	140.00171.000	Sedeo Original footplate for SW 42 cm (size 15) right
8	1	140.00121.000	Sedeo Original footplate for SW 46 cm (size 17) right
8	1	140.00181.000	Sedeo Original footplate for SW 50 cm (size 19) right
8	1	1003220	Sedeo Original footplate for SW 54 cm (size 21) right
8	1	140.00201.000	Sedeo Original footplate for SW 42 cm (size 15) left
8	1	140.00131.000	Sedeo Original footplate for SW 46 cm (size 17) left
8	1	140.00211.000	Sedeo Original footplate for SW 50 cm (size 19) left
8	1	1003219	Sedeo Original footplate for SW 54 cm (size 21) left
9-11	2	1004189	Calf plate hinge standard
9-11	2	1004188	Calf plate hinge small
10	2	00000.2601	Toothed spring washer M6
11	2	00000.3725	Button head M6X10
12	4	00000.3206	Screw M5x16
12,13	1	220.00121.505	Calf plate standard
12,13	1	220.00141.505	Calf plate small
-	1	01503.8400	Sedeo Original electrically elevating legrest top complete right
-	1	01503.8500	Sedeo Original electrically elevating legrest top complete left
-	1	1015149	Sedeo Original electrically elevating legrests with calf support for SW 42 cm, set (mounting brackets included)

Continued from **Electrically elevating legrest**


Pos	Qty	Article number	Description
-	1	1015150	Sedeo Original electrically elevating legrests with calf support for SW 46 cm,set (mounting brackets included)
-	1	1015151	Sedeo Original electrically elevating legrests with calf support for SW 50 cm,set (mounting brackets included)
-	1	1015152	Sedeo Original electrically elevating legrests with calf support for SW 54 cm,set (mounting brackets included)
-	1	9006539	Torx set Sedeo Original legrest
-	2	1002599	Legrest inner extrusion 10 cm extended
-	2	1005119	Depth adjustable calf plate bracket (max. 4,5 cm)
-	1	1005160	Depth adjustable calf plate bracket (max. 4,5 cm), set
-	2	9001233	Legrest actuator harness

9 Central legrest

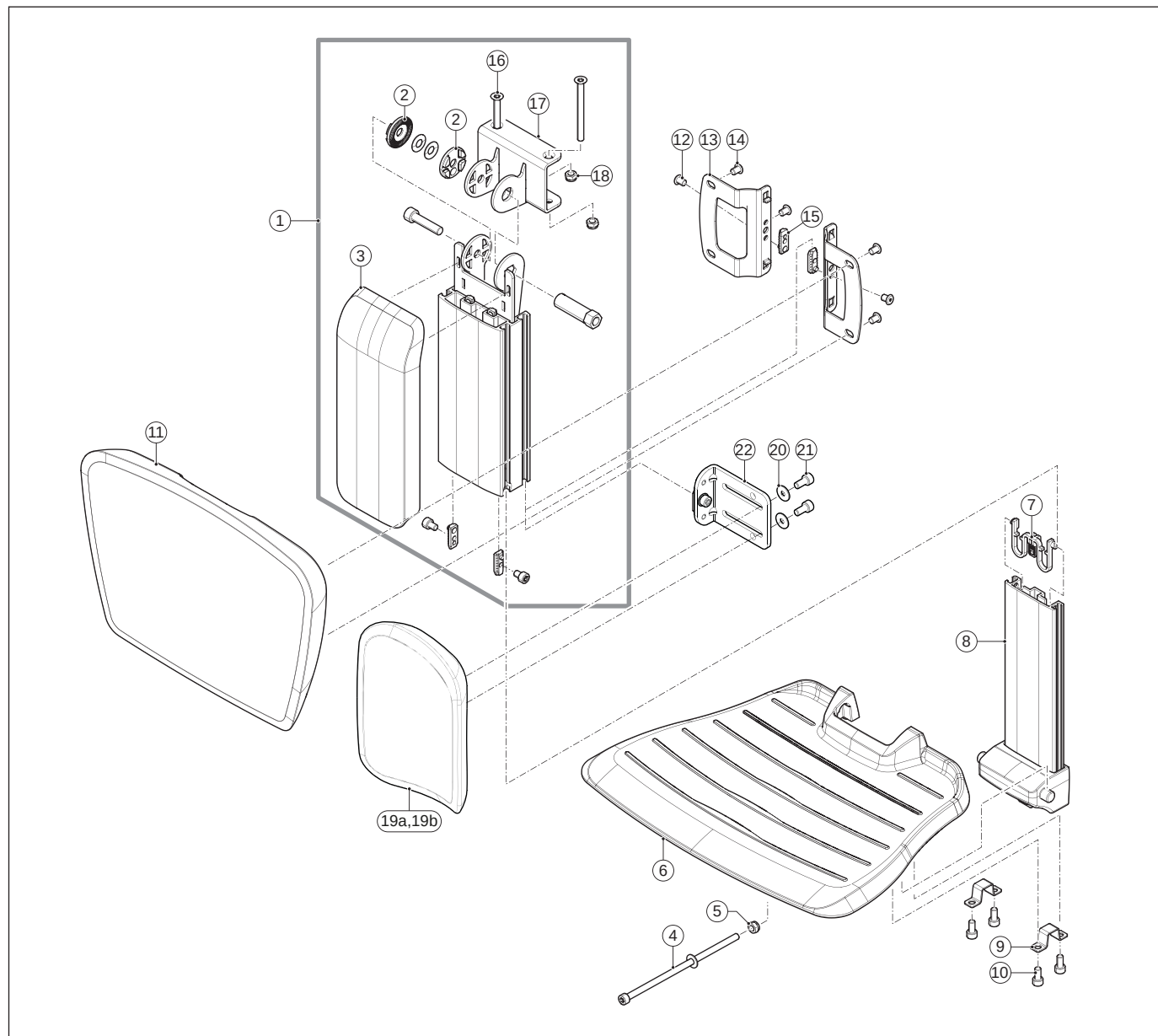

Pos	Qty	Article number	Description
1	1	1010941	Standard central legrest upper part
1-10,16-18	1	1011156	Standard central legrest
1-18	1	1011157	Standard central legrest with calf plate
2	1	1010943	Adjustment wheel, set
3	1	1010942	Central legrest cover
4	1	9003894	Angle adjustment bolt
5	1	-	Lock nut M6 D985/8/VZ
4-10	1	1010939	Central legrest lower part
4,5,6,9,10	1	1010940	Footplate
10	4	1015332	Torx socket head screw M6x12 iso14579 Precote 85 (100 pieces)
11-15	1	1010948	Central legrest calf plate with mouning set
11,14	1	1010947	Central legrest calf plate
12,13,15	1	1010946	Central legrest calf plate mounting, set
12	2	1015330	Torx socket head screw M6x8 din 7984T Precote 85 (100 pieces)
14	4	1015330	Torx socket head screw M6x8 din 7984T Precote 85 (100 pieces)
16-18	1	1016446	Standard central legrest mounting bracket
19a,20,21	2	1016775	Separate calf plate for central legrest (L=R)
19b,20,21	2	1016776	Separate calf plate side support for central legrest (L=R)
19a,20-22	2	9010234	Separate calf plate for central legrest with mounting set (L=R)

10 Electrically elevating central legrest

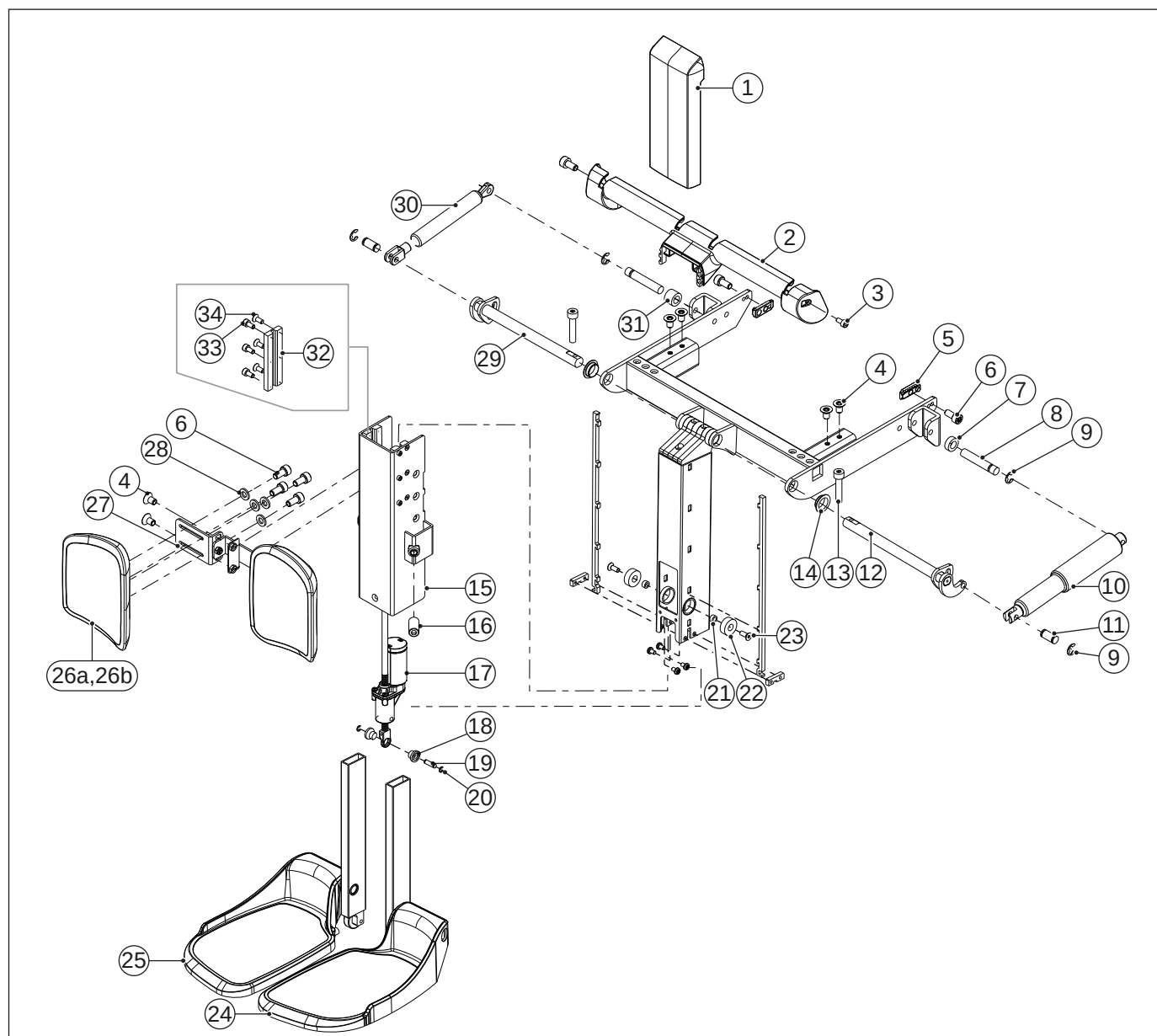


Pos	Qty	Article number	Description
1	1	1016763	Front cover
2,3	1	1016764	Axle cover
4	8	1016778	Torx countersunk socket head screw M6x10 ISO10642 (100 pieces)
4,6,26a,27,28	1	1016774	Calf plate bracket complete with calf plates
4,6,26b,27,28	1	1016782	Calf plate bracket complete with calf plates side support
4,27	1	1016777	Calf plate mounting bracket
5	2	1010881	Clamp nut profiles (10 pieces)
6	10	1015332	Torx socket head screw M6x12 ISO14579 (100 pieces)
6,26a,28	2	1016775	Separate calf plate for central legrest (L=R)
6,26b,28	2	1016776	Separate calf plate side support for central legrest (L=R)
7-11	1	1016768	Tilt actuator
8,9,11,30,31	1	1016767	Gas spring
12-14	1	1016766	Tilt actuator axle
13	2	1016780	Torx socket head screw M6x30 ISO14579 (100 pieces)
13,14,29	1	1016765	Gas spring axle
15,16,32-34	1	1016770	Outside profile
16	2	1016779	Torx set screw M10x 10 ISO4026 (100 pieces)
17-23	1	1016769	Height adjustment actuator
24	1	1016773	Footplate left
25	1	1016772	Footplate right

Continued from **Electrically elevating central legrest**


Pos	Qty	Article number	Description
32-34	1	1016771	Slider bearing
1-25,29-34	1	1016781	Electrical elevating central legrest without calfplates
1-25,29-34	1	1017106	Electrical elevating central legrest with lenght adjustment without calfplates

11 Flip up armrest

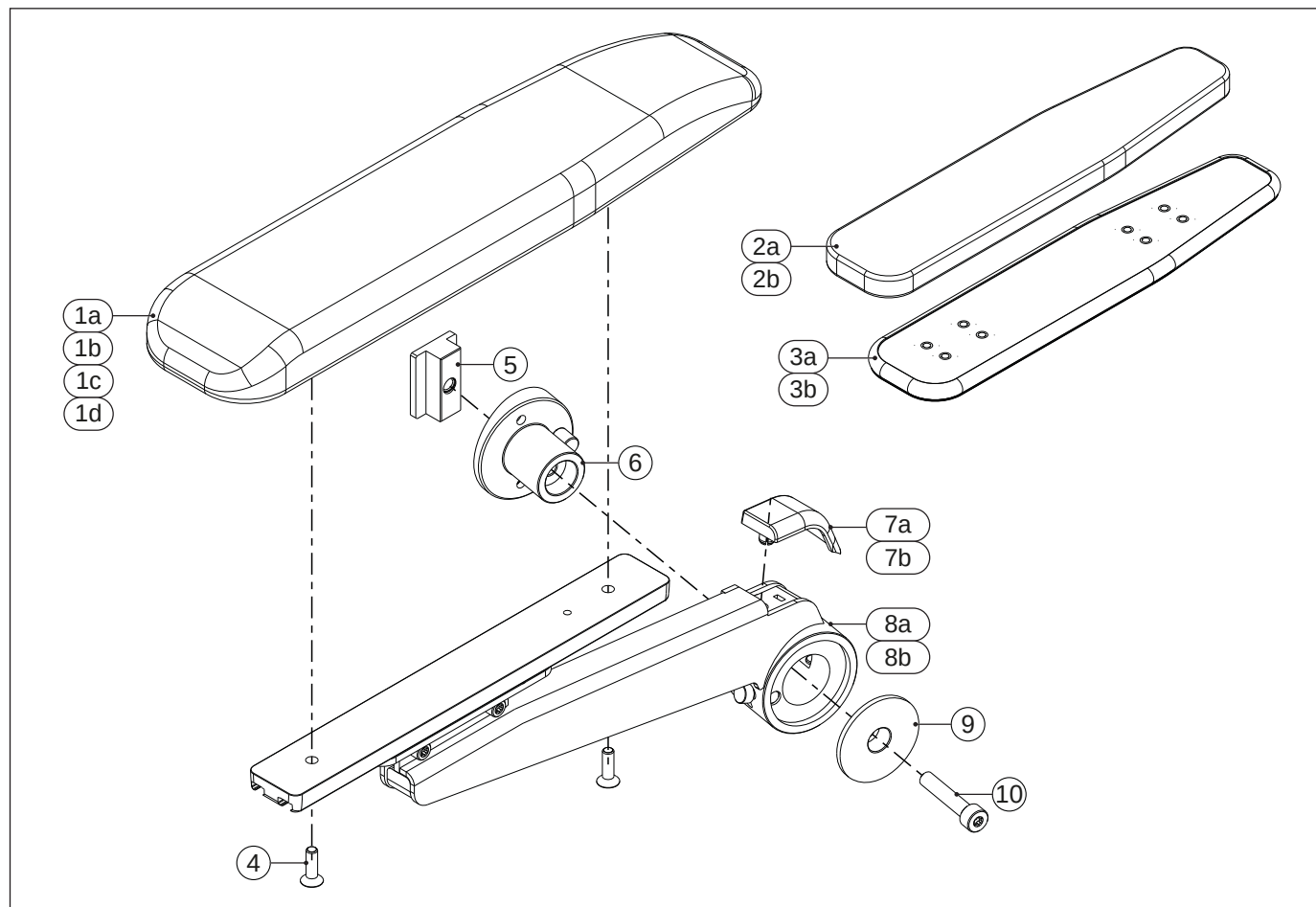


Pos	Qty	Article number	Description
1a,4	2	1017052	Armrest pad 30 cm upholstered, stretch rubber, right = left
1b,4	2	1017053	Armrest pad 40 cm upholstered, stretch rubber, right = left
1c,4	1	1010908	Armrest pad 30 cm, PU, right
1d,4	1	1010909	Armrest pad 30 cm, PU, left
2a	2	1017056	Armrest pad top 30 cm removeable upholstery, fabric, right = left
2b	2	1017057	Armrest pad top 40 cm removeable upholstery, fabric, right = left
2a,3a,4	2	1017054	Armrest pad 30 cm removeable upholstery, fabric, right = left
2b,3b,4	2	1017055	Armrest pad 40 cm removeable upholstery, fabric, right = left
4	4	1015340	Torx countersunk socket head screw M6x20 iso10642 Precote 85 (100 pieces)
5,6,7a,8a,9,10	1	1016791	Flip up armrest, right
5,6,7b,8b,9,10	1	1016792	Flip up armrest, left
7a,7b	1	1016793	Armrest cover, set

12 Headrest

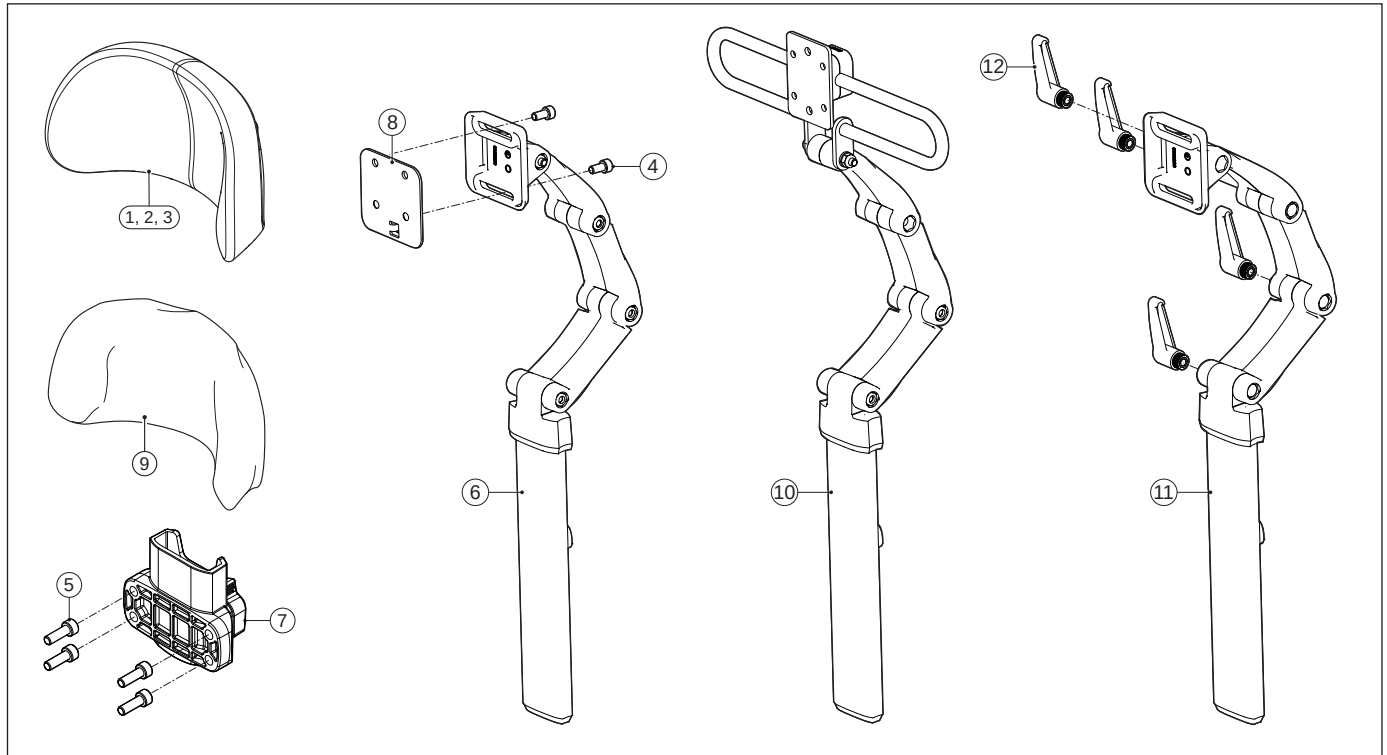


Pos	Qty	Article number	Description
1,4	1	1010957	Headrest small
2,4	1	1010996	Headrest large
3,4	1	1010997	Headrest flexible
4	2	1015332	Torx socket head screw M6x12 iso14579 Precote 85 (100 pieces)
5	4	1015333	Torx socket head screw M6x16 iso14579 Precote 85 (100 pieces)
5,7	1	1017703	Headrest socket set
6	1	1010954	Headrest metal parts
8	1	9004116	Headrest adaptor for Sedeo Original
9	1	3009220	Headrest cover small
10	1	1010994	Headrest metal parts width adjustable
11,12	1	1010955	Headrest parts quick release
12	1	1017947	Quick release lever for headrest metal parts, set
-	1	1015403	Headrest hinge set

13 Joystick bracket

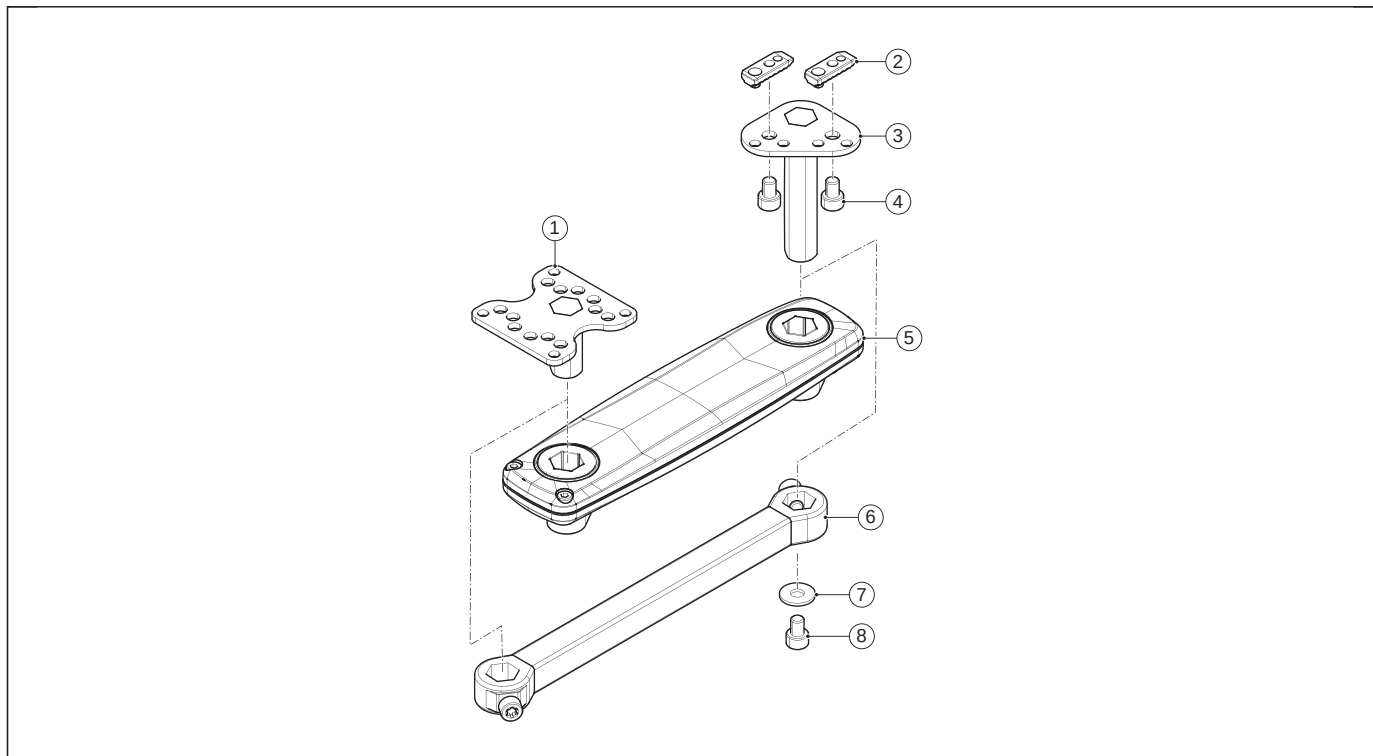


Pos	Qty	Article number	Description
1	1	1010949	Remote interface bracket
2	4	1010881	Clamp nut (10 pieces)
2-4, 7, 8	1	1010950	Remote Armrest bracket
4	2	1015330	Torx socket head screw M6x8 din7984T Precote 85 (100 pieces)
5	1	1010951	Swing-away remote holder
6	1	1010952	Standard remote holder
1-5, 7, 8	1	1011164	Swing-away remote holder complete
1-4, 6-8	1	1011165	Standard remote holder complete
7	1	-	Washer 3xD M6 D9021
8	1	1015331	Torx socket head screw M6x10 iso14579 Precote 85 (100 pieces)

14 Lighting

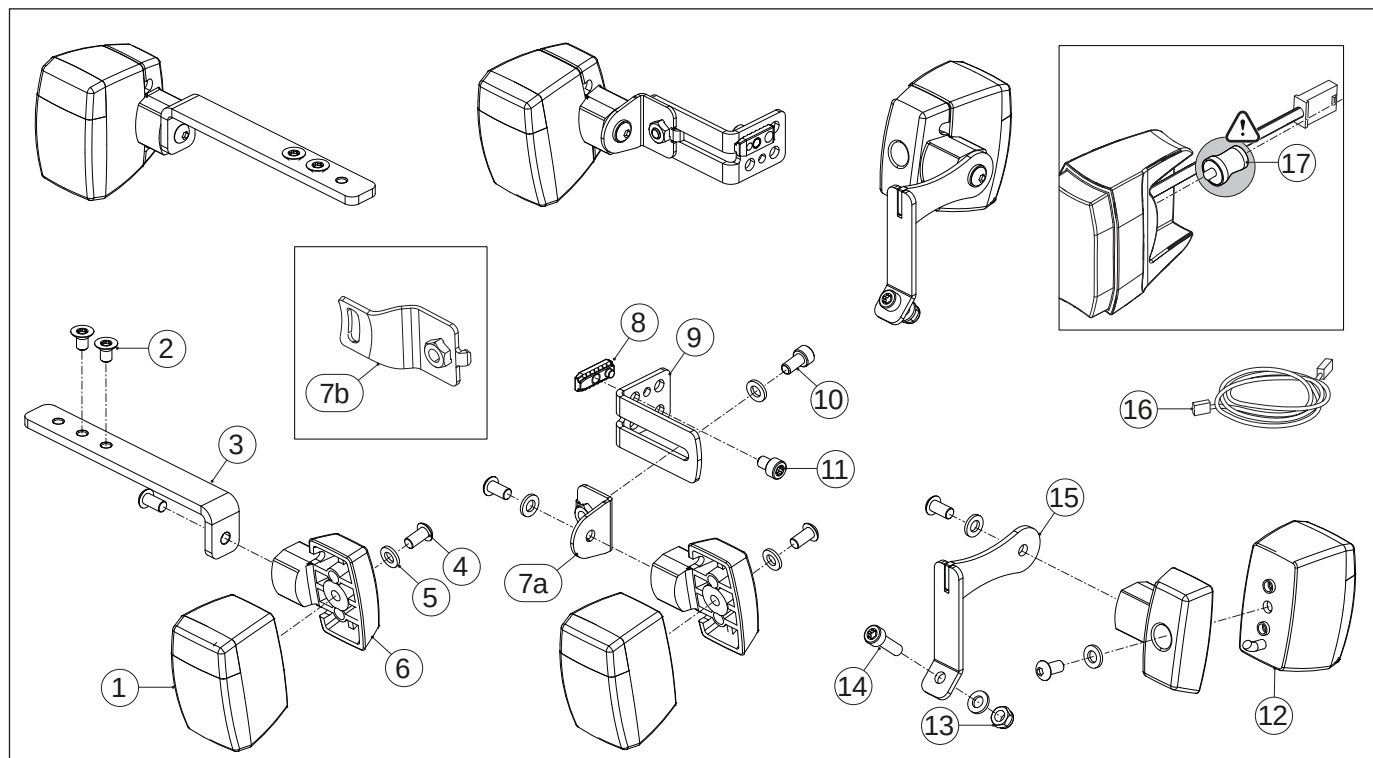


Pos	Qty	Article number	Description
1	2	6002661	Front light module standard
1	2	6000038	Front light module StVZO
1,4-6,7a,8-11	2	1016799	Front light module for mounting on seat frame, left = right
1,4-6,7a,8-16	1	1016796	Lighting set R-net for mounting on seat frame
1-6	2	1016797	Front light module for mounting with central electrical legrest, left = right
1-6,12-16	1	1016795	Lighting set R-net for mounting with central electrical legrest
2	4	1016778	Torx countersunk socket head screw M6x10 iso10642 (100 pieces)
4	8	-	Buttonhead M6x12 ISO7380
4-6,12-15	2	1016801	Rear light module incl mounting bracket, left = right
5,7a,8-11	2	1018044	Lighting bracket front on seat frame, left = right, version > June 2015
5,7b,8-11	2	1016804	Lighting bracket front on seat frame, left = right, version < July 2015
6	4	6000056	Light rubber mounting unit, left = right
10	2	1015331	Torx socket head screw M6x10 iso14579 Precote 85 (100 pieces)
12	2	6000039	Rear light module
14	2	1016898	Torx socket head screw M6x20 iso14579 (100 pieces)
16	1	9006331	Lighting harness R-net lighting
16	1	9006332	Lighting harness R-net lighting and horn
17	4	1016517	Rubber damper (10 pieces) for version < July 2015

15 Abduction support

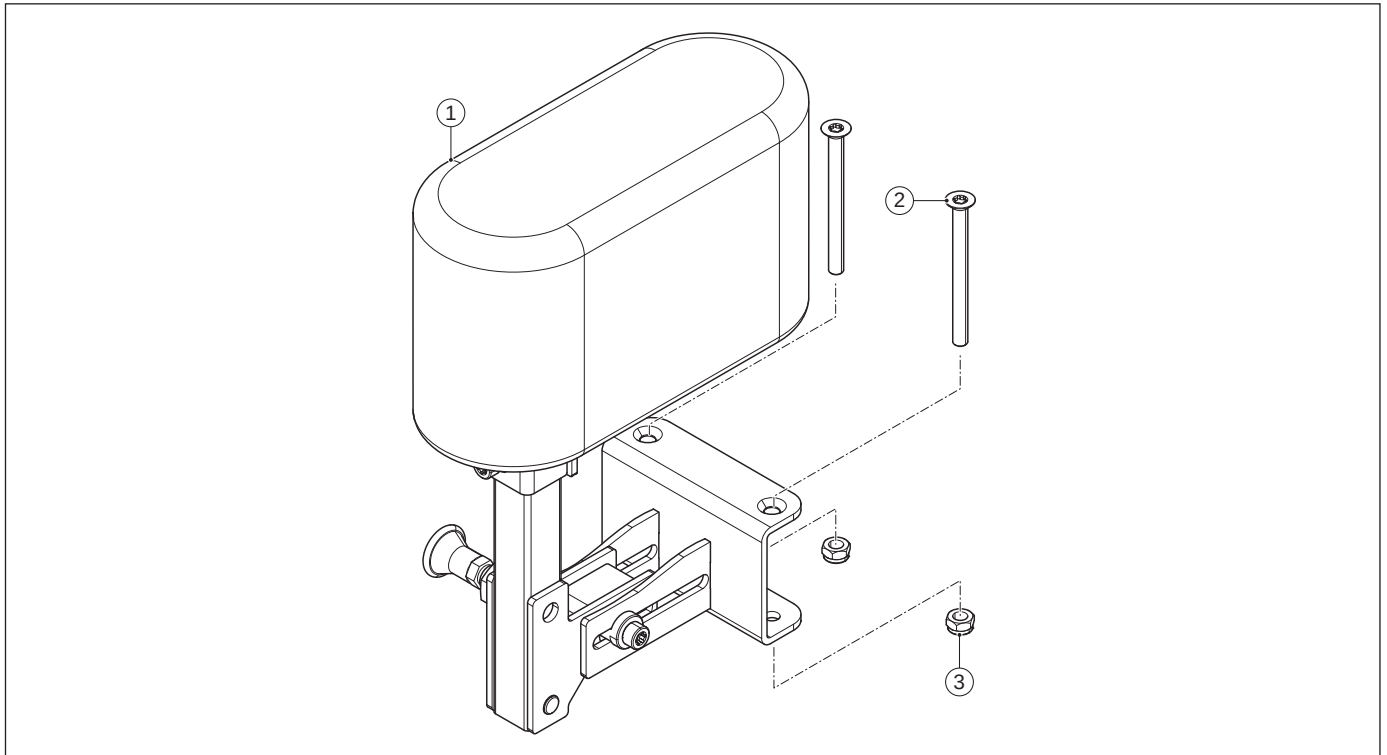


Pos	Qty	Article number	Description
1-3	1	1011002	Abduction support
2	2	1015342	Torx countersunk socket head screw M6x70 iso10642 (100 pieces)
3	2	-	Lock nut M6 D985/8/VZ

16 Mirror - Side plate

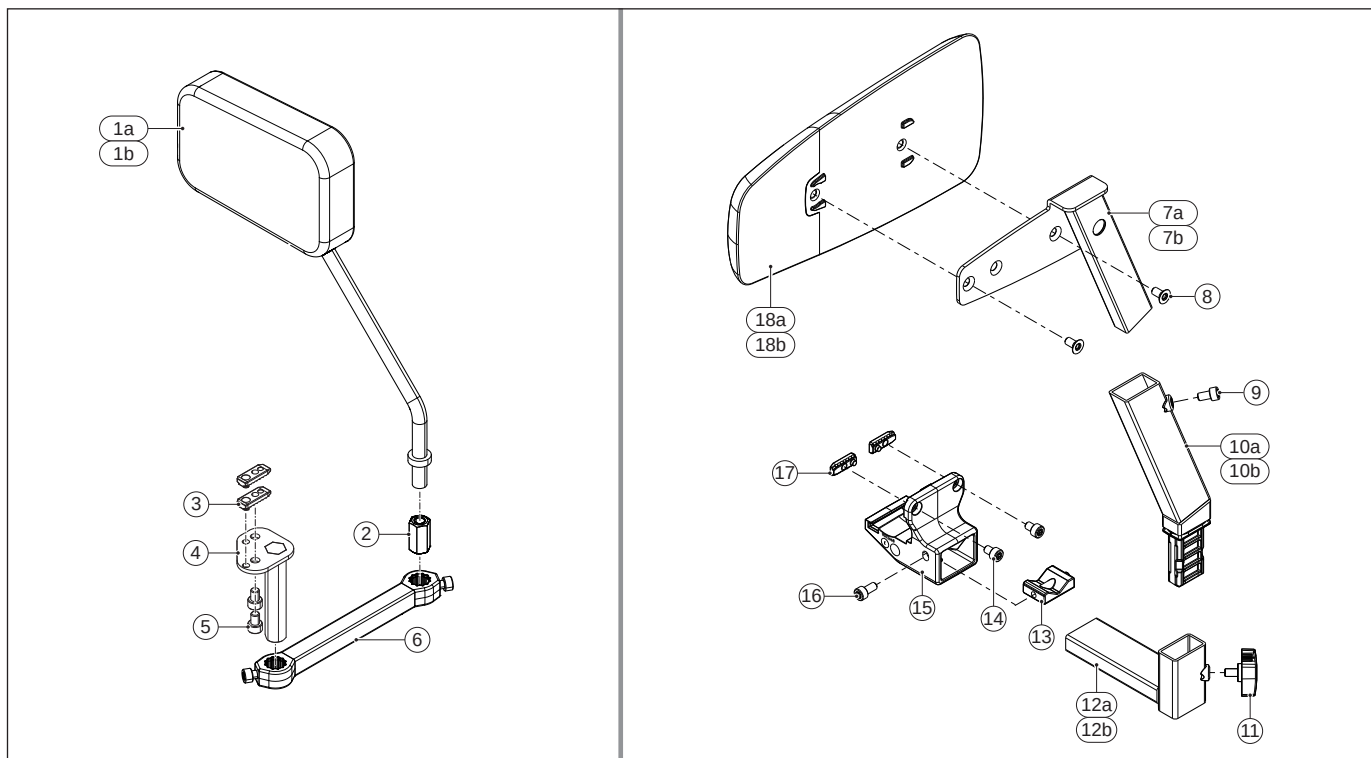


Pos	Qty	Article number	Description
1a,2-6	1	1012319	Mirror not joystick side, right
1b,2-6	1	1012321	Mirror not joystick side, left
1a,2-6	1	1012269	Mirror joystick side, right
1b,2-6	1	1012320	Mirror joystick side, left
3	4	1010881	Clamp nut (10 pieces)
7a,9,10a	1	1016902	Side support bracket, right
7b,9,10b	1	1016903	Side support bracket, left
8	4	1015339	Torx countersunk socket head screw M6x12 iso10642 (100 pieces)
8,18a	2	1016900	Side support PU, right = left
8,18b	2	1016901	Side support PA, right = left
9	2	1015331	Torx socket head screw M6x10 iso14579 Precote 85 (100 pieces)
11,12a	1	1012936	Armrest / side support socket, right
11,12b	1	1012937	Armrest / side support socket, left
13-17	1	1010880	Armrest / side support holder, set
14	4	1015330	Torx socket head screw M6x8 din7984T Precote 85 (100 pieces)
16	2	1015332	Torx socket head screw M6x12 iso14579 Precote 85 (100 pieces)
17	4	1010881	Clamp nut (10 pieces)
7a,9,10a,11,12a,13-17	1	1017219	Side support complete, right
7b,9,10b,11,12b,13-17	1	1017220	Side support complete, left

17 Armrest bag - Bag hooks - Hemiplegia armrest

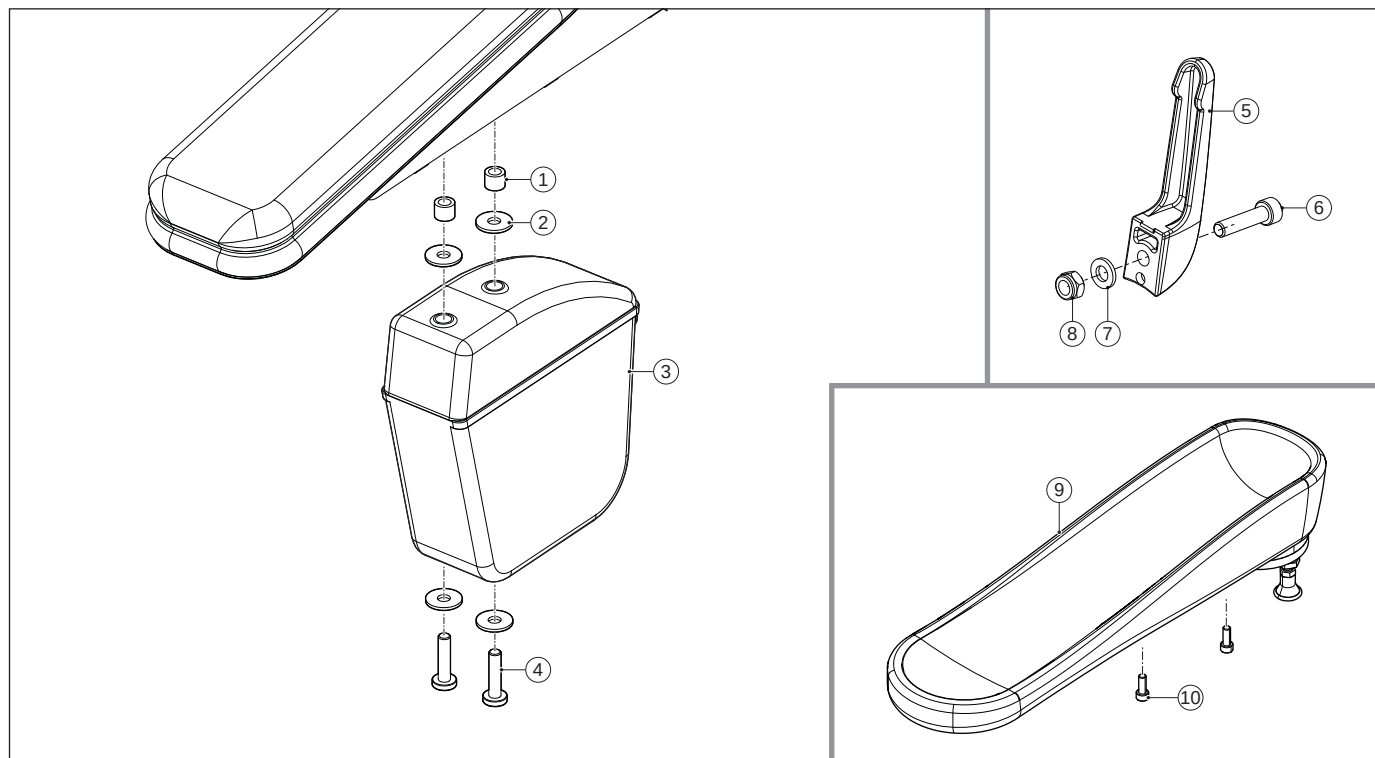


Pos	Qty	Article number	Description
1-4	1	1017051	Armrest bag including mounting bracket
4	2	-	Screw Pan Head Torx M6x25 Taptite
5-8	1	1016794	Bag hooks, set
6	1	1015334	Torx socket head screw M6x25 iso14579 (100 pieces)
7	1	-	Washer flat M6 DIN125A/VZ
8	1	-	Lock nut M6 DIN985
9-10	1	1010977	Hemiplegia armrest
10	2	1015333	Torx socket head screw M6x16 iso14579 Precote 85 (100 pieces)

18 Positioning belt

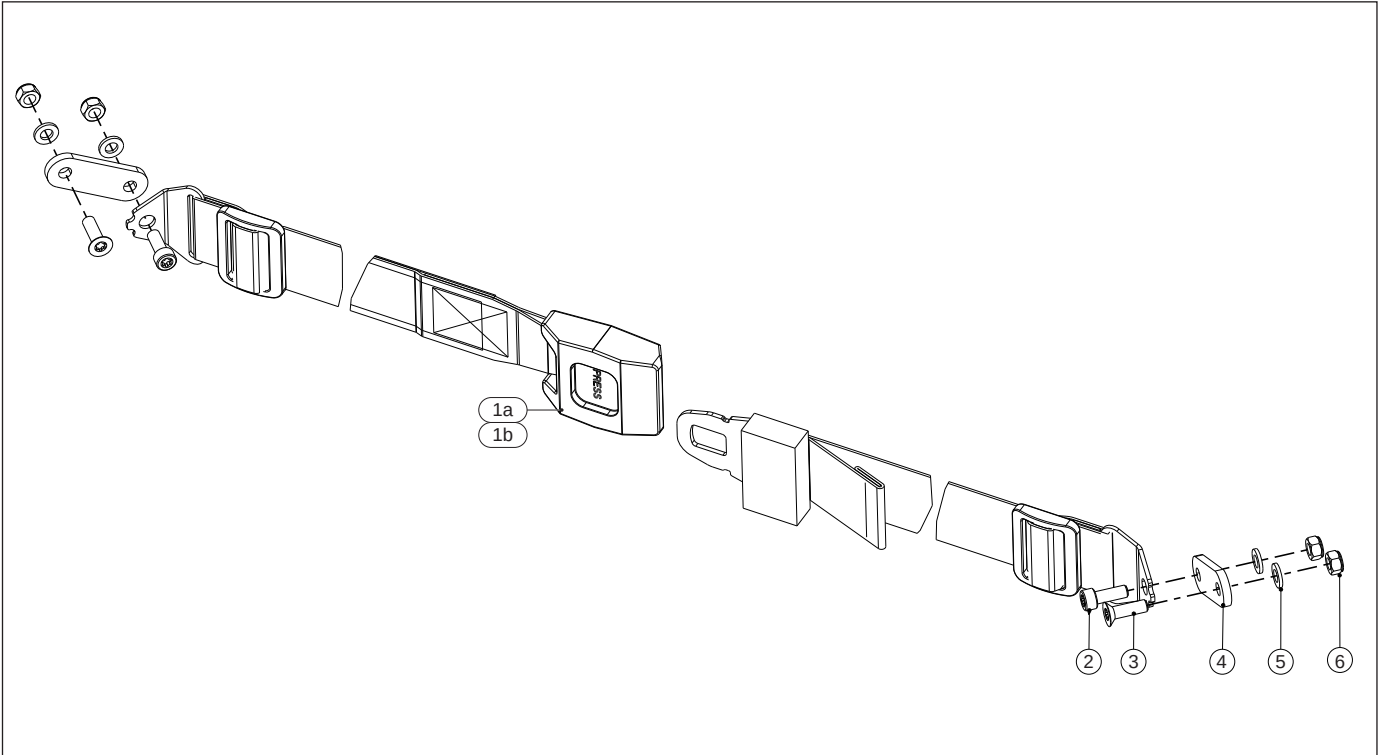


Pos	Qty	Article number	Description
1a,2-6	1	1016819	Positioning belt short
1b,2-6	1	1016820	Positioning belt long
2	2	1015333	Torx socket head screw M6x16 iso14579 Precote 85 (100 pieces)
3	2	1015340	Torx countersunk socket head screw M6x20 iso10642 Precote 85 (100 pieces)
5	4	-	Flat washer M6 DIN125
6	4	-	Lock nut M6 DIN985

19 Cane holder - Thoracic support

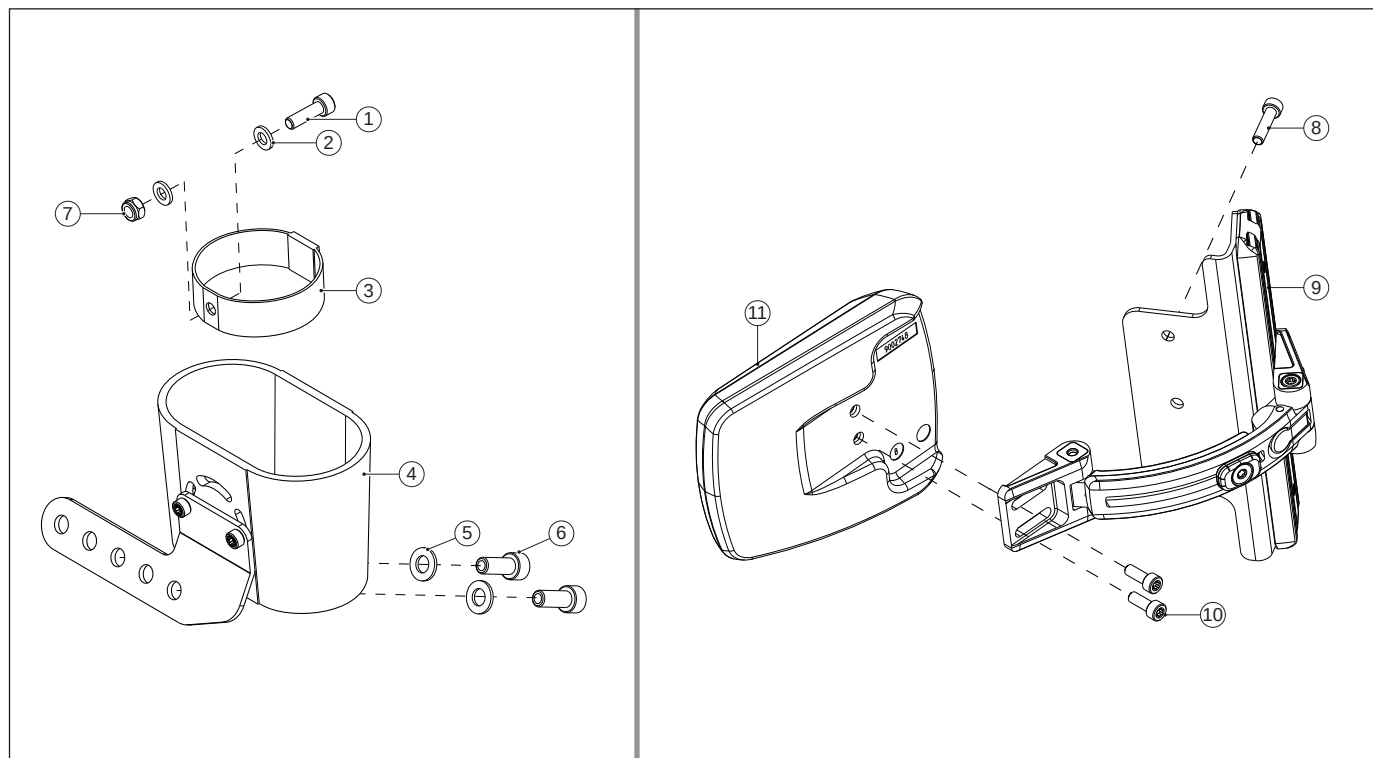

Pos	Qty	Article number	Description
1	1	1016898	Torx socket head screw M6x20 iso14579 (100 pieces)
1-7	1	1016897	Cane holder
6	1	1016850	Torx socket head screw M8x16 iso14579 (100 pieces)
8	1	1015344	Torx socket head screw m5x25 iso14579 (100 pieces)
8,9	1	1016905	Thoracic support mounting, right
8,9	1	1016906	Thoracic support mounting, left
8,9	1	1016907	Thoracic support mounting, set
8-11	1	1016904	Thoracic support complete, set
10	1	1015332	Torx socket head screw M6x12 iso14579 Precote 85 (100 pieces)
10,11	2	1010980	Thoracic support cushion, right = left

20 Worktops power wheelchair

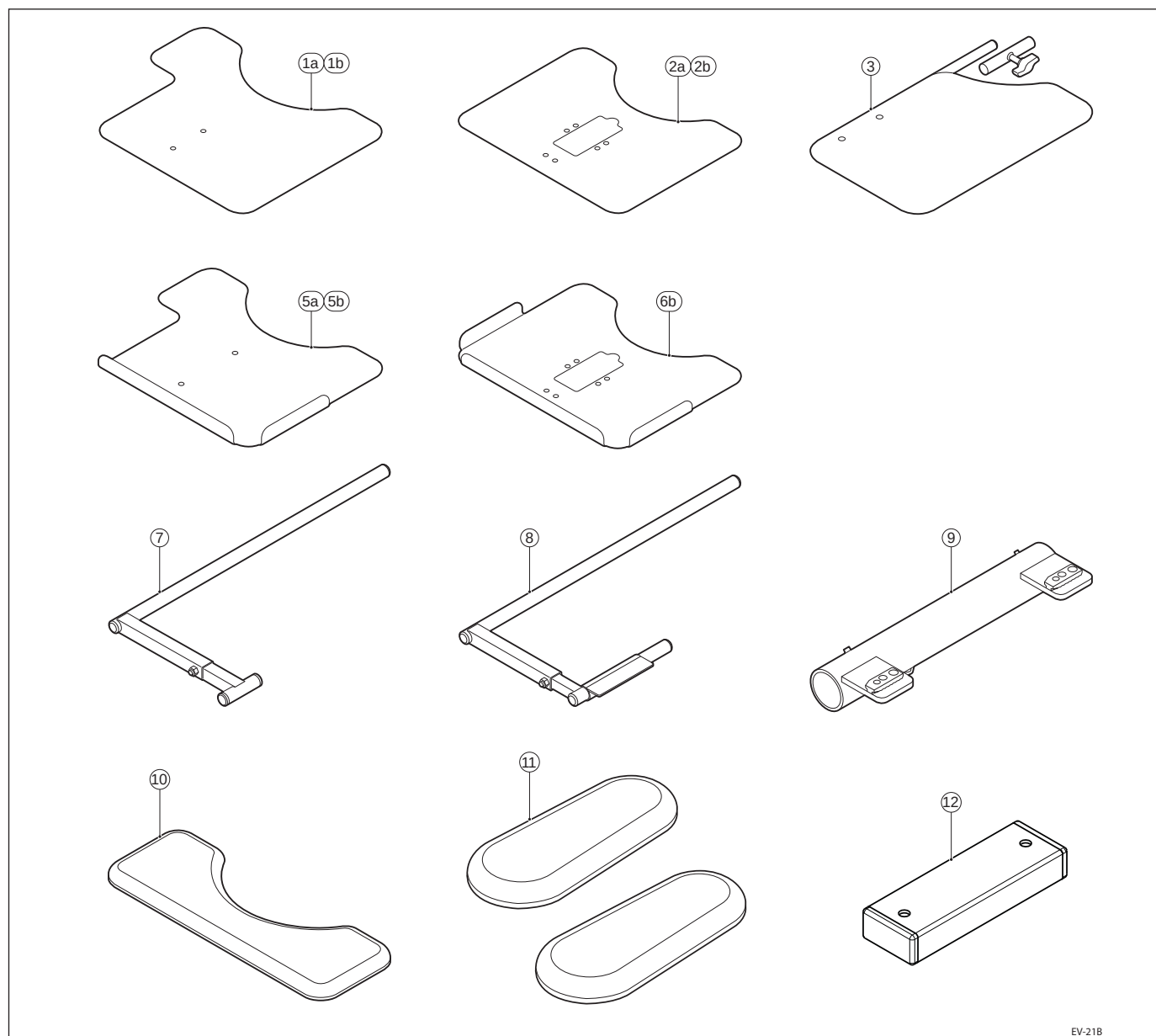

EV-21B

Pos	Qty	Article number	Description
1a,8,9	1	1015167	Worktop, folding, acrylic + joystick recess right hand side
1b,8,9	1	1015168	Worktop, folding, acrylic + joystick recess left hand side
2a,7,9	1	1015169	Worktop with center remote recess for DX
2b,7,9	1	1015455	Worktop with center remote recess for R-net
3	1	1015170	Half worktop, acrylic
5a,8,9	1	1015162	Worktop, folding, acrylic with raised edges, suitable for upholstery + joystick recess right hand side
5b,8,9	1	1015164	Worktop, folding, acrylic with raised edges, suitable for upholstery + joystick recess left hand side
6b,7,9	1	1015165	Worktop, folding, acrylic with raised edges, suitable for upholstery, with Center Remote Recess
9	1	1015157	Worktop receiving tube assy
10	1	223.00041.505	Partially padded, for worktop, black Sellaskin
11	2	221.00211.505	Elbow pad
12	1	1015664	Worktop holder spacer
-	1	1015171	Receiving tube for hemi worktop

5 Service instructions

5.1 Maintenance plan / Reimbursement plan

Below, we have indicated what needs to be checked, how often this should be done and by whom.

Time	Description	To be carried out by	
		User	Supplier
Monthly	• Cleaning the seating system	x	-
	• Cleaning the upholstery (if necessary)	x	-
Annually	• Inspecting the mechanical parts	-	x
	• Checking all fastenings and bolts: tighten if necessary	-	x

For disinfecting spray the chair with alcohol.

This document is a reference book to be used to order parts for the Sedeo seating system.

Assembly, replacement and adjustment instructions

The next paragraph cover illustrated assembly, replacement and adjustment instructions.

Every separate visual instruction is prefaced with:

- Preparation remarks
- Instructions
- Relevant article numbers
- Tools used
- Notice (if applicable): extra notes covering information you need when executing the specific task)

5.2 Assembly, replacement and adjustment instructions Sedeo Pro+

5.2.1 Adjusting or mounting the seat cushion fixation set

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

Follow the illustrated instructions for adjusting or mounting the seat cushion fixation.



Note! Ex-factory the fixation brackets on the cushion or seat are always mounted in the middle position! In case you want to reposition the seat cushion, place the fixation brackets at the front and backside always in the same holes! The seat board Sedeo SD40 has only two positions. Ex-factory the fixation brackets on the seat board are mounted in the back position.

Relevant article numbers

1010953 Seat cushion fixation set

1012267 Interface board Sedeo Original

Tools used

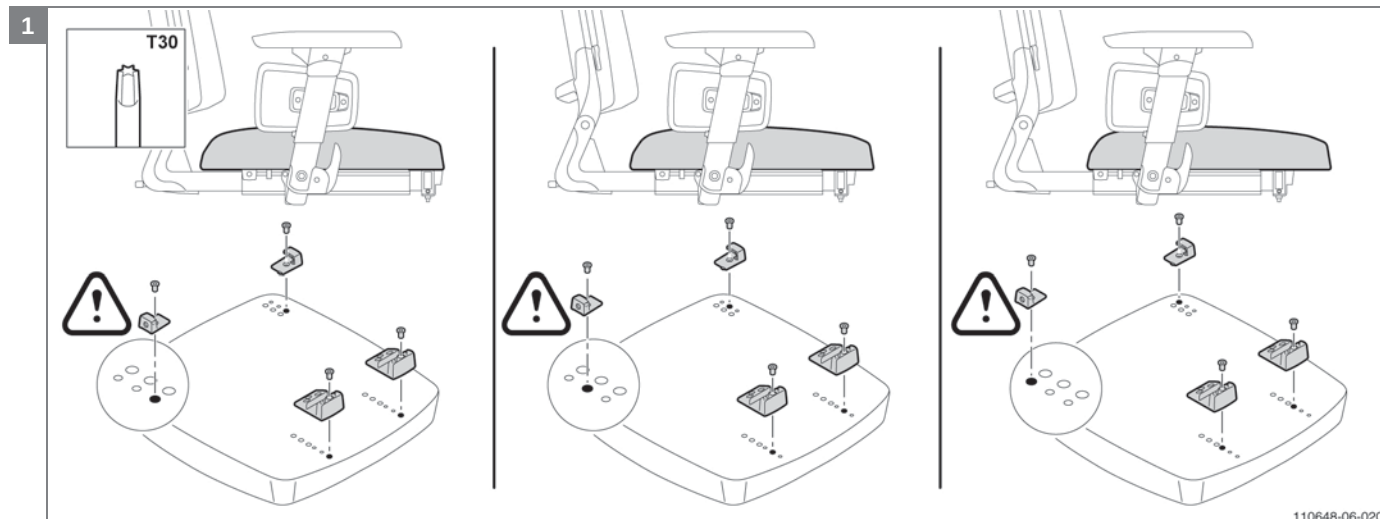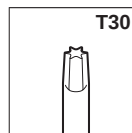

5.2.2 Replacing the seat cushion cover

Preparation

Note! Switch off the wheelchair via remote control.

Instructions






Note! Suitable work area required!


- Pull the seat up and forwards to remove it from the chair.
- Unscrew the bolts (4x) from the mounting brackets at the bottom of the seat.
- Remove the mounting brackets (4x) from bottom plate.
- Unzip the seat cushion cover.
- Replace the seat cushion cover and mount the seat in reversed order.

Relevant article numbers

All seat cushions, see service manual page 14

Tools used

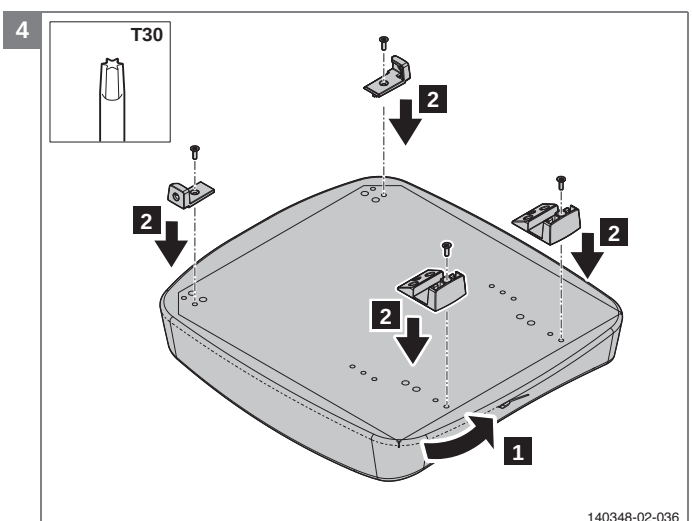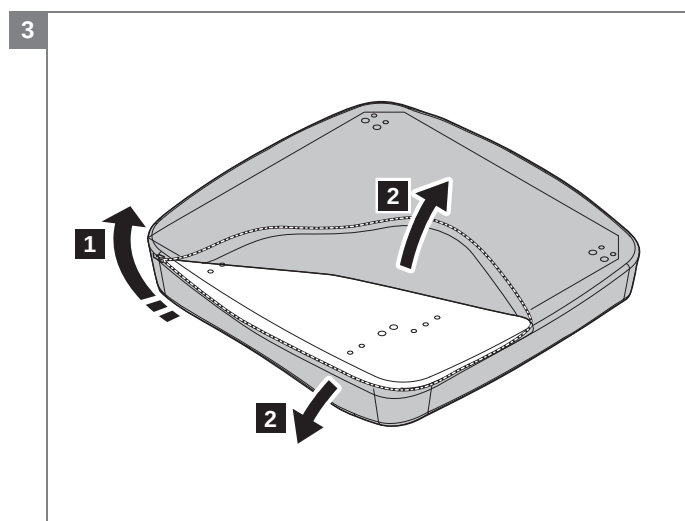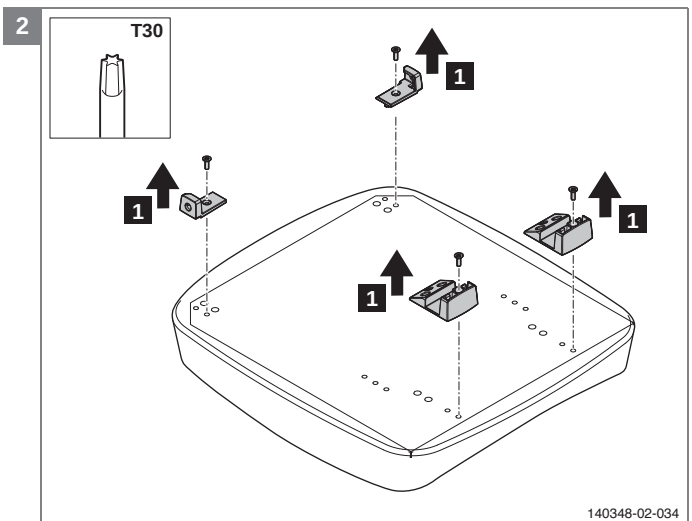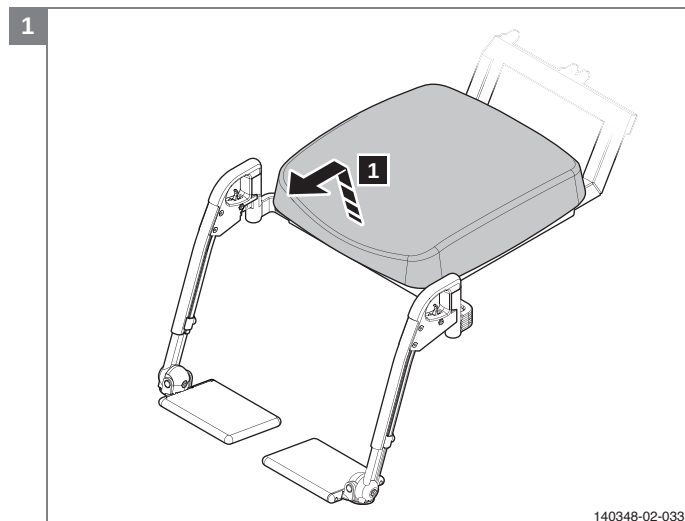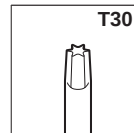

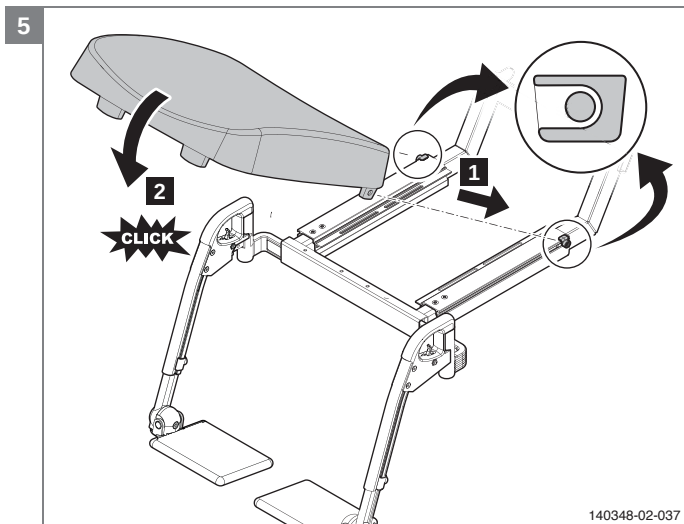

5.2.3 Mounting the backrest cushion fixation set

Preparation

Note! Switch off the wheelchair via remote control.



Instructions

Note! Suitable work area required!

- **Note!** Define by feeling the lowest hole at the top of the frame and the lowest hole at the bottom of the frame.
- Rebuild afterwards the wheelchair in reversed order.

Relevant article numbers

1016821 Backrest mounting set

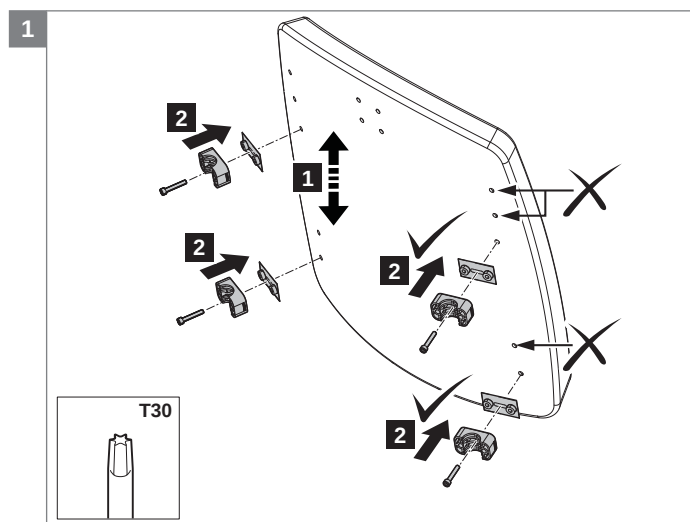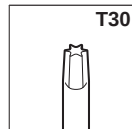Tools used

5.2.4 Replacing the backrest cushion cover

Preparation

Note! Switch off the wheelchair via remote control.

Instructions





Note! Suitable work area required!

- Unscrew the bolts (8x) from the mounting brackets, which attach the backrest to the frame.
- Remove the mounting brackets (4x).
- Pull the backseat upwards to remove it from the chair.
- Unscrew the bolts (4x) from the counterparts of the mounting brackets.
- Remove the counterparts of the mounting brackets from the backrest.
- Unzip the backrest cushion cover.
- Replace the backseat cushion cover and mount the backrest in reversed order.

Relevant article numbers

All backrest cushions, see service manual page 23

Tools used

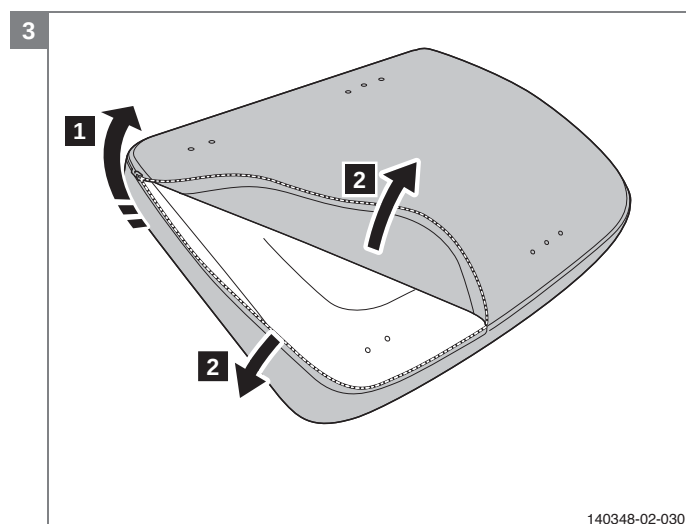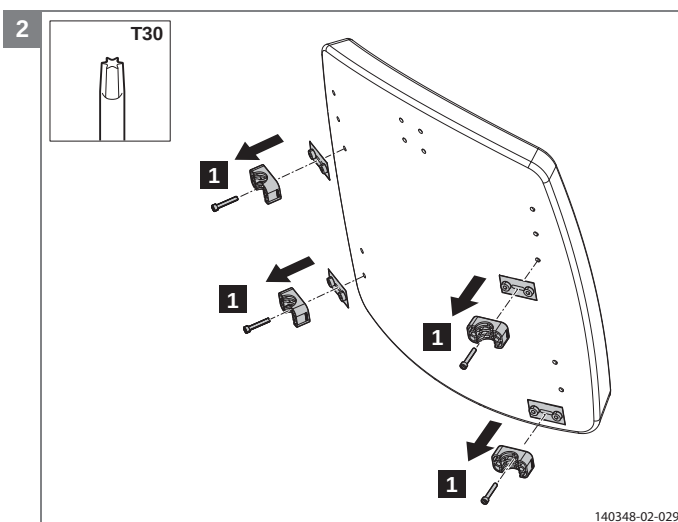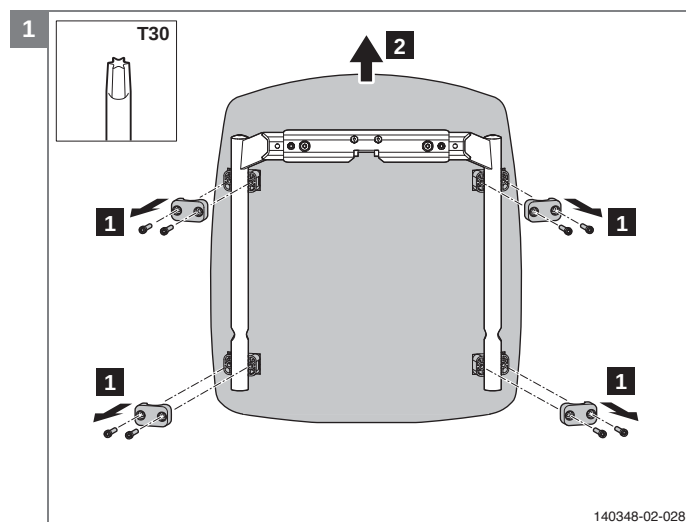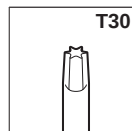

5.2.5 Mounting the thoracic supports

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

- Follow the illustrated instructions for mounting the thoracic supports.
- Rebuild afterwards the wheelchair in reversed order.



Notice

Some parts needs to be replaced; the replacement is indicated by a  symbol.

Relevant article numbers

1016905 Thoracic support mounting, right
 1016906 Thoracic support mounting, left
 1016907 Thoracic support mounting, set
 1010980 Thoracic support cushion, right = left
 1016904 Thoracic support complete, set

Tools used

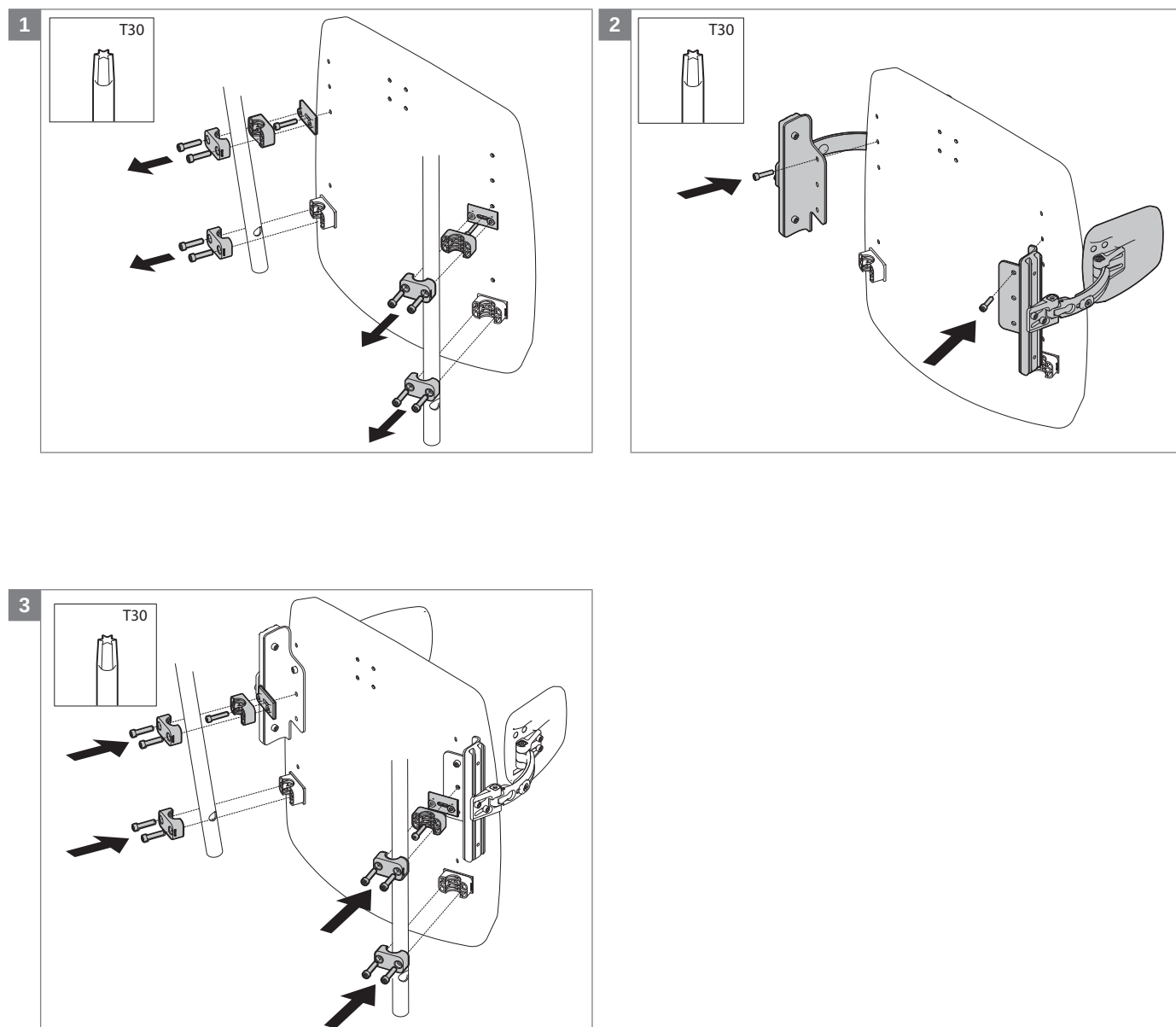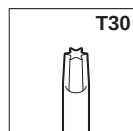

5.2.6 Mounting the standard legrest

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

- Mount the legrest mounting bracket to the seat support frame.
- Release the footplate fixing screw and lower the footplate to remove it from the legrest.
- Position the clamp nut in the legrest profile and attach the calf plate to the clamp nut.
- Adjust the calf plate into the desired position, tighten the clamp nut
- Reposition the footplate and tighten the fixing screw.
- Lower the legrest into the tube of the mounting bracket, holding the footplate outwards. Swing the legrest inwards to lock it in position.




Note! Make sure the locking pin falls in place (with a 'click') and secures the legrest from becoming detached.





Note! If locking fails: Adjust the screw from the locking mechanism.

Relevant article numbers

1003563 Standard legrest for SW 42 cm (size 15), right
 1003562 Standard legrest for SW 42 cm (size 15), left
 01500.9019 Standard legrest for SW 46 cm (size 17), right
 01500.9119 Standard legrest for SW 46 cm (size 17), left
 1003565 Standard legrest for SW 50 cm (size 19), right
 1003564 Standard legrest for SW 50 cm (size 19), left
 1003567 Standard legrest for SW 54 cm (size 21), right
 1003566 Standard legrest for SW 54 cm (size 21), left

Tools used

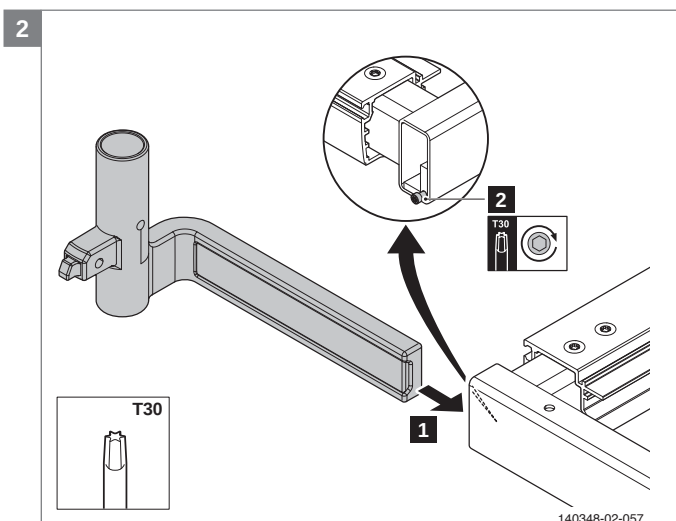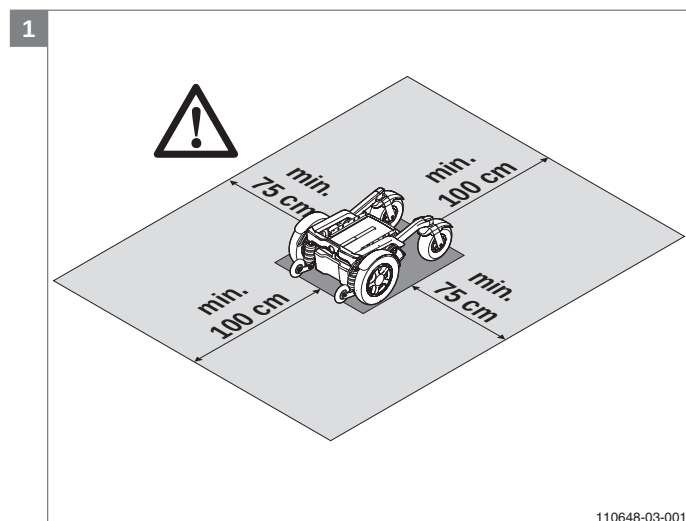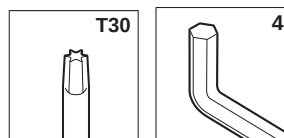

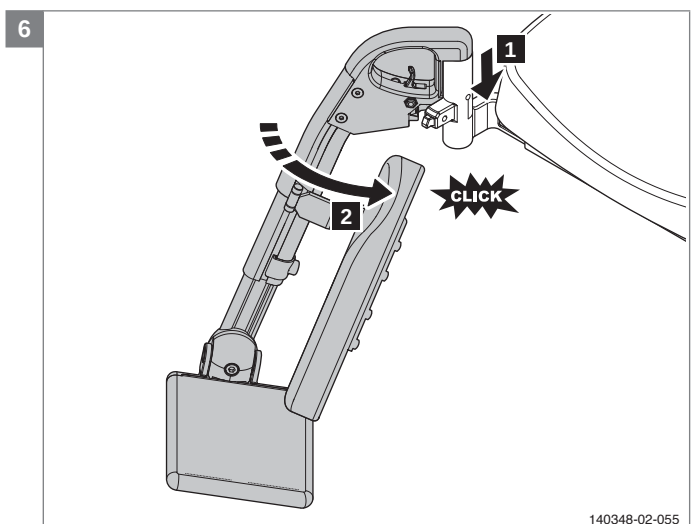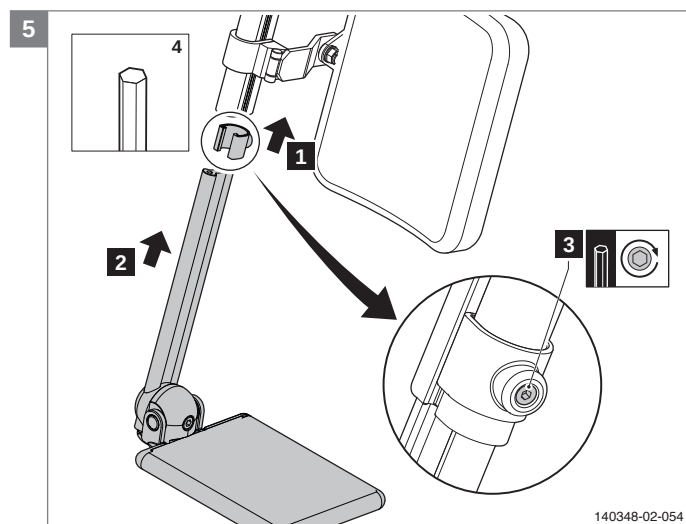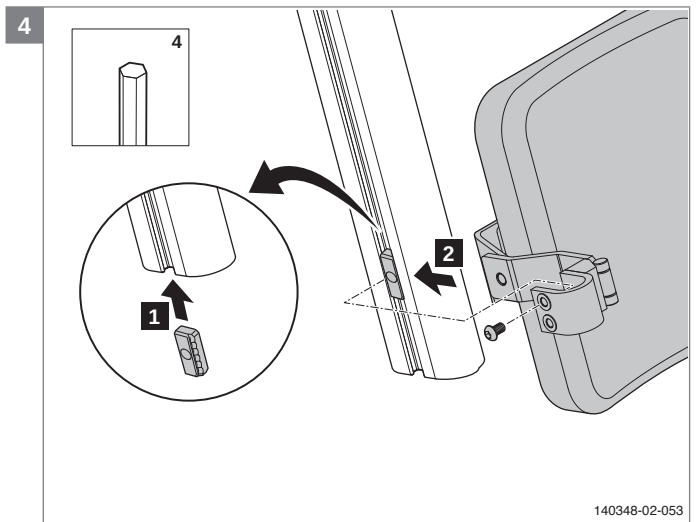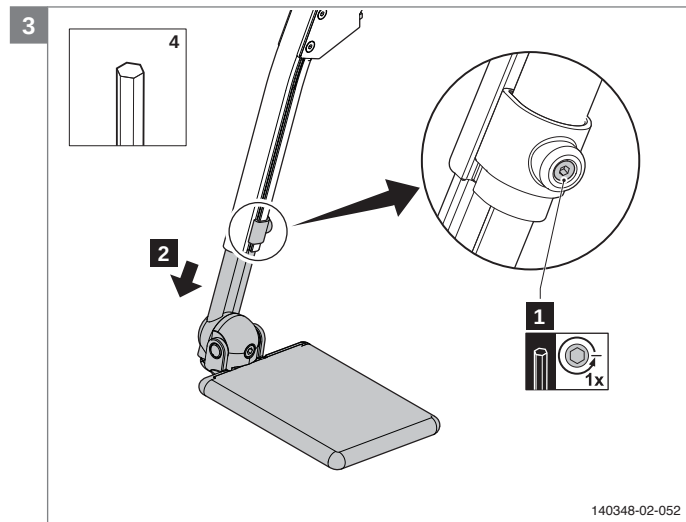

5.2.7 Mounting the manually elevating legrest

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

- Mount the legrest mounting bracket to the seat support frame.
- Secure the legrest in its upper position.
- Release the footplate fixing screw and lower the footplate to remove it from the legrest.
- Unscrew the lower clamp of the damper unit.
- Position the clamp nut in the legrest profile and attach the calf plate to the clamp nut.
- Adjust the calf plate into the desired position, tighten the clamp nut.
- Remount the lower clamp of the damper unit.
- Reposition the footplate and tighten the fixing screw.
- Lower the legrest into the tube of the mounting bracket, holding the footplate outwards. Swing the legrest inwards to lock it in position.




Note! Make sure the locking pin falls in place (with a 'click') and secures the legrest from becoming detached.




Note! If locking fails: Adjust the screw from the locking mechanism.

Relevant article numbers

- 140.00351.000 Manually elevating legrest for SW 42 cm (size 15), right
- 140.00321.000 Manually elevating legrest for SW 42 cm (size 15), left
- 01500.9219 Manually elevating legrest for SW 46 cm (size 17), right
- 01500.9319 Manually elevating legrest for SW 46 cm (size 17), left
- 140.00361.000 Manually elevating legrest for SW 50 cm (size 19), right
- 140.00331.000 Manually elevating legrest for SW 50 cm (size 19), left
- 1003601 Manually elevating legrest for SW 54 cm (size 21), right
- 1003598 Manually elevating legrest for SW 54 cm (size 21), left

Tools used

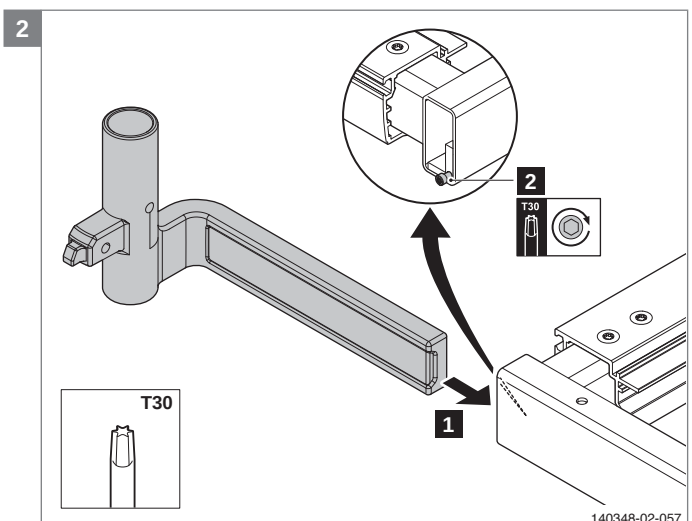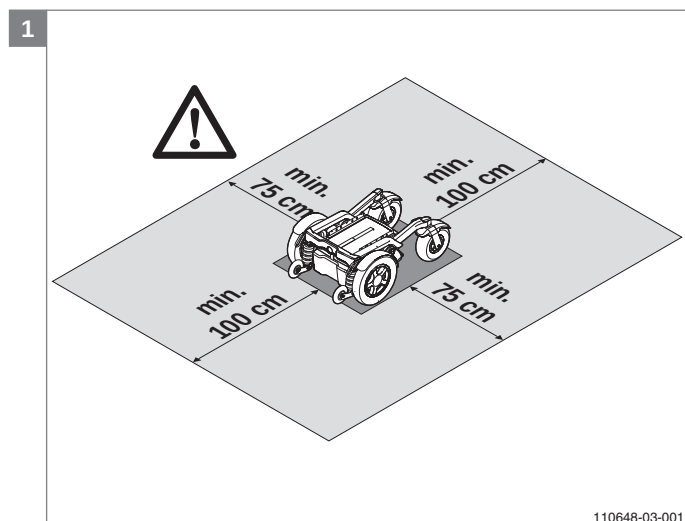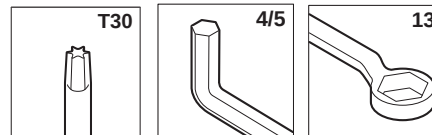

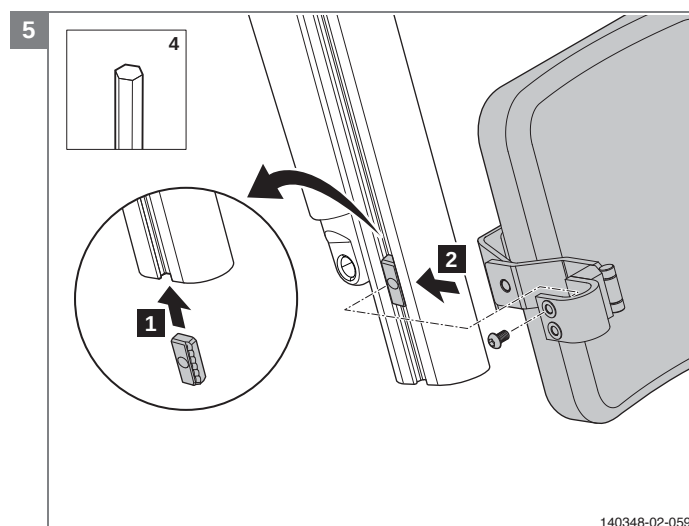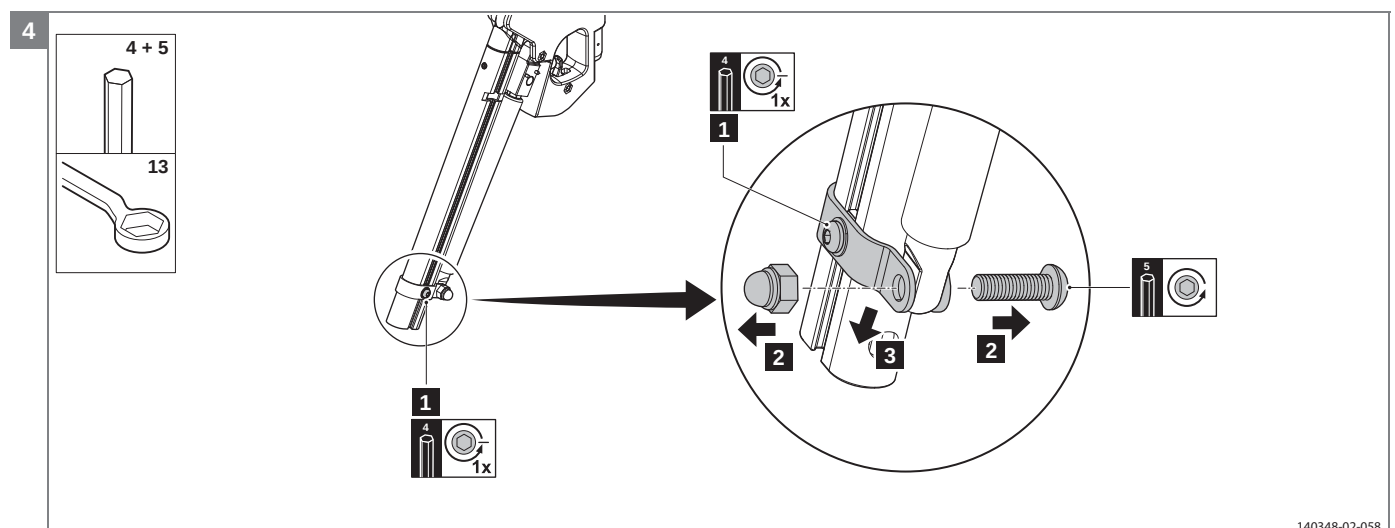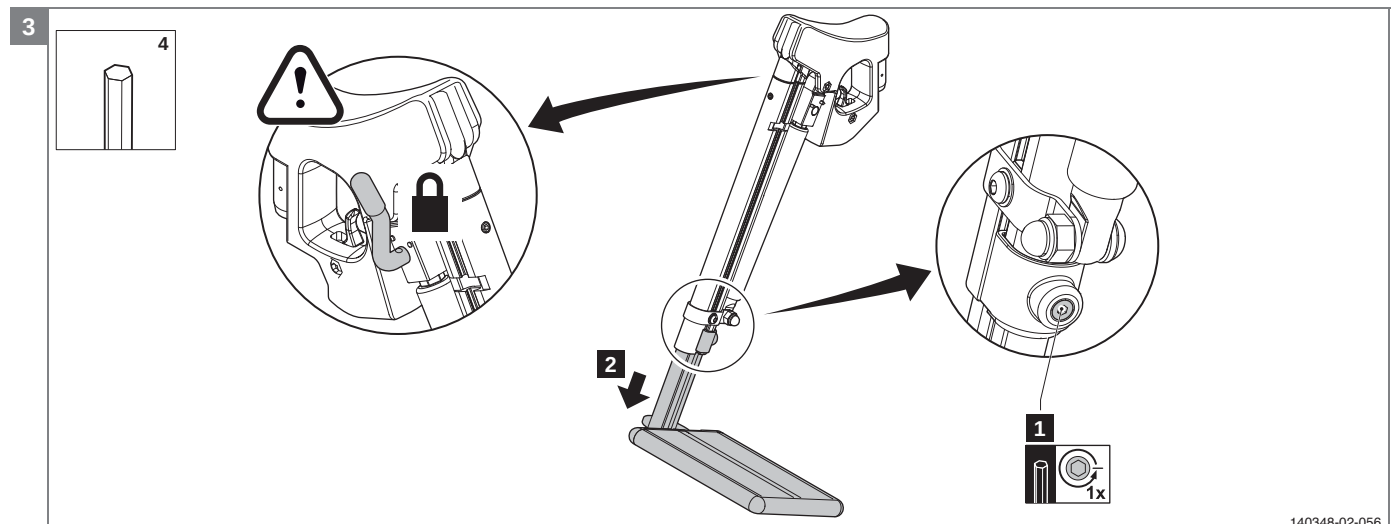

5.2.8 Mounting the electrically elevating legrest

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

- Remove seat and legrests.
- Place cable with plug from the top downwards into the tube from the legrest support. Secure the plug.
- Attach the connectors to the cables.
Note! Attach the brown cable to connector position with the (dot-) marking.
- Attach the cable extensions to the connectors and control box.
Note! See the backside of the control box cover for a wiring diagram.
- Release the footplate fixing screw and lower the footplate to remove it from legrest.
- Unscrew the lower clamp of the damper unit.
- Position the clamp nut in legrest profile and attach the calf plate to the clamp nut.
- Adjust the calf plate into the desired position, tighten the clamp nut.
- Remount the lower clamp of the damper unit.
- Reposition the footplate and tighten the fixing screw.
- Lower the legrest into the tube of the mounting bracket, holding the footplate outwards. Swing the legrest inwards to lock it in position.




Note! Make sure the locking pin falls in place (with a 'click') and secures the legrest from becoming detached.




Note! If locking fails: Adjust the screw from the locking mechanism.




Relevant article numbers

- 1003618 Sedeo Original electrically elevating legrest for SW 42 cm (size 15) right
- 1003613 Sedeo Original electrically elevating legrest for SW 42 cm (size 15) left
- 01502.9013 Sedeo Original electrically elevating legrest for SW 46 cm (size 17) right
- 01502.9113 Sedeo Original electrically elevating legrest for SW 46 cm (size 17) left
- 1003619 Sedeo Original electrically elevating legrest for SW 50 cm (size 19) right
- 1003614 Sedeo Original electrically elevating legrest for SW 50 cm (size 19) left
- 1003620 Sedeo Original electrically elevating legrest for SW 54 cm (size 21) right
- 1003615 Sedeo Original electrically elevating legrest for SW 54 cm (size 21) left

Tools used

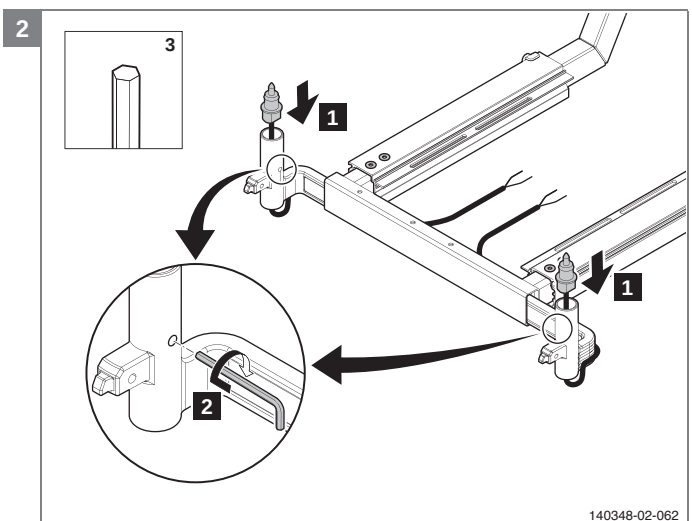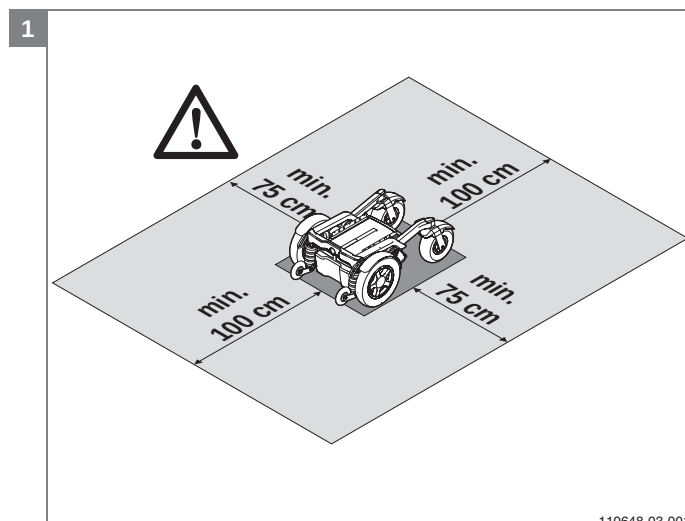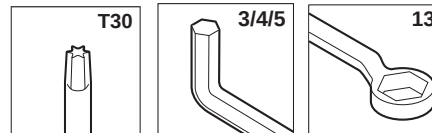

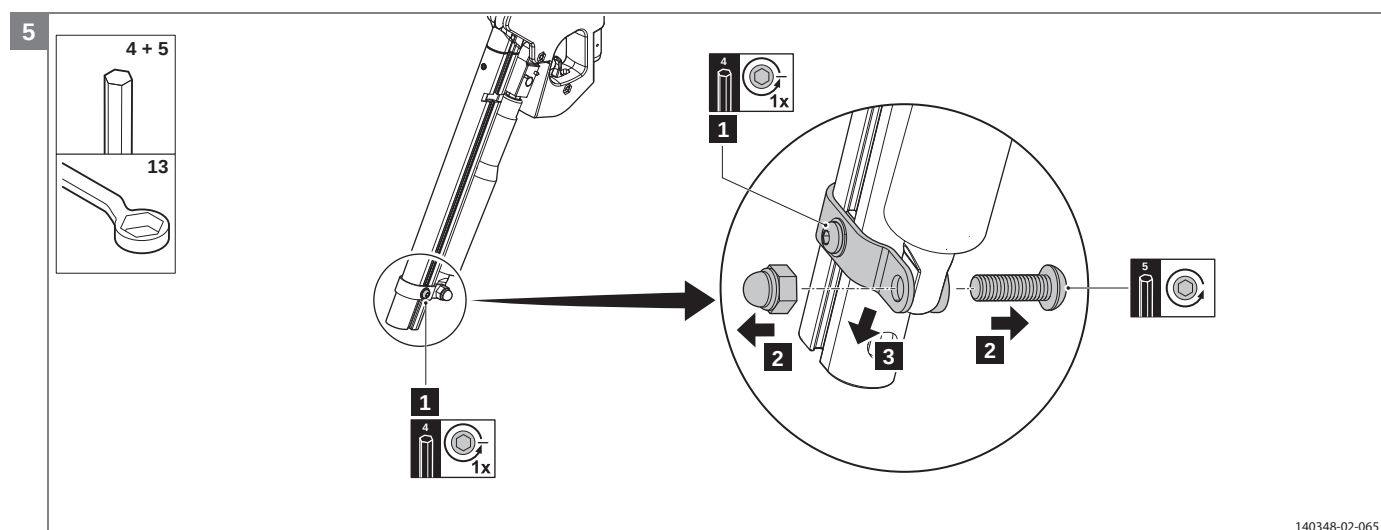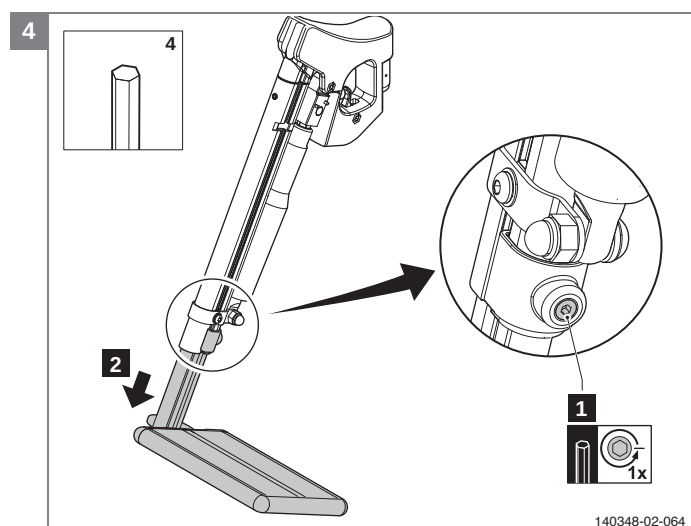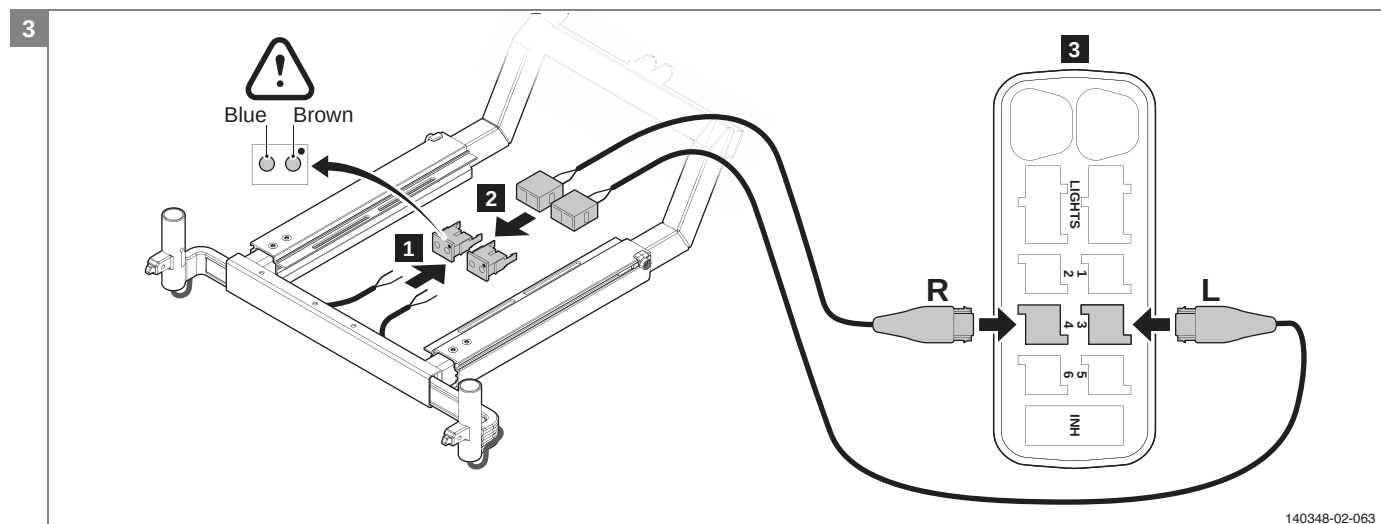

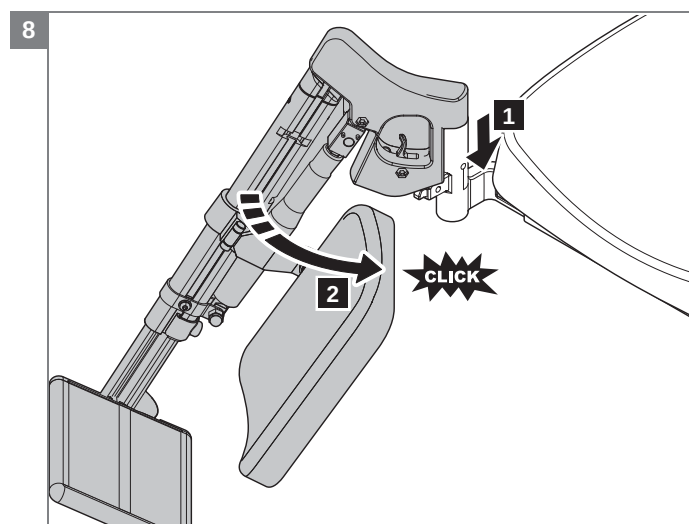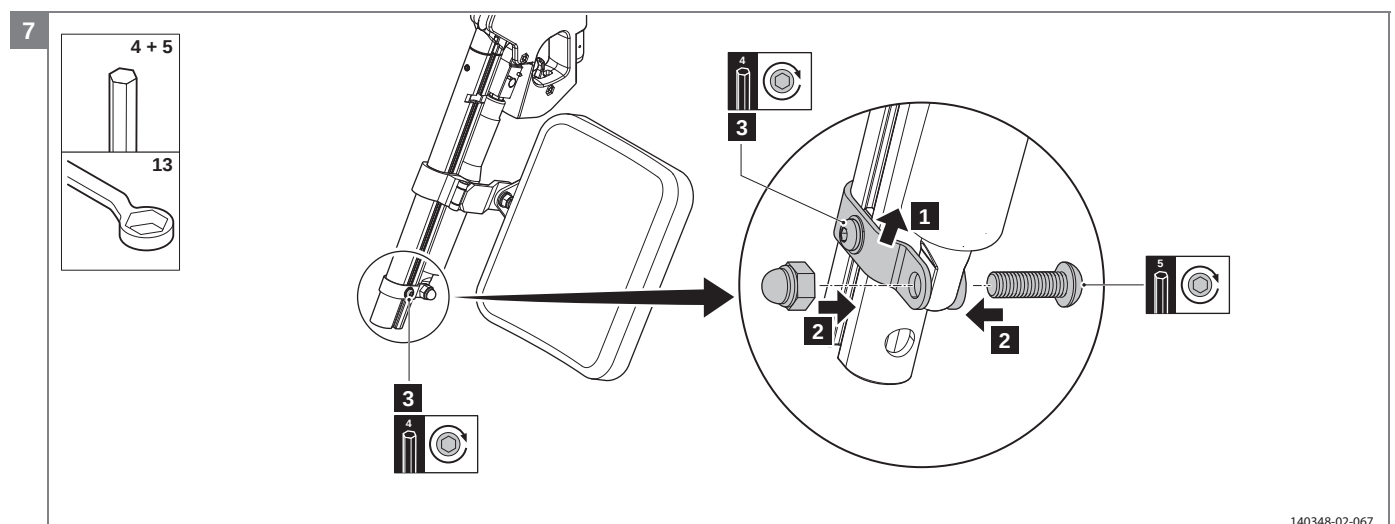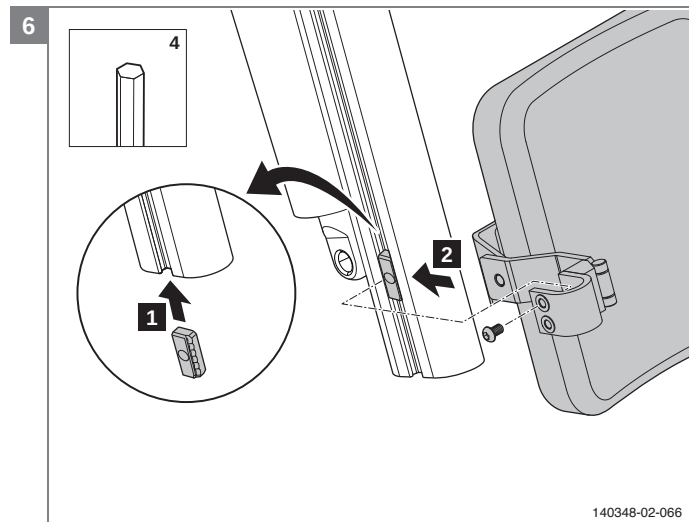

5.2.9 Mounting the central legrest

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

Adjust the level of friction of the footplate by tightening / unscrewing the Torx on the bottom of the footplate

Relevant article numbers





1011156 Central legrest

1011157 Central legrest with calf plate

1011158 Manually elevating central legrest

1011159 Manually elevating central legrest with calf plate

Tools used

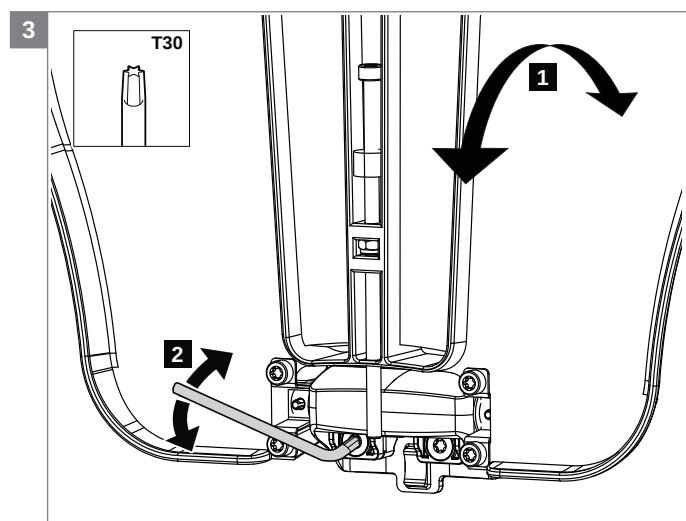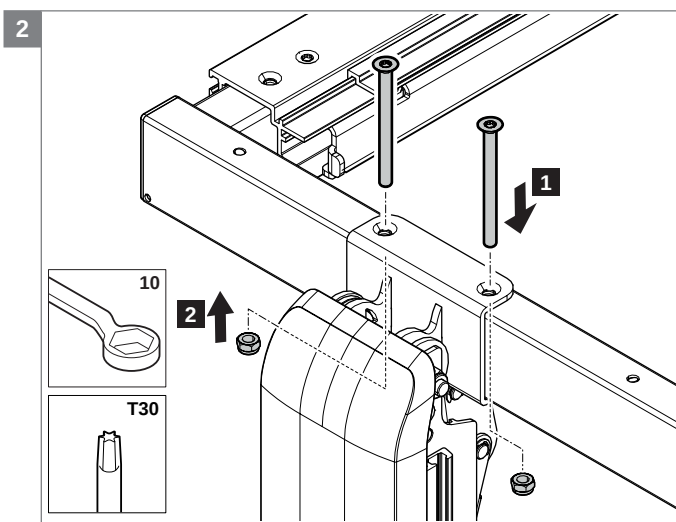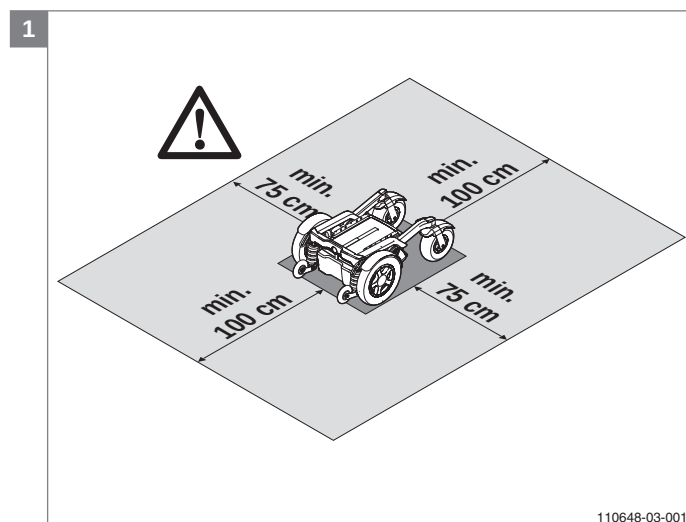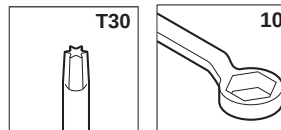

5.2.10 Mounting the electrically central legrest

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

- Remove seat and legrests.
- Unscrew the bolts (4x) that attach the legrest support frame to the seat support frame.
- Release the footplate fixing screw and lower the footplates to remove it from the legrest.
- Screw the calf plate mounting bracket to the central column.
- Mount the calf plates to the mounting bracket.
- Slide the central legrest unit into the front of the seat support frame and secure the central legrest unit with the bolts (4x) on the top and one bolt on either side of the seat support frame.
- Attach cable extensions to connectors and control box

Note! See backside of control box cover for wiring diagram.






- Hook the seat in place and push it downwards into position.
- **Caution!** Reprogram the controls of the chair.





Relevant article numbers

1016781 Electrically elevating central legrest without calfplates

1017106 Electrically elevating central legrest with length adjustment without calfplates

Tools used

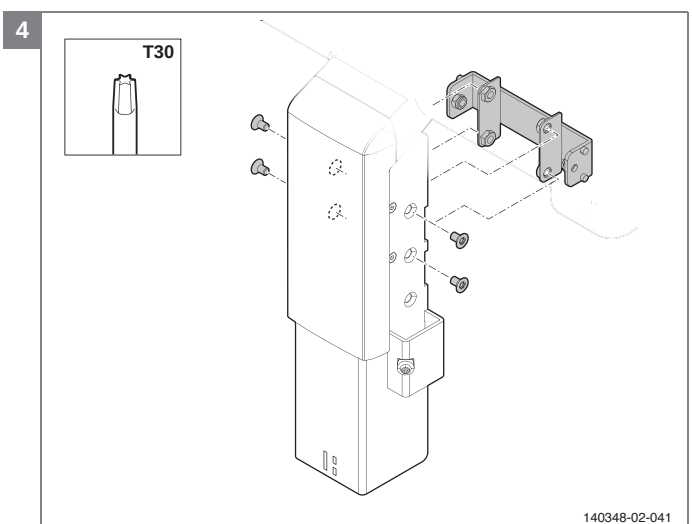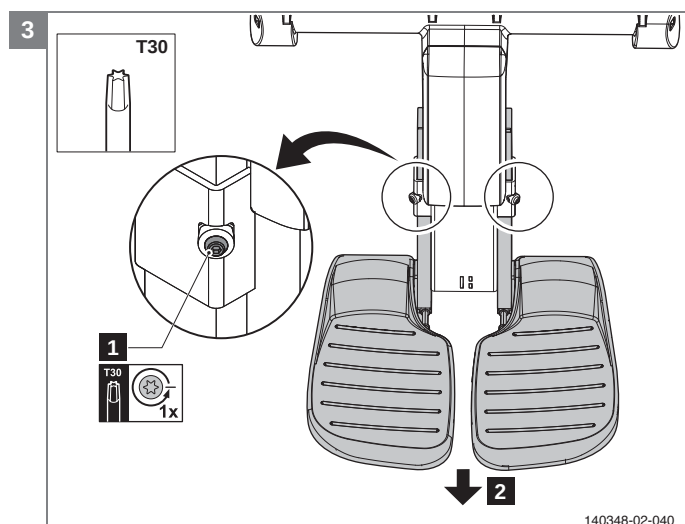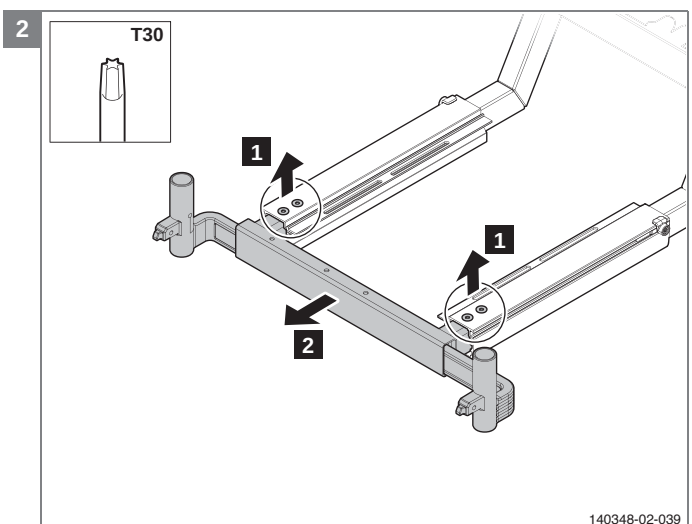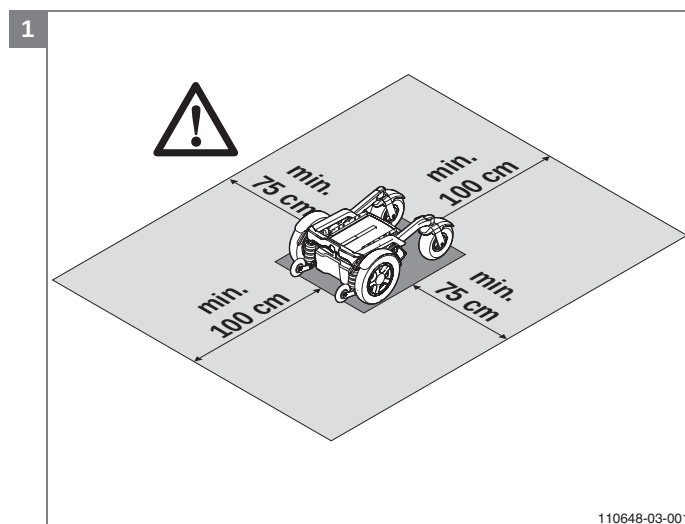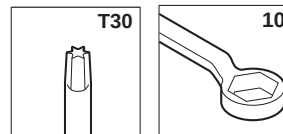

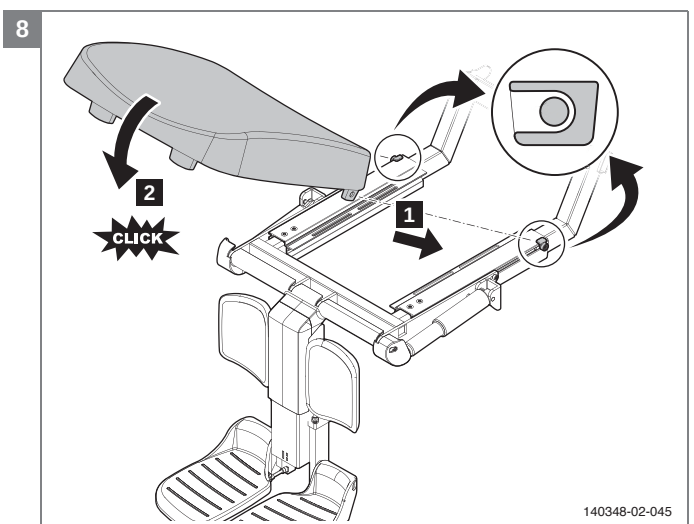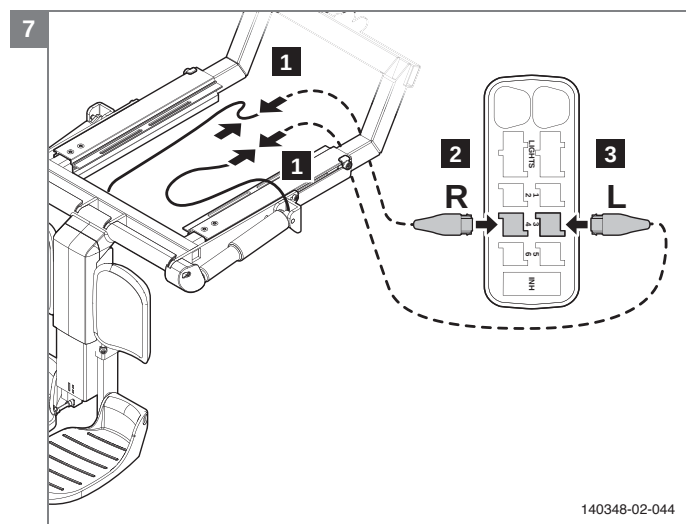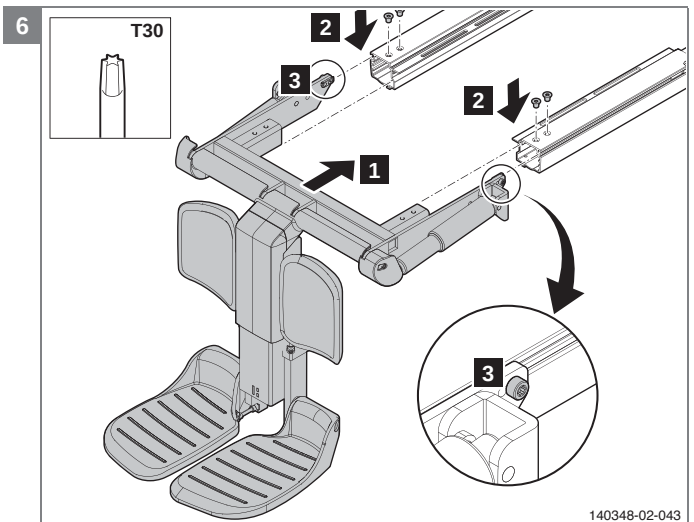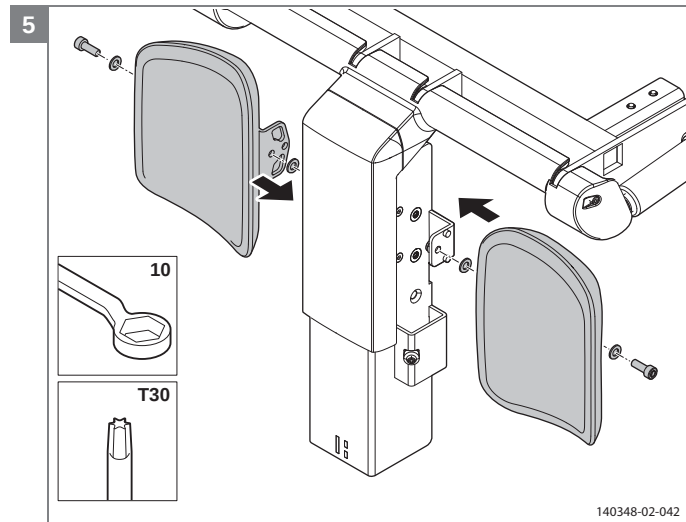

5.2.11 Mounting the standard or swing-away controller bracket

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

Note! Notice the position of the hexagon at any bracket. The hexagon should point with an edge forward!

Note! Replace the bolt in step 6 by a new precoted bolt

Note! Check the cable laying after you have mounted the swing-away controller bracket, connected the controller and adjusted the height, width and depth of the armrest. The controller cable needs enough clearance to swing away the controller freely inwards and outwards.

Relevant article numbers



1010951 Swing-away remote holder

1010952 Standard remote holder

1011164 Swing-away remote holder complete

1011165 Standard remote holder complete

Tools used

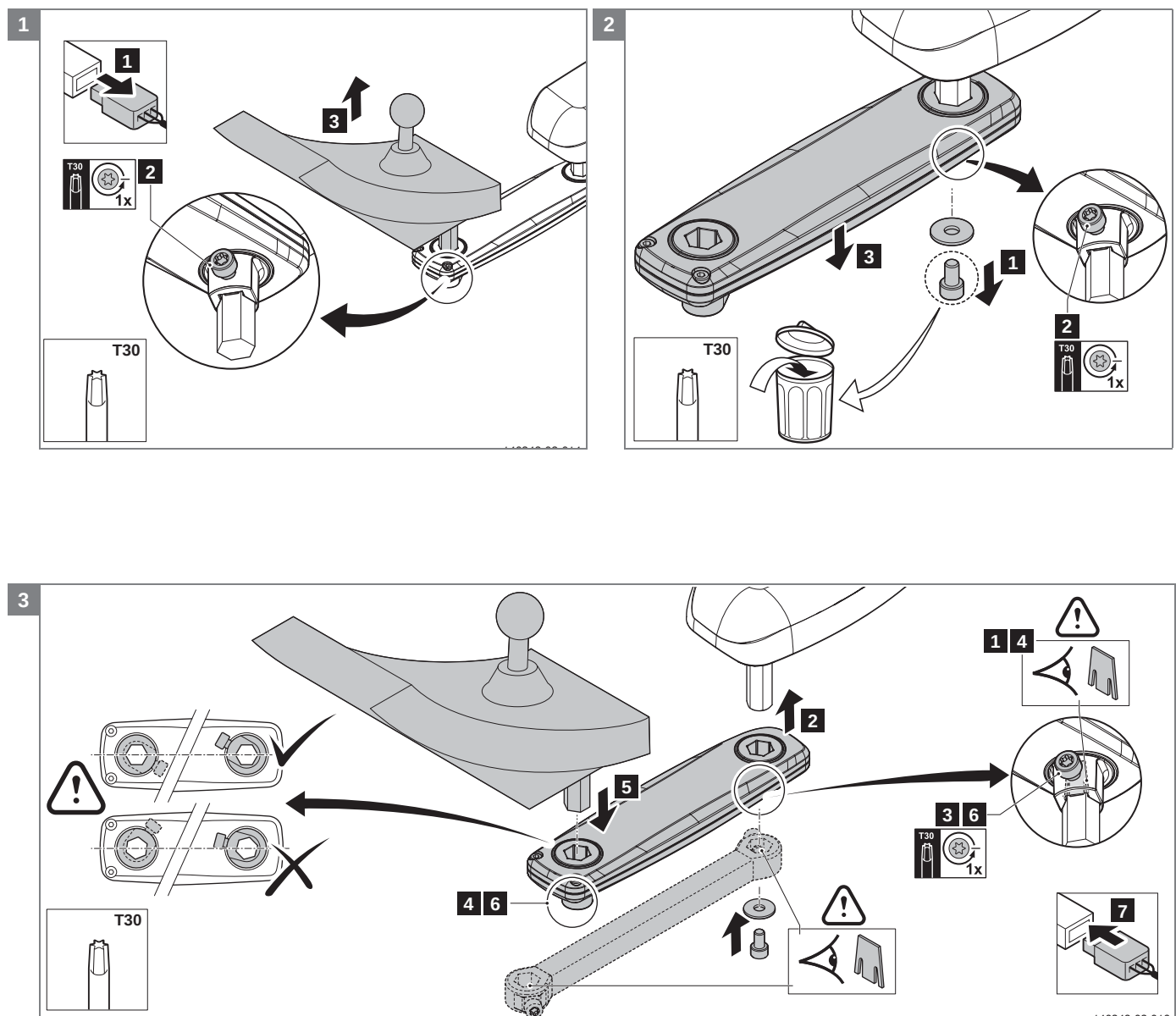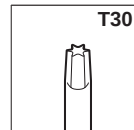

5.2.12 Mounting the lights on the seat frame

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!
Follow the illustrated instructions

Relevant article numbers

1016796 Lighting set R-net for mounting on seat frame




1016799 Front light module for mounting on seat frame, left
= right

1016801 Rear light module incl mounting bracket, left = right

9006331 Lighting harness R-net lighting

9006332 Lighting harness R-net lighting and horn

Tools used

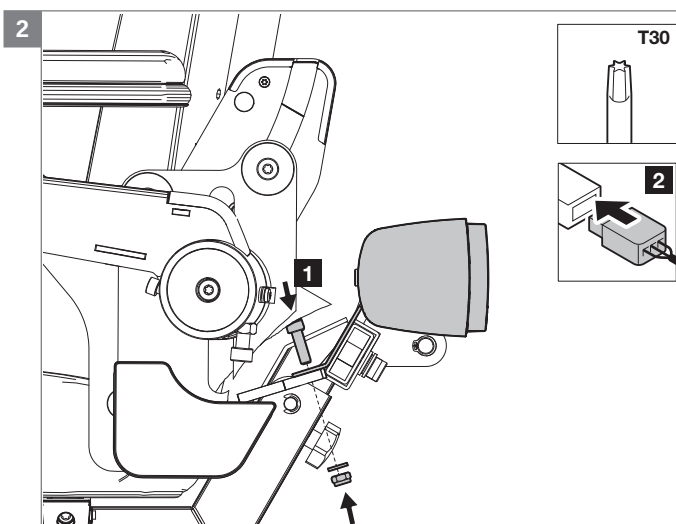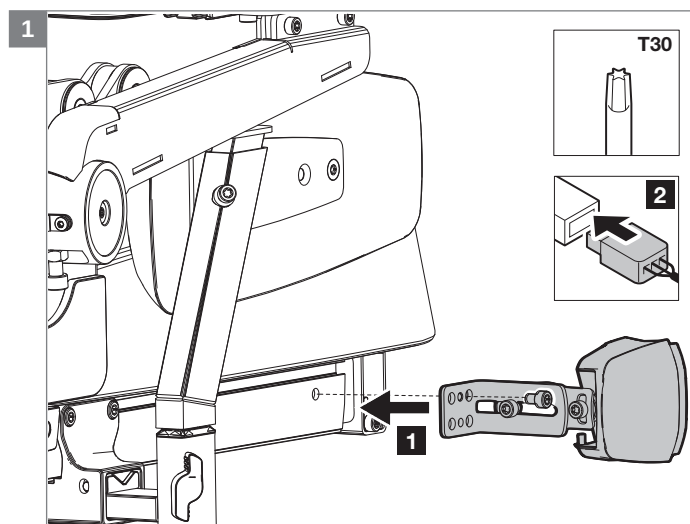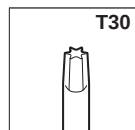

5.2.13 Mounting the lights on the central legrest frame

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

Follow the illustrated instructions

Relevant article numbers



1016795 Lighting set R-net for mounting with central electrical legrest

1016797 Front light module for mounting with central electrical legrest, left = right

1016801 Rear light module incl mounting bracket, left = right

9006331 Lighting harness R-net lighting

9006332 Lighting harness R-net lighting and horn

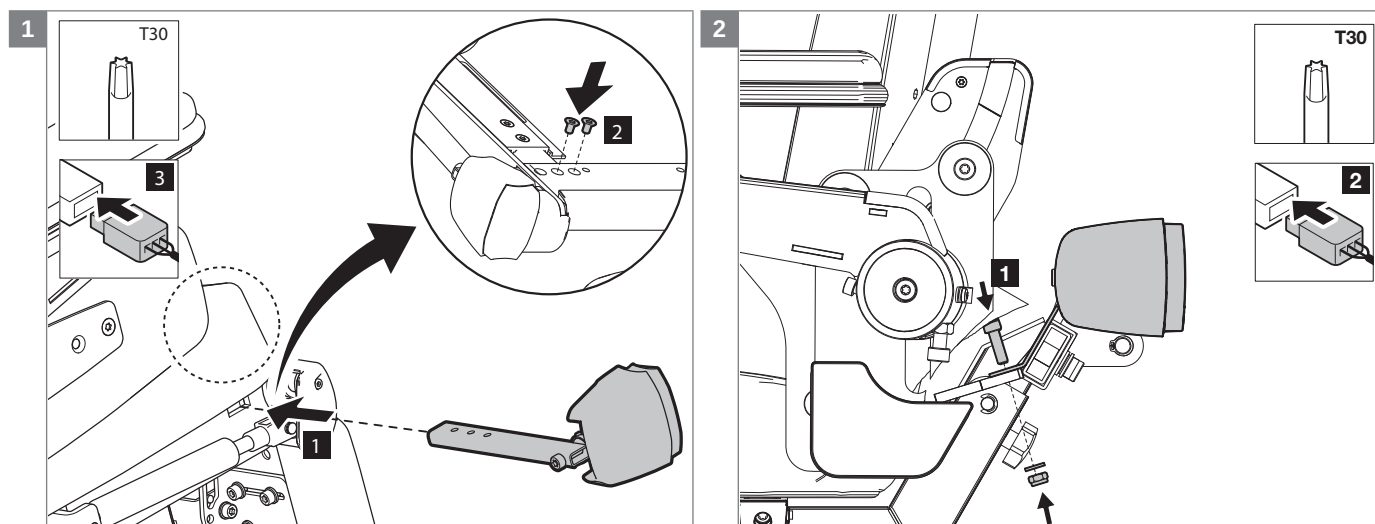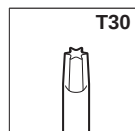Tools used

5.2.14 Mounting the side support

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!



Follow the illustrated instructions

Relevant article numbers

1017219 Side support complete, right

1017220 Side support complete, left

Tools used

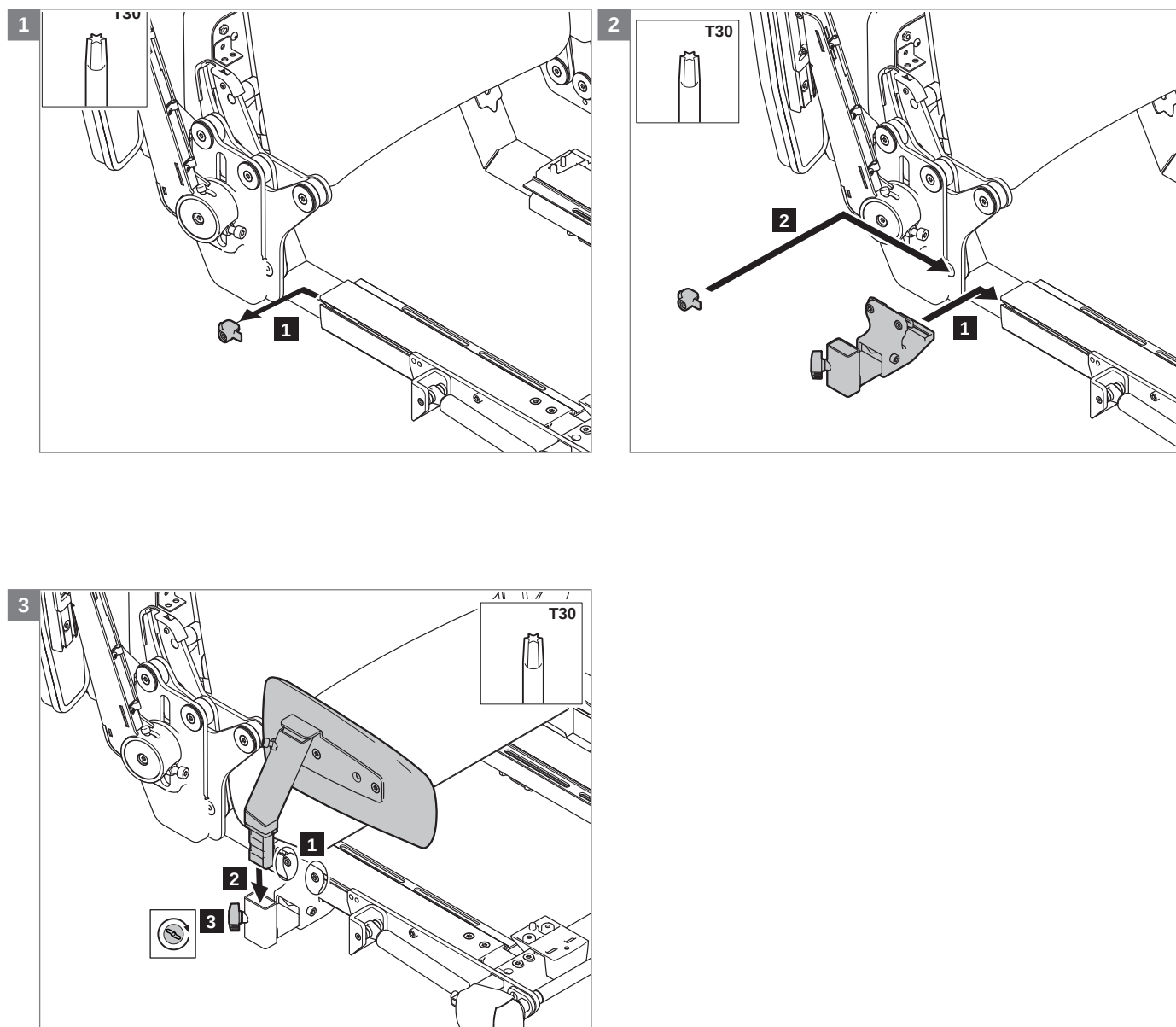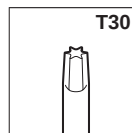

5.2.15 Mounting the abduction support

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!






Follow the illustrated instructions


- **Note!** Adjust the abduction support to the desired position, see the Sedeo Pro+ user manual.

Relevant article numbers

1011002 Abduction support

Tools used

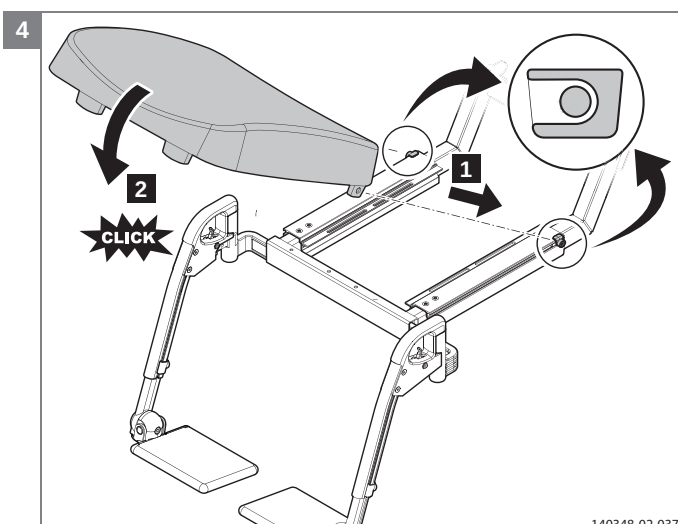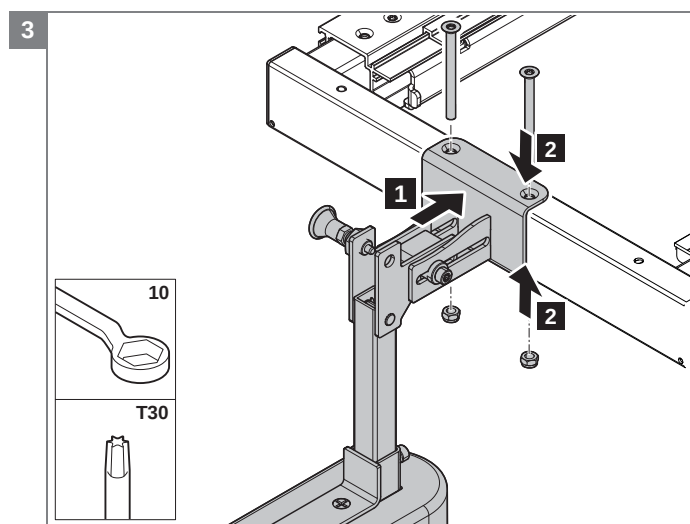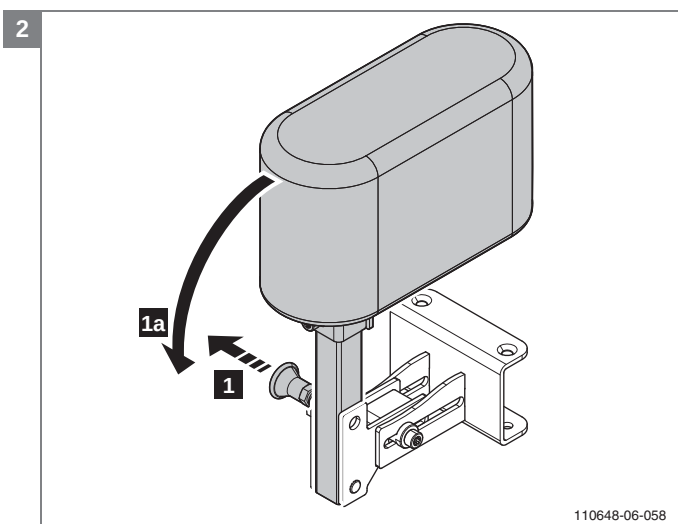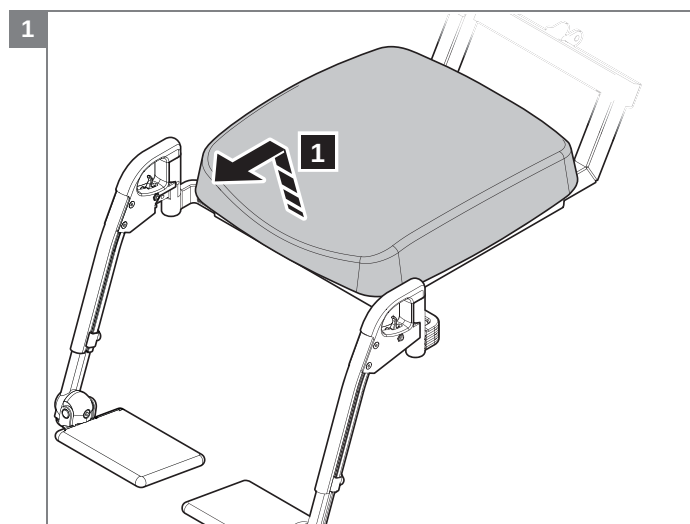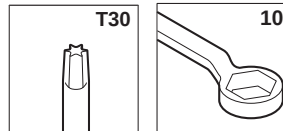

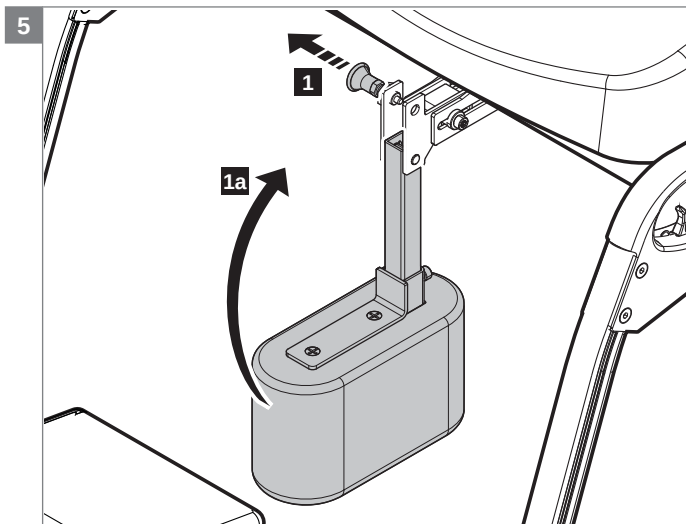

5.2.16 Mounting the armrest pads

Preparation

Note! Switch off the wheelchair via remote control

Instructions

Note! Suitable work area required!

- Place two clamp nuts on the top each armrest frame
Note! Do not fully tighten the clamp nuts.
- Slide mounting profile on the clamp nuts.
- Screw the armrest pads on the mounting profile.
- Adjust the armrest pads into the desired position.
- Tighten the clamp nuts.

Relevant article numbers

1017052 Armrest pad 30 cm upholstered, stretch rubber, right = left

1017053 Armrest pad 40 cm upholstered, stretch rubber, right = left





1010908 Armrest pad 30 cm, PU, right

1010909 Armrest pad 30 cm, PU, left

1017054 Armrest pad 30 cm removeable upholstery, fabric, right = left

1017055 Armrest pad 40 cm removeable upholstery, fabric, right = left

Tools used

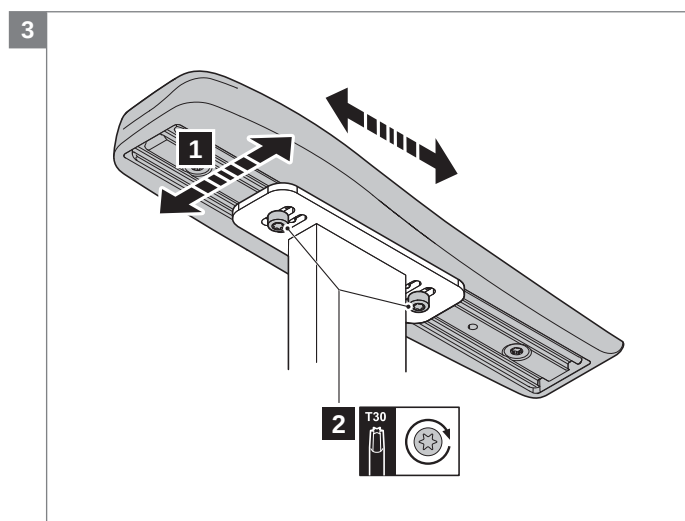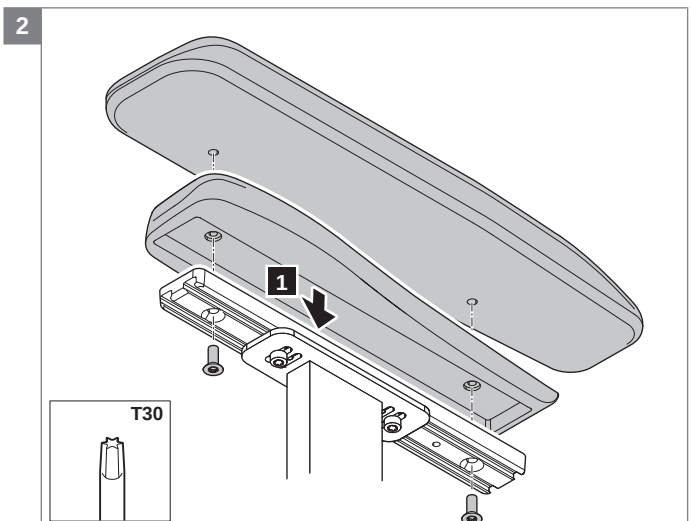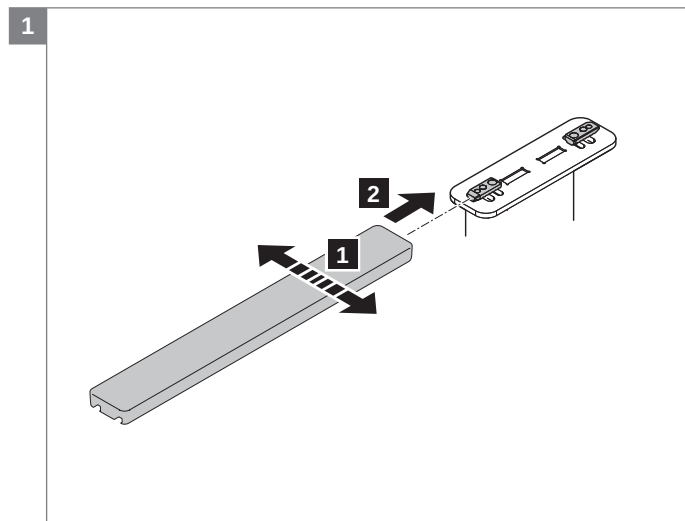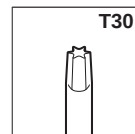

5.2.17 Mounting the hemiplegia armrest

Preparation

Note! Switch off the wheelchair via remote control.

Instructions






Note! Suitable work area required!




- Unscrew the armrest which is to be replaced from its mounting profile.
- Release the clamp nuts holding the mounting profile.
- Slide the mounting profile from the clamp nuts.
- Release the locking pin on the Hemiplegia rest and turn the pad sideways to expose the mounting bracket.
- Screw the mounting profile to the Hemiplegia mounting bracket.
- Slide the Hemiplegia armrest to the clamp nuts.
- Fasten the clamp nuts at the desired position.
- Swing back the Hemiplegia armrest and lock it in position.

Relevant article numbers

1010977 Hemiplegia armrest

Tools used

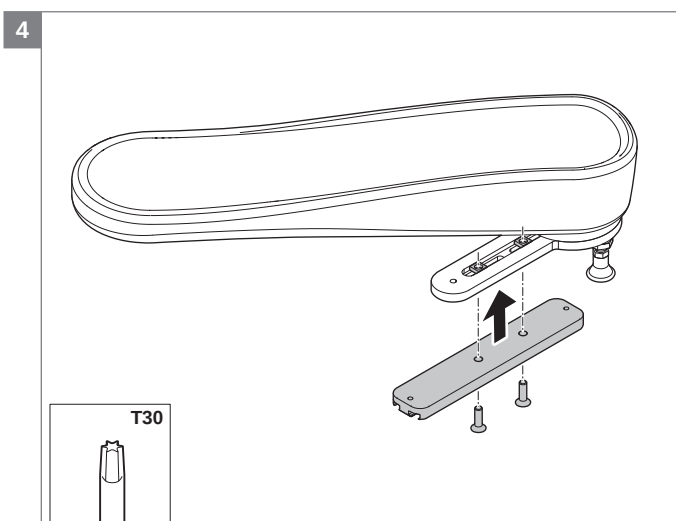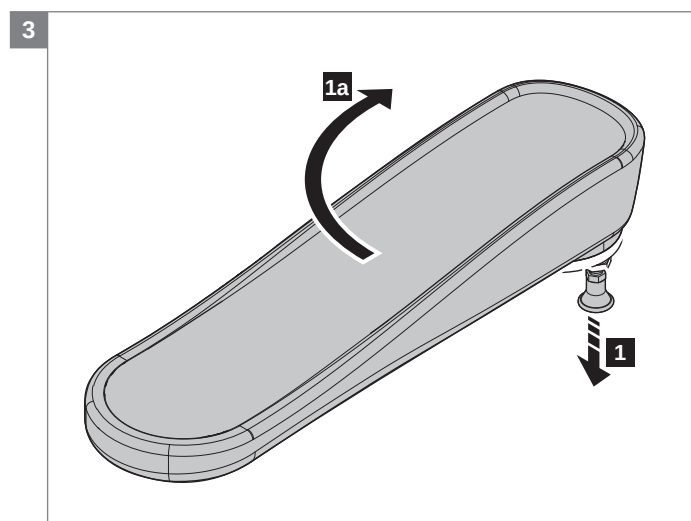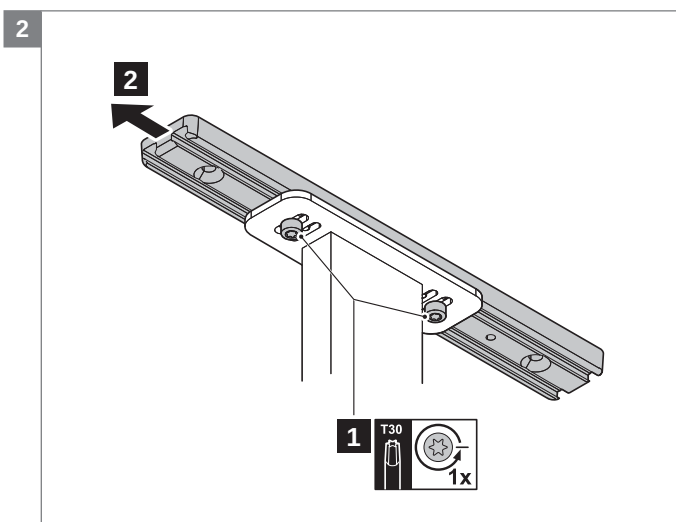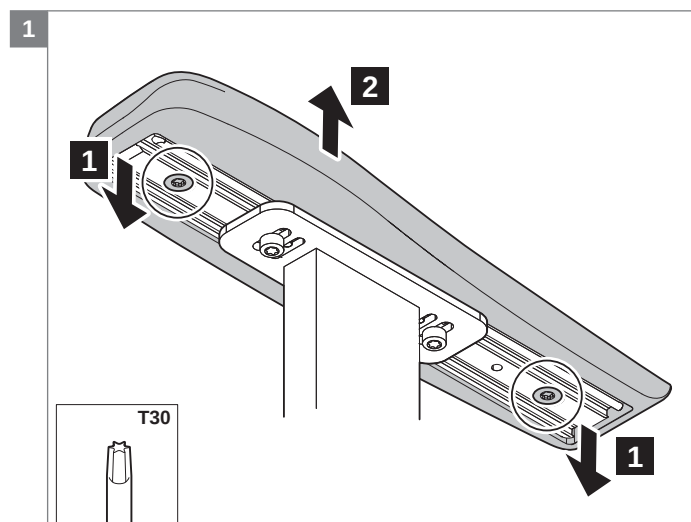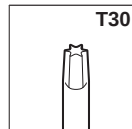

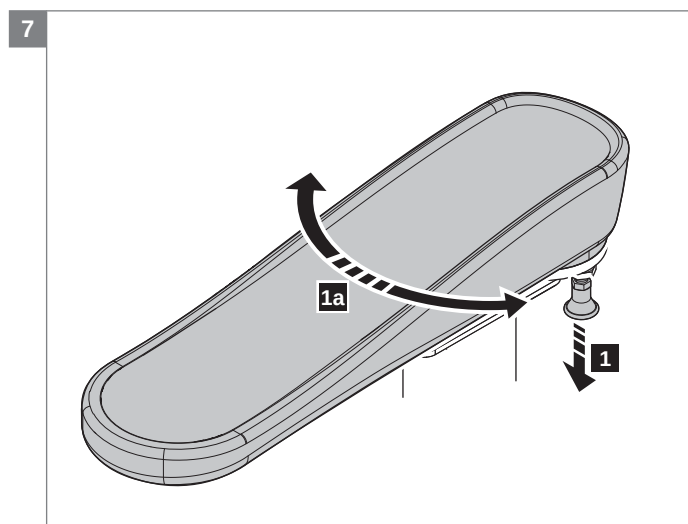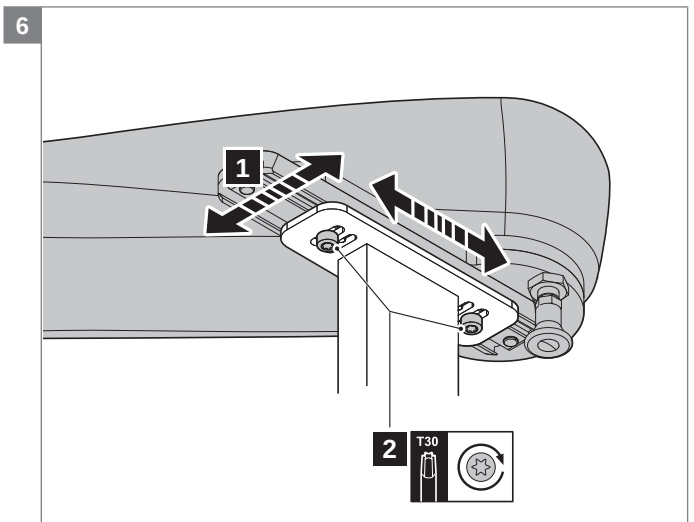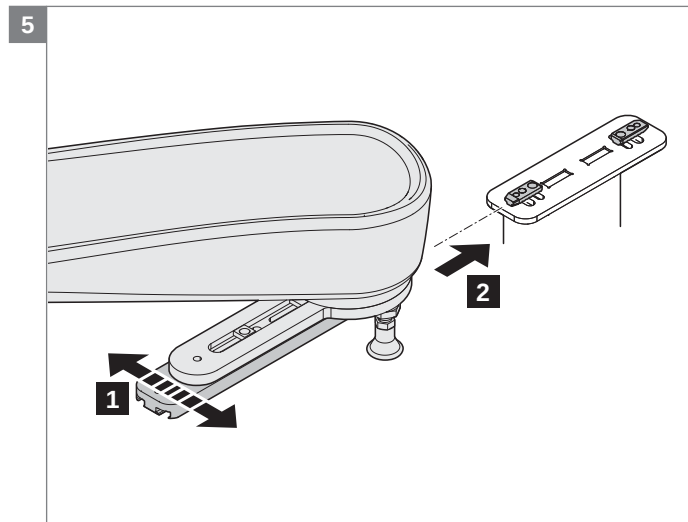

5.2.18 Mounting the mirror

Preparation

Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!

Follow the illustrated instructions

1-4 joystick side

5-8 not joystick side

- **Note!** Adjust mirror bracket when charger plug not in reach
- **Note!** Replace the bolt in step 7.6 by a new precoted bolt

Relevant article numbers



1012269 Mirror joystick side right





1012320 Mirror joystick side left

1012319 Mirror not joystick side right

1012321 Mirror not joystick side left

Tools used

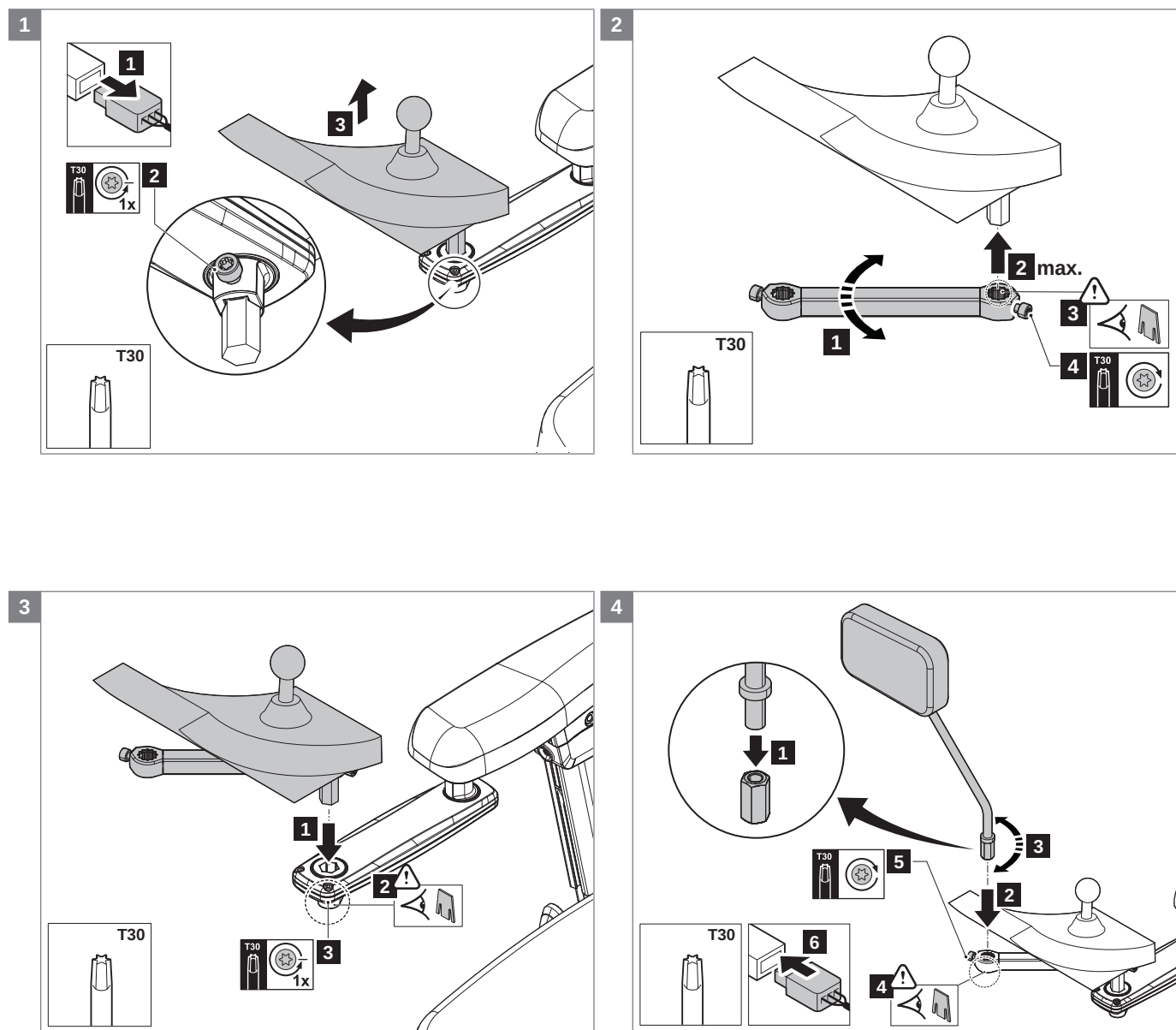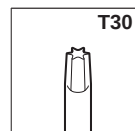

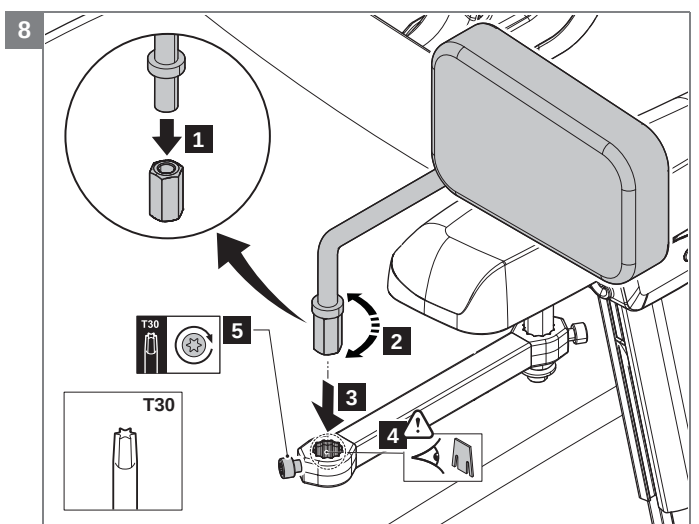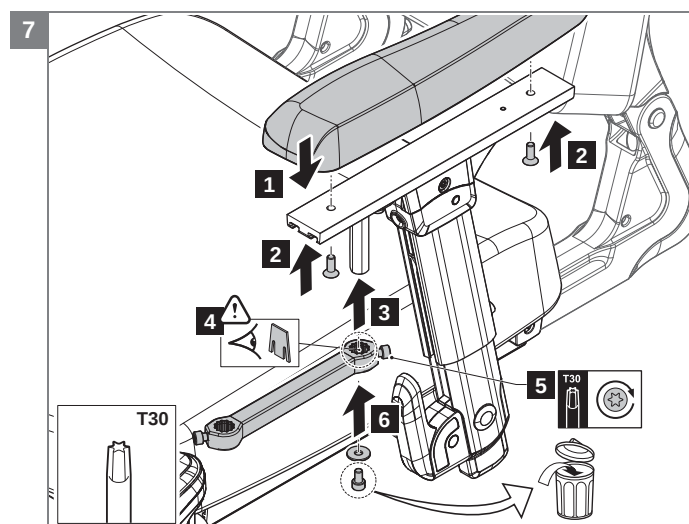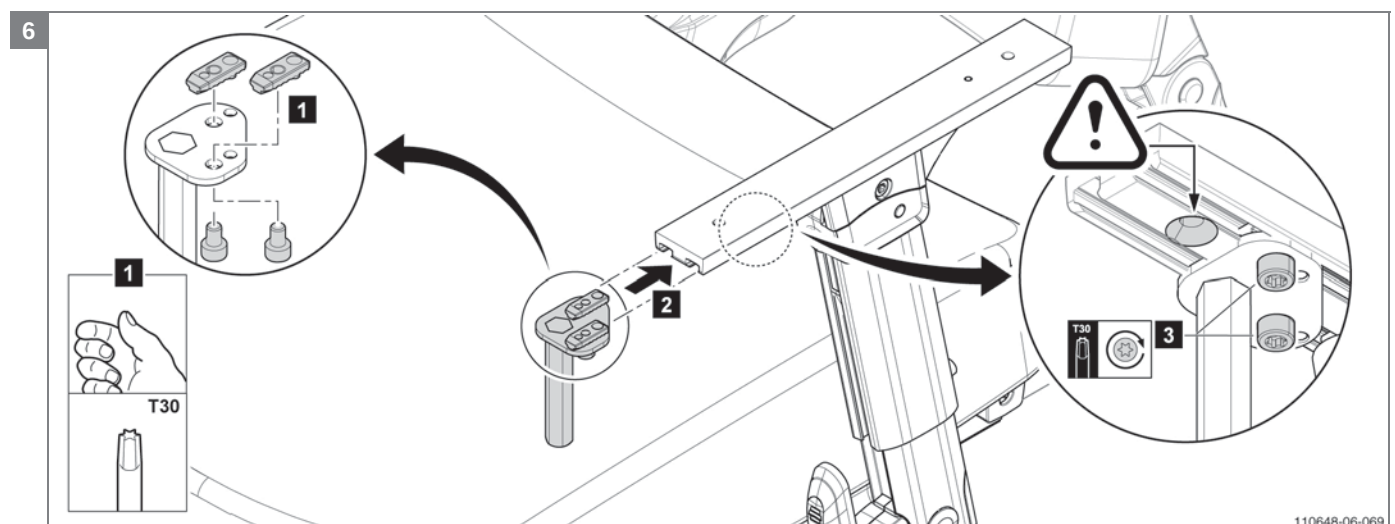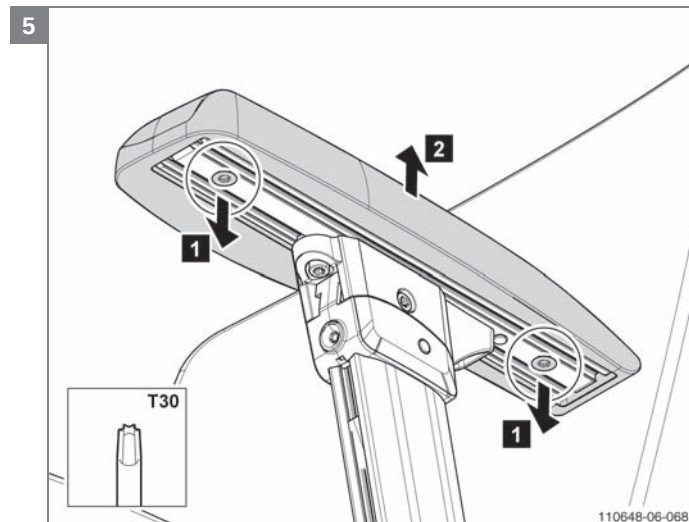

5.2.19 Mounting the positioning belt

Preparation

Note! Switch off the wheelchair via remote control.



Instructions

Note! Suitable work area required!
Follow the illustrated instructions

Relevant article numbers

1016819 Positioning belt short

1016820 Positioning belt long

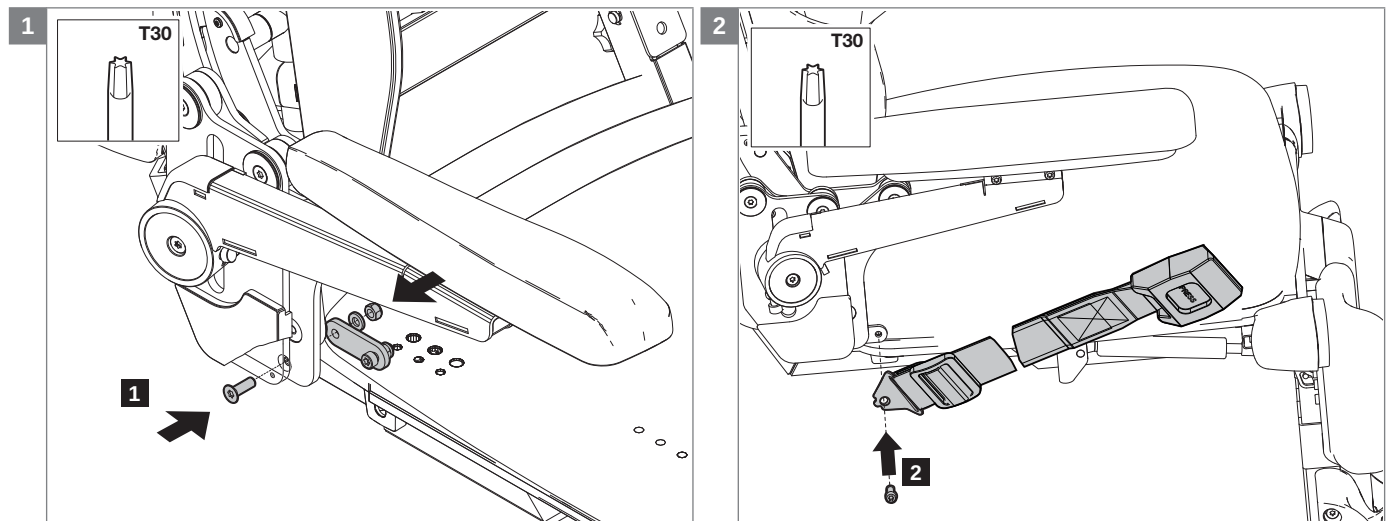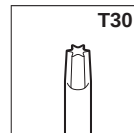Tools used

5.2.20 Mounting the cane holder

Preparation




Note! Switch off the wheelchair via remote control.

Instructions

Note! Suitable work area required!
Follow the illustrated instructions

Relevant article numbers

1016897 Cane holder

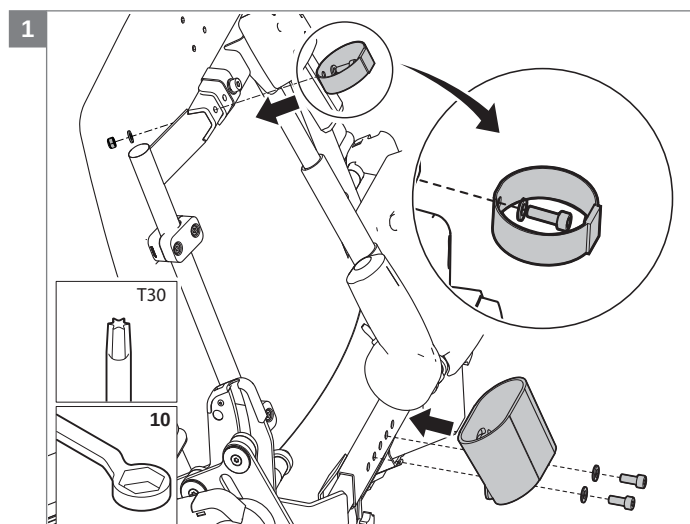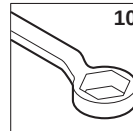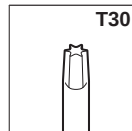Tools used

5.2.21 Mounting the table


Preparation

Note! Switch off the wheelchair via remote control.

Instructions




Note! Suitable work area required!
Follow the illustrated instructions

Tools used

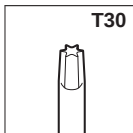

Relevant article numbers

- 1015167 Worktop, folding, acrylic + joystick recess right hand side
- 1015168 Worktop, folding, acrylic + joystick recess left hand side
- 1015169 Worktop with center remote recess for DX
- 1015455 Worktop with center remote recess for R-net
- 1015170 Half worktop, acrylic
- 1015162 Worktop, folding, acrylic with raised edges, suitable for upholstery + joystick recess right hand side
- 1015164 Worktop, folding, acrylic with raised edges, suitable for upholstery + joystick recess right left side
- 1015165 Worktop, folding, acrylic with raised edges, suitable for upholstery with centre remote recess

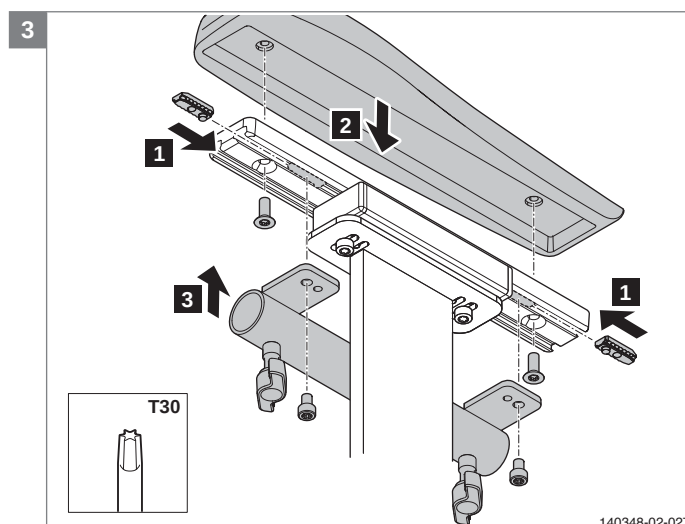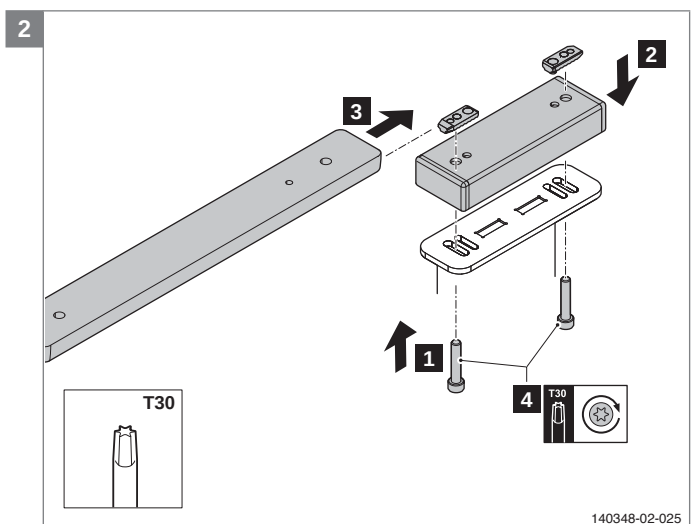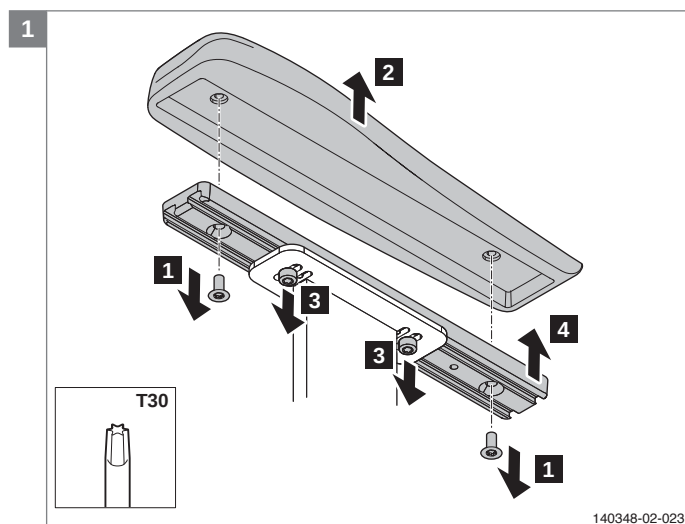

6 Trouble shooting

If the Sedeo seating system is not working or not working as it should, contact your supplier.

7 Specifications

7.1 CE declaration of conformity



The product is in conformity with the provisions of the Medical Devices Directive and carries the CE mark.

The product also meets the following requirements and standards. This has been verified by independent test organizations.

1. EN 12183 (2009) Manual seating systems or EN 12184 (2009) Electrically powered seating systems.
2. ISO 7176-8 Requirements for impact, static and fatigue strength.
3. ISO 7176-9 Climatic tests for electric seating systems and scooters.
4. ISO 7176-16 Requirements for resistance to ignition of upholstered parts.
5. ISO 7176-19 Crash test.
6. EN-ISO 10993-5 (1999-11) Biological evaluation of medical devices.
7. EN 1021 1,2 Assessment of the ignitability of upholstered furniture. Ignition source: smouldering cigarette/match flame equivalent.

7.2 Technical information Sedeo Pro+

Seating system	
User weight	160 kg
Seat width	380 - 600 mm
Seat depth	440 - 600 mm
Point of gravity adjustment range	95 mm
Backrest	
Backrest height	480 - 670 mm
Backrest recline	90 - 135°
Armrest	
Armrest pad width	-40 till +20 mm
Armrest pad height (from seat frame)	240 - 360 mm
Armrest pad depth controller bracket side	95 mm
Armrest pad depth non controller bracket side	140 mm
Side support	
Side cushion height	100 mm
Central legrest	
Lower leg length (with compressed seat cushion)	330 - 495 mm
Knee angle	30 - 110°
Ankle angle	+20° tot -40°
Central calf support height (option)	120 mm
Electrically elevating central legrest	
Lower leg length (with compressed seat cushion)	365 - 530 mm
Knee angle	94 - 170°
Max length compensation	180 mm
Ankle angle	+20 till -20°
Central calf support height	90 mm
Width calf supports	350 - 430 mm
Standard legrest (set)	
Lower leg length (with compressed seat cushion)	340 - 510 mm
Knee angle	70°
Ankle angle	0 - 180°
Elevating and electrical elevating legrest (set)	
Lower leg length (with compressed seat cushion)	340 - 530 mm
Knee angle	70° - 0°
Ankle angle	0 - 180°
Calf support height	120 mm
Headrest	
Headrest height (memory)	170 mm
Headrest depth	180 mm
Headrest vertical angle	-45° - 45°
Headrest width (option) standard headrest	-20 - 20 mm
Headrest width (option) width adjustable headrest	-110 - 110 mm

Controller bracket	
Controller bracket height	60 mm
Controller bracket depth	35 mm
Controller bracket swing-away	
Controller bracket height	50 mm
Controller bracket depth	35 mm
Controller position angle	0 - 45°
Hemiplegia	
Hemiplegia armrest angle (option)	0 - 45°
Abduction support	
Abduction support height	60 - 120 mm
Abduction support depth	60 - 110 mm
Thoracic support	
Thoracic support width	-60 - 0 mm
Thoracic support height	200 - 400 mm
Thoracic support depth	20 - 70 mm

8 Warranty

8.1 Definitions of terms

Definitions of terms used in this warranty:

- After sales service part: Part purchased after the initial product that is durable and may be subjected to natural wear and tear or natural contamination during normal operation within the lifetime of the product.;
- Consumable part: Part that is subjected to natural wear and tear or natural contamination during normal operation within the lifetime of the product (section 9 of Handicare's general terms and conditions of sale);
- Client: Those who purchase the product directly from Handicare;
- Corrective action: Repair, replace or refund of the product;
- Dealer: Those who re-sell the product to the User;
- Defect: Any circumstance due to which the product is not sound or fit to use, caused by a lack of quality of the material used to manufacture the product as well as the quality of the manufacturing process;
- Option: An accessory delivered with the initial product by Handicare to extend the standard product model;
- Product: Product that is delivered according to brochure or contract (e.g. wheelchair, scooter, battery-charger etc.);
- Part: Part of product that can be exchanged or replaced. This can be an option, accessory, service part or consumable part;
- Returns: Product or part that needs to be returned;
- RMA-process: Process to return goods, contact Handicare's Customer Service;
- User: Those who use the product;
- Warranty: The rights and obligations set forth in this document;
- Warranty period: The period of time during which the warranty is valid;
- Warranty provider: Handicare B.V., Vossenbeemd 104, 5705 CL Helmond, The Netherlands.

Notwithstanding the rights and obligations of Handicare, Client and User set forth in Handicare's general terms and conditions of sale, the rights of the Client and/or User towards Handicare in case of defects are limited to the provisions set forth in this warranty. For the duration of the warranty period Handicare guarantees that the product is without defects. In case of any defects the User is required –within two weeks after discovery of the defect- to contact the dealer. He has to complete a return form and return the product or part via the RMA-process. Handicare will, at its sole discretion, take the corrective action it seems fit under the given circumstances within a reasonable period of time (depends on nature of claim) from receipt of the completed return form. The warranty period will not be extended after a corrective action.

8.2 Warranty period table

Seating

Description	Warranty period	Examples include, but are not limited to the parts mentioned below
Frame	2 years	Weldment/frame
Electronics	1 year	Electronic components
After sales service parts	New: 1 year after invoice Repaired: 90 days after invoice	Metal parts
Consumable parts	40 days after invoice	Upholstery etc.
Options/Accessories	2 years	Lap strap, bag brackets etc. Delivered with the initial product

Handicare will only accept shipment costs and corrective costs related to warranty on equipment during the warranty period.

This warranty will void in case of:

- The product and/or its parts being modified;
- Changes in cosmetic appearance by use;
- Failure to observe the instructions for use and maintenance, use other than normal use, wear and tear, negligence, collateral damage by neglect of earlier symptoms, overloading, third-party accidents, non-original parts used and defects not caused by the product;
- Circumstances beyond our control (flood, fire, etc.).

This warranty does not cover:

- Tyres and inner tubes
- Batteries (covered by the battery manufacturer's warranty).

Clients and/or Users have legal (statutory) rights under applicable national laws relating to the sale of consumer products. This warranty does not affect statutory rights you may have nor those rights that cannot be excluded or limited, nor rights against the entity from whom the product was purchased. Clients may assert any rights they have at their sole discretion.



Handicare B.V.
Vossenbeemd 104
5705 CL Helmond
The Netherlands

T +31 (0)492 593 888
F +31 (0)492 537 931
www.handicare.com